

Electrostatic extraction of cold molecules from a cryogenic reservoir

L.D. van Buuren, C. Sommer, M. Motsch, S. Pohle, M. Schenk, J. Bayerl, P.W.H. Pinkse, and G. Rempe
Max-Planck-Institut für Quantenoptik, Hans-Kopfermann-Str. 1, D-85748 Garching, Germany

(Dated: May 9, 2019, PREPRINT)

We present a source which delivers a continuous, high-density beam of slow and internally cold polar molecules. In the source, warm molecules are injected into a cryogenic cell where their external and internal degrees of freedom are cooled by collisions with a helium buffer gas. Cold molecules are extracted out of the cryogenic environment by means of an electrostatic quadrupole guide. For ND_3 the source produces fluxes up to $(7 \pm 1) \times 10^{10}$ molecules/s with peak densities up to $(1.0 \pm 0.6) \times 10^9$ molecules/cm³. The cooling of the rotational state distribution has been monitored by depletion spectroscopy of H_2CO within the guide.

PACS numbers: 37.10.Mn, 37.20.+j, 33.80.-b

Keywords: cold molecules

Presently, there is a wide interest in cold molecules (see [1, 2, 3] and references therein). Cold polar molecules are employed in high-precision experiments and they are promising candidates for the implementation of quantum computation. Dense samples will allow the study of collisions and chemistry at low temperatures (see [4] and references therein). This would enable us, for example, to mimic conditions similar to the ones observed in interstellar clouds and explore formation rates of molecules there. By adding external fields, collisions and chemical rates can be manipulated, due to the Stark (Zeeman) interaction between electric (magnetic) fields and the molecular dipole (magnetic) moment. At sub-kelvin temperatures quantum effects will start dominating and resonances are expected to occur in the collision cross section. In general, these applications benefit from or demand large high-density samples of state-selected molecules. This letter reports on a general method to fulfill these requirements.

The complex internal structure of molecules precludes, in general, the technique of laser cooling. Therefore, several approaches are being pursued to produce large and dense samples of cold molecules. Indirectly, translational ultracold alkali and alkali-like molecules have been created from laser-cooled atoms using photo-association and Feshbach resonances (see [5] and references therein). Other, more general methods have been developed to bring molecules directly into the cold regime, where they can be trapped. Most of these methods are described in review articles [1, 2, 3]. Here we combine two separately developed techniques, buffer-gas cooling [6] and electric filtering of cold polar molecules [7], to produce a dense and state-purified beam of slow molecules. The benefit of combining buffer-gas cooling with magnetic filtering has already been demonstrated [8]. We extend this method to a wider variety of molecules by electrically extracting a continuous beam of polar molecules out of a buffer gas. Due to collisions with the buffer gas, the rotational temperature of the molecules is reduced. The resulting state population of the extracted molecular beam is studied by applying an optical depletion technique within our guide [9].

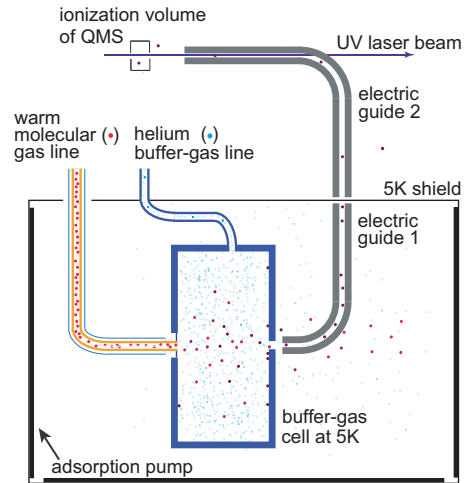


FIG. 1: (color online). Schematic view of the experiment. Warm molecules are cooled by collisions with a cold helium gas. Behind the cell an electric quadrupole extracts slow molecules in low-field-seeking states out of the cryogenic environment into the detection chamber, where they are detected by a QMS. Depletion measurements are performed by counter-propagating a UV laser beam through the guide.

In the source, warm molecules [10] are injected into a buffer-gas cell (see Fig. 1), where they are cooled by collisions with a helium gas at cryogenic temperatures. Due to these collisions not only the velocity distribution is compressed and shifted to lower values, but also the internal (rotational) temperature of the molecules is reduced as shown in a previous buffer-gas experiment [11]. Most of the molecules hit the cell wall where they stick. A fraction of the molecules exits the cell through a small aperture ($\varnothing = 1$ mm), after which the transversely slow molecules in low-field-seeking states are captured in an electric quadrupole guide similar to the one described in [7]. By bending the quadrupole guide (bending radius of 25 mm) a selection in longitudinal velocity is made. The guide transports the molecules to a separate vacuum system at room temperature where they are detected by a quadrupole mass spectrometer (QMS). The geometry of

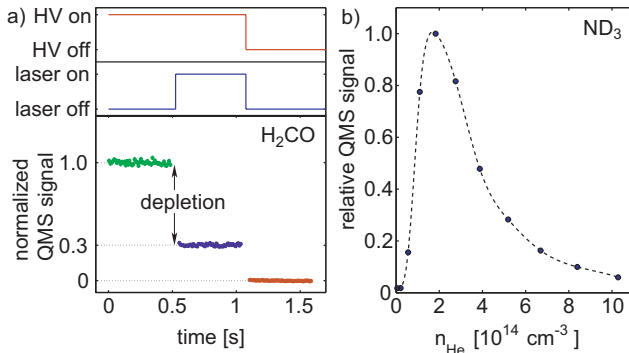


FIG. 2: (color online). a) Recording cycle to measure guided molecules and the depletion signal. The data are background subtracted and normalized to the average count rate without laser light. b) QMS signal (normalized to the maximal value) as function of n_{He} . For ND_3 , the highest signal is obtained with $n_{\text{He}} = 1.8 \times 10^{14} \text{ cm}^{-3}$. The distribution looks similar for H_2CO . The dashed curve is a guide to the eye.

this cross-beam-type QMS is such that it allows for free passage of a laser beam through the straight section of the guide. Information on the purity of the guided beam is obtained from depletion measurements, where we overlap the molecular beam in the guide with an UV laser beam tuned to a transition from a specific state [9]. The relative decrease in count rate observed with the QMS indicates the population of that state as shown in Fig. 2 a). Due to the internal cooling only a few low-lying states are populated causing large depletion signals.

The buffer-gas cell is cooled to cryogenic temperatures with a pulse-tube refrigerator (Cryomech PT 410) within 2 hours, starting from room temperature. To reach temperatures of $T \approx 5 \text{ K}$, the heat load of the system has to be limited to $P \approx 1 \text{ W}$. Two gold-coated copper shields are installed to block thermal radiation. The shields are bolted to the cold stages of the pulse-tube cooler. The cooler itself and all electric connections (including the guide and detection chamber) are mounted on the top flange of the vacuum vessel. This allowed us to easily adjust the experiment by lowering the vacuum vessel and radiation shields. The helium buffer gas is brought into the cell by a capillary, which is heat sunk to the 2 stages of the refrigerator ($T_1 \approx 35 \text{ K}$ and $T_2 \approx 5 \text{ K}$), whereas the molecular capillary is heated to avoid clogging. A 4 cm diameter and $12.5 \mu\text{m}$ thin Kapton foil thermally disconnects the warm injection line from the cold cell. Both helium and molecule flows are controlled by regulator valves. We typically flow $(2 - 30) \times 10^{15} \text{ He/s}$ and $(1 - 10) \times 10^{17} \text{ molecules/s}$ into the cell. The quadrupole guide is positioned at 1 mm distance from the exit aperture of the cell. Collisions between molecules trapped in the guide and background helium gas are minimized by a large area ($A \approx 10^3 \text{ cm}^2$) charcoal adsorption pump epoxied on the inner side of the second stage radiation shield. The quadrupole guide brings the molecules through openings in the top plates of both ra-

diation shields and through differential pumping sections into the detection chamber, where it ends $\sim 2.5 \text{ cm}$ before the QMS. The signal of guided molecules is the difference in count rate when switching on and off the high voltage (HV) in electric guide 1 (see Fig. 2 a)). The QMS response is calibrated by monitoring the increase in count rate and pressure when injecting small flows of molecules from the gas-handling system directly into the detection region.

The technique presented here works for all polar species which possess reasonably large Stark shifts, and which can be introduced into the buffer-gas cell. We have used deuterated ammonia (ND_3) to optimize the flux from the source. Previous reports showed successfully the buffer-gas cooling of NH_3 [12]. ND_3 has also been employed in electric guiding and filtering experiments [7]. Besides, it is easy to obtain and handle. To optimize the source we have varied the gas flows of both the helium and molecular gas and the cell length ($l_{\text{cell}} = 1, 2, \text{ and } 3 \text{ cm}$) at guide voltages of $\pm 5 \text{ kV}$, corresponding to an electric field of 93 kV/cm . The highest flux is expected with a high density of helium gas (n_{He}) in a short cell, such that the mean free path λ inside the buffer-gas cell is of order of the diameter of the exit aperture ($\varnothing = 1 \text{ mm}$) [8]. For $\lambda \ll \varnothing$ the molecules are accelerated by collisions with helium in close proximity of the exit aperture causing a smaller fraction to remain guided through the bend in the quadrupole, while for $\lambda \gg \varnothing$ only a very long cell will allow for enough collisions to slow down the molecules. But in that case most molecules will diffuse to the cell walls and stick there. Best results are obtained for $L_{\text{cell}} = 2 \text{ cm}$ and n_{He} of order 10^{14} cm^{-3} at which the guided ND_3 molecules have a measured average velocity of $\langle v \rangle \approx 65 \text{ m/s}$. A typical scan of the helium density in the buffer-gas cell is shown in Fig. 2 b). For optimal settings, the highest flux corresponds to $(7 \pm 4) \times 10^{10} \text{ molecules/s}$ in the guide, which includes corrections for the limited spatial acceptance of the QMS. This results in a peak density in the guide up to $(1.0 \pm 1.0) \times 10^9 \text{ molecules/cm}^3$ assuming a Gaussian beam profile with $\sigma = 400 \mu\text{m}$ [13], where σ represents the square root of the variance.

The flux coming out of the source decreases in time, due to the formation of ammonia ice near the inlet and/or exit apertures. Since the inlet is heated continuously, ammonia is not directly formed on the inlet capillary, but on the Kapton foil which connects the capillary to the cryogenic cell. In time an ice layer grows slowly on the foil which can finally obstruct the incoming molecular beam. The exit aperture is closing with ice more rapidly. To remove the ice, we have made the aperture in a $50 \mu\text{m}$ thick copper foil on which a Nichrome heater wire and Pt100 temperature sensor are glued. By heating, the frozen ammonia evaporates within minutes from the copper foil. We repeat this heating procedure every ~ 3 hours to keep fluxes above $\sim 50\%$ of the flux obtained in the beginning of the measurement. Typically, we stop cooling after a day of measurements (up to $\sim 10 \text{ h}$) to re-

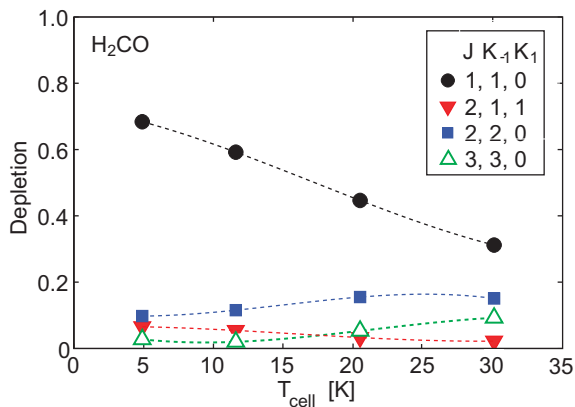


FIG. 3: (color online). Depletion as function of T_{cell} for different states. The dashed lines are guides to the eye.

move the helium gas from the charcoal and the molecules from the cell walls. By closing a valve between two guide sections separated by 1 mm, contamination of the ultra-high vacuum in the detection chamber is avoided.

A useful feature of our set up is that we can gradually adjust the cooling of the internal degrees of freedom. The amount of cooling depends on both the temperature of the buffer-gas cell (T_{cell}) and n_{He} . This allows us not only to follow the cooling process, but also to tune the source to other preferred internal temperatures. The lowest guidable state of ND_3 is expected to have a population of 96%, if we assume complete thermalization of the molecules to 5 K. However, it is a priori unclear how much the internal degrees have cooled inside the buffer-gas cell, since the source is optimized for maximum flux. The questions of state population and purity of guided molecules can be answered by laser depletion measurements as described for a similar setup in [9]. Since ND_3 has no easily accessible optical transition, we performed these measurements with formaldehyde (H_2CO). Switching these gases required only some minor changes in the gas handling system. The light at a wavelength of 330 nm is produced by second-harmonic generation of light from a ring dye laser. Upon excitation the molecules dissociate and are lost from the guide. The amount of depletion is given by the relative decrease in count rate observed by the QMS when the laser is switched on (see Fig. 2 a)). For these measurements the guide is switched to a voltage of $\pm 3\text{kV}$, which compromises between high flux and small line broadening due to the inhomogeneous electric field in the guide. The data shown here are taken at high laser power ($P_{\text{laser}} \approx 100\text{ mW}$) for which the depletion values reach $\sim 99\%$ of the high-power limit, determined by fitting the measured power dependence.

The depletion signal changes for all measured states as function of T_{cell} as shown in Fig. 3. Here, the states are labeled by $|J, K_{-1}, K_1\rangle$ in which J represents the rotational quantum number and K_{-1} (K_1) its projection along the molecular axis for the prolate (oblate) limiting case [14]. Due to varying electric field strength and

orientation along the laser beam, all guidable substates with different M -values, where M is the projection of the angular momentum on the electric field direction, contribute to the measured depletion signal. Highest populations are found for the $|1, 1, 0\rangle$ state. At $T_{\text{cell}} \approx 5\text{ K}$ no other contributing states than the ones shown were found. For the data shown in Fig. 3, the pressure in the buffer-gas line was kept constant at room temperature at the optimal settings found in Fig. 2, causing variations in n_{He} depending on T_{cell} . Still the increase of the $|2, 2, 0\rangle$ and $|3, 3, 0\rangle$ data (rotational energies of 40 cm^{-1} and 88 cm^{-1} , resp.) and the decrease of the $|1, 1, 0\rangle$ and $|2, 1, 1\rangle$ data (rotational energies of 11 cm^{-1} and 16 cm^{-1} , resp.) indicate clearly a dependence of the rotational temperature on the temperature of the buffer gas. Since the data do not indicate thermalization inside the buffer-gas cell, an increase in the $|1, 1, 0\rangle$ population might be feasible by increasing n_{He} at minimal T_{cell} , as discussed now.

Fig. 4 presents depletion measurements as function of n_{He} at $T_{\text{cell}} \approx 5\text{ K}$. Again, highest populations are found for the $|1, 1, 0\rangle$ state. At low densities, the population increases strongly with n_{He} , caused by the onset of rotational cooling of the molecules by collisions with cold helium gas. Also a decrease in the relatively high-energetic $|3, 3, 0\rangle$ state is visible. For $n_{\text{He}} \gtrsim 3 \times 10^{14}\text{ cm}^{-3}$ the populations converge to constant values. In this region, the total depletion signal sums up to 80-90%, illustrating the low number of contributing states to the guided beam. The remaining 10-20% is attributed to molecular orbits which do not interact with the laser beam and to molecules which do not dissociate, but decay back into guidable states. A similar fraction was found in our previous experiments performed at room temperature [9]. Taking this into account the population in the $|1, 1, 0\rangle$ state of H_2CO is $(82 \pm 10)\%$ for the optimal flux setting of n_{He} . From the depletion data we can also extract an entropic temperature, T_e , defined as the temperature of a thermal gas with an entropy S equal to the entropy in our beam ($S = -\sum p_i \ln p_i$), where p_i represents the measured population in state i . Due to the low number of contributing states, we find $T_e < 2\text{ K}$ for the optimal flux setting of n_{He} .

Note that we did not observe any population in the $|1, 0, 1\rangle$ state (rotational energy of 2.4 cm^{-1}) for which a population of $\sim 15\%$ is expected, assuming the guide extracts the molecules from a thermal ensemble at 5 K without any boosting effects at the exit aperture. A possible explanation for the absence of this state in the beam is given by its relatively low Stark shift. We know from velocity distributions [15] at high n_{He} , that the molecules can be accelerated when leaving the cell. Therefore, this Stark shift might be too low to keep the molecules guided in the quadrupole bend. At low n_{He} the velocity boost is small, but in this case data indicate higher rotational temperatures for which the population in $|1, 0, 1\rangle$ is decreasing rapidly.

In general, the addition of buffer-gas cooling to the electric filtering method purifies the extracted molec-

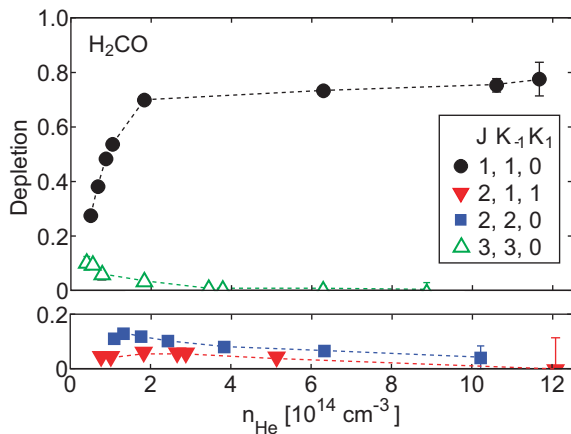


FIG. 4: (color online). Depletion as function of n_{He} for different states. The large statistical errorbars at high n_{He} are caused by the relatively low count rates at these settings. The lines are guides to the eye.

ular beam. This is especially beneficial for molecules with linear Stark shifts, because these molecules populate many guidable states at room temperature. For low-mass molecules with quadratic Stark shifts, such as water, the number of states contributing to an electrically filtered beam is much less than for the linear case [16]. Still the described source can enhance the purity of the extracted beam by adjusting the temperature of the buffer-gas cell to the rotational temperature of a low-lying state with large Stark shift. Currently, the source is being employed to produce beams consisting of polar molecules with masses up to ~ 100 amu, showing the generality of our method.

In summary, we have developed a source which extracts a continuous beam of slow molecules out of a cryogenic buffer gas and delivers them into a separate ultra-

high vacuum chamber. The source can operate for many hours with peak densities inside the quadrupole guide up to $(1.0 \pm_{0.6}^{1.0}) \times 10^9$ molecules/cm³ for ND₃. The internal temperature of the molecules is reduced by collisions with the helium gas. This leads to a strong increase of the purity of the beam. Comprising high flux and high purity, it is for H₂CO most beneficial to run the source just above the optimal flux setting of $n_{\text{He}} = 2 \times 10^{14}$ cm⁻³ (see Fig. 2 b)). If required, even higher purities can be obtained at the cost of lower flux or by depleting other contributing states. The extracted beam can be employed to load traps [17] in which collisions and chemistry at low temperature can be observed. Such experiments might also be performed by crossing separate beams of cold molecules. The high purity allows the study of state-dependent processes and facilitates optical detection schemes. The source can also be used to increase the resolution in high-precision experiments. Enhancement of the flux might be reached by pumping the (high-field seeking) ground-state population to a low-field seeking state in the region between the buffer-gas cell and the guide. Similar population transfer schemes in a molecular beam have been demonstrated (e.g. [18]). Finally, extraction of high-field-seeking states could be achieved by AC-switching the voltages of the guide electrodes [19, 20]. This would, for example, enable the production of a high flux, guided beam of slow ground-state molecules.

Acknowledgments

Support by EuroQUAM (Cavity-Mediated Molecular Cooling) and by the Deutsche Forschungsgemeinschaft through the excellence cluster "Munich Centre for Advanced Photonics" is acknowledged.

-
- [1] J. Doyle *et al.*, Eur. Phys. J. D **31**, 149 (2004).
 - [2] H. L. Bethlem and G. Meijer, Int. Rev. Phys. Chem. **22**, 73 (2003).
 - [3] O. Dulieu, M. Raoult, and E. Tiemann, J. Phys. B **39** (2006).
 - [4] R. V. Krems, Int. Rev. Phys. Chem. **24**, 99 (2005).
 - [5] J. M. Hutson and P. Soldan, Int. Rev. Phys. Chem. **25**, 497 (2006).
 - [6] J. Weinstein *et al.*, Nature (London) **395**, 148 (1998).
 - [7] T. Junglen *et al.*, Eur. Phys. J. D **31**, 365 (2004).
 - [8] D. Patterson and J. M. Doyle, J. Chem. Phys. **126**, 154307 (2007).
 - [9] M. Motsch *et al.*, Phys. Rev. A **76**, 061402(R) (2007).
 - [10] The optimal temperature of the capillary depends on the vapour pressure of the specific molecules; here we used temperatures in the range of $T \approx 180$ -295 K.
 - [11] D. Egorov *et al.*, Eur. Phys. J. D **31**, 307 (2004).
 - [12] D. R. Willey *et al.*, J. Chem. Phys. **120**, 129 (2004).
 - [13] T. Junglen, Ph.D. thesis, Max-Planck-Institut für Quantenoptik (2006).
 - [14] C. H. Townes and A. L. Schawlow, *Microwave Spectroscopy* (Dover Publications, Inc., 1975).
 - [15] Velocity distributions are obtained by differentiating the signal rise when the guide is switched on [7].
 - [16] T. Rieger *et al.*, Phys. Rev. A **73**, 061402(R) (2006).
 - [17] T. Rieger *et al.*, Phys. Rev. Lett. **95**, 173002 (2005).
 - [18] S. van de Meerakker *et al.*, Phys. Rev. A **68**, 32508 (2003).
 - [19] T. Junglen *et al.*, Phys. Rev. Lett. **92**, 223001 (2004).
 - [20] F. Filsinger *et al.*, Phys. Rev. Lett. **100**, 133003 (2008).