

# A Xenon Condenser with a Remote Liquid Storage Vessel

S. Slutsky<sup>1a</sup>, Y.-R. Yen<sup>a</sup>, H. Breuer<sup>a</sup>, A. Dobi<sup>a</sup>, C. Hall<sup>a</sup>, T. Langford<sup>a</sup>,  
D. S. Leonard<sup>a</sup>, L. J. Kaufman<sup>a</sup>, V. Strickland<sup>b</sup>, N. Voskanian<sup>a</sup>

<sup>a</sup>*Department of Physics, University of Maryland, College Park, MD 20742 USA*

<sup>b</sup>*TRIUMF, 4004 Wesbrook Mall, Vancouver, BC, V6T 2A3 Canada*

*P. Eng. at Carleton University, Ottawa, Canada*

---

## Abstract

We describe the design and operation of a system for xenon liquefaction in which the condenser is separated from the liquid storage vessel. The condenser is cooled by a pulse tube cryocooler, while the vessel is cooled only by the liquid xenon itself. This arrangement facilitates liquid particle detector research by allowing easy access to the upper and lower flanges of the vessel. We find that an external xenon gas pump is useful for increasing the rate at which cooling power is delivered to the vessel, and we present measurements of the power and efficiency of the apparatus.

*Key words:*

Xenon, Condenser, Recirculation pump

---

## 1. Introduction

Particle detectors based on condensed noble gases have found wide application in high energy physics, astro-particle physics, and nuclear physics. Noble liquids are attractive candidates for particle detectors due to their ease of purification, good charge transport properties, high scintillation efficiency, and in the case of xenon, high density and short radiation length. Examples of recent rare-event searches based on cryogenic noble gases include CLEAN/DEAP [1], XMASS [2], ZEPLIN [3], XENON [4], LUX [5], WARP [6], ArDM [7], EXO [8], ICARUS [9], and MEG [10].

The experimental methods of liquid noble gas detectors have been developed in small prototype instruments whose liquid volumes range from a few cubic centimeters to tens of liters. In these prototypes, the system is

---

<sup>1</sup>Corresponding author. E-mail: [slutsky@umd.edu](mailto:slutsky@umd.edu)

31 typically cooled either by liquid nitrogen or by directly coupling the storage  
32 vessel to the cold head of a refrigerator. These techniques have several attrac-  
33 tive features, including simplicity, robustness, and the availability of great  
34 cooling power. However, both of these cooling strategies typically require a  
35 significant amount of space above or below the storage vessel be devoted to  
36 the cryogenics. In contrast, room temperature detector technologies, such  
37 as gas proportional counters or plastic scintillator, are not burdened by a  
38 cryogenic system, which is particularly advantageous for prototyping work  
39 where the freedom to make maximal use of the space around the detector  
40 provides valuable flexibility.

41 In this article we describe a system for condensing and storing a noble gas  
42 (xenon) where the cooling system is spatially separate from the liquid storage  
43 vessel, leaving only the cryostat in the space surrounding the vessel. This  
44 configuration facilitates many common laboratory operations, particularly  
45 those which require access to the upper or lower flanges of the storage vessel.  
46 Examples include the installation of detector structures and the insertion  
47 and removal of material samples and radioactive sources. This setup also  
48 provides direct optical access to the interior of the vessel for viewing or for  
49 laser injection, and it simplifies the construction of a lead shield.

50 Our primary motivation for pursuing the remote cooling method discussed  
51 here is to allow the space above the liquid xenon vessel to be used for ion  
52 tagging and retrieval experiments in the context of the EXO double beta  
53 decay search. EXO has proposed to eliminate radioactive backgrounds by  
54 identifying the barium ion produced in the double beta decay of  $^{136}\text{Xe}$  [11].  
55 The identification method may require that a device be inserted into the  
56 active volume of the double beta decay detector to retrieve the final state  
57 nucleus. The condenser and liquid xenon vessel described in this article will  
58 allow both a barium ion calibration source and an insertion and retrieval  
59 device to be coupled to the xenon vessel from above.

60 In the last decade pulse tube cryocoolers have attracted attention as con-  
61 venient and reliable means to liquefy noble gases. For example, technology  
62 has been developed for the MEG and XENON experiments using a modified  
63 cryocooler integrated into a liquid xenon storage vessel [12]. Several other  
64 recent articles have reported on the development of a small-scale helium con-  
65 denser based on a pulse tube cryocooler where the condenser is located di-  
66 rectly in the neck of a liquid helium storage dewar [13][14][15]. In our system,  
67 we also employ a pulse tube cryocooler, and our condenser is conceptually  
68 similar to the helium condenser described in references [13]-[15], although

69 in our case the remote storage vessel represents an additional complication.  
70 Note that we report here the results of cooling and liquefaction tests carried  
71 out with xenon. Similar results could likely be obtained with other heavy  
72 noble gases having lower saturation temperature, such as argon or krypton,  
73 provided that the heat leaks in the system are minimized.

## 74 **2. Apparatus**

75 The system consists of two units: a helical copper condenser and a stain-  
76 less steel liquid xenon (LXe) storage vessel. A system drawing and plumbing  
77 schematic are shown in Figures 1 and 2, respectively. The condenser is lo-  
78 cated above the vessel, and condensed liquid flows downward through a 1/4”  
79 stainless steel (SS) tube to the top of the vessel. The condenser is cooled by  
80 the cold head of a cryocooler, while the vessel is cooled only by the xenon.  
81 The entire system is enclosed in a vacuum-insulated cryostat composed of  
82 two SS cylinders connected by a bellows. The two cylinders and bellows  
83 form one vacuum volume for pump-out purposes. The radiative heat leak is  
84 reduced by wrapping the system components in super-insulation consisting  
85 of 10-15 alternating layers of aluminized mylar and fabric.

### 86 *2.1. Condenser and Temperature Control*

87 The condenser is a helical coil of 1/4” diameter Oxygen Free High Con-  
88 ductivity copper (OFHC) tube, which was chosen for its purity and thermal  
89 conductivity. The tubing is partially annealed; non-annealed tube was found  
90 to kink when coiled. The coil is brazed to a 2” diameter cylindrical shank  
91 of OFHC copper, which is mechanically attached to the coldhead of a pulse  
92 tube cryocooler [16]. The total length of the cooled portion of the tube is 30  
93 inches (five turns).

94 As shown in Figure 3, the cold head has two stages, with base temper-  
95 atures of 22 K and 8 K for the first and second stages, respectively. These  
96 base temperatures are far below what is required to liquefy xenon (triple  
97 point of 161 K), but should allow the system to condense the lighter noble  
98 gases. To ensure temperature uniformity in the condenser, and to achieve  
99 good thermal control, each end of the condenser is coupled to one of the two  
100 stages of the cold head. We refer to the upper half of the condenser (cou-  
101 pled to the first stage of the cold head) as the “pre-cooler”, and the lower  
102 half (second stage) as the “post-cooler”. The function of the pre-cooler is to  
103 cool the room temperature xenon gas to the saturation temperature, while

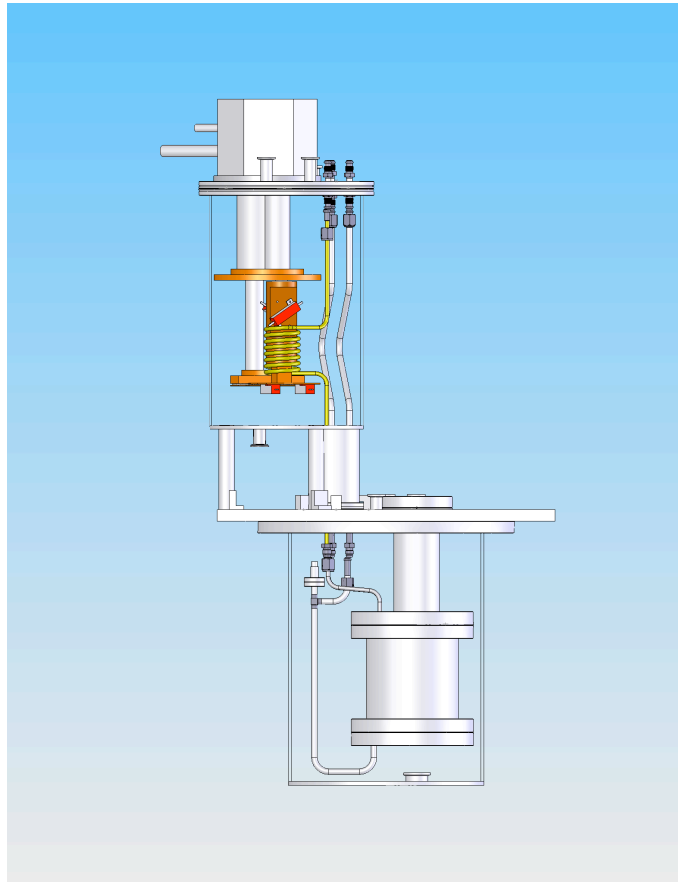


Figure 1: A drawing of the full xenon condenser system. The upper part of the cryostat holds the condenser, and the lower part holds the liquid xenon vessel.

104 the function of the post-cooler is to remove the latent heat of vaporization,  
105 thereby effecting the phase change.

106 The pre-cooler and the post-cooler are independently temperature con-  
107 trolled by trim heaters. The heaters are driven by a common 77 W regulated  
108 DC power supply, and each heater circuit is controlled by a PID temperature  
109 feedback unit [17]. The controllers adjust the current in each heater via the  
110 gate voltage on two power FETs.

111 The pre-cooler trim heaters are mounted on flats in the cylindrical shank  
112 of the condenser. These flats make a narrow “neck” between the coil and  
113 the upper cold stage, thus reducing the cooling power and allowing for more  
114 precise temperature control. A stainless steel shim between the condenser

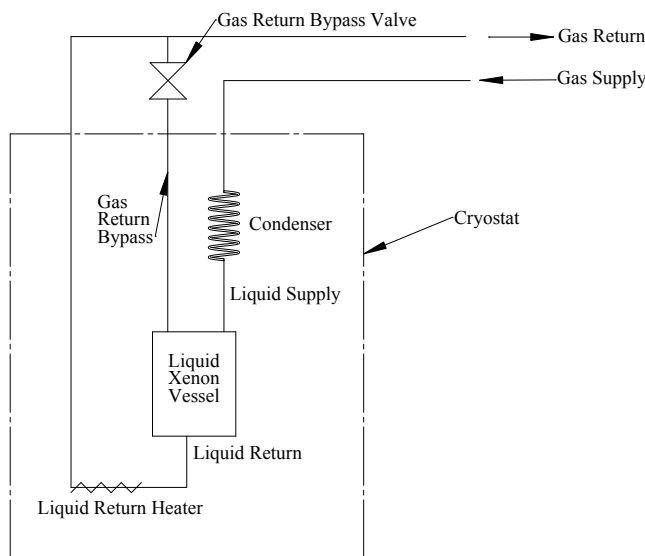


Figure 2: A schematic of the condenser system. Room temperature gas flows in the supply to the condenser. Liquid drips into the vessel, and cool gas returns through the gas return line. The gas can also return through the liquid return line when no liquid is present.

115 shank and the cold head further reduces the cooling power of the pre-cooler,  
 116 which was found to be excessive for our purposes. The post-cooler trim  
 117 heaters are attached to a plate on the bottom of the shank, which is itself  
 118 attached to the second stage of the cold head via a flexible copper braid. The  
 119 braid was chosen to provide a flexible thermal bridge so that the condenser  
 120 is mechanically constrained at its upper end only. The control temperatures  
 121 for the two PID feedback loops are measured by thermocouples located on  
 122 the neck and on the plate for the pre-cooler and post-cooler, respectively.

123 In typical operation, the pre-cooler set point temperature is chosen to be  
 124 179.5 K, and the post-cooler temperature is chosen to be a few degrees cooler.  
 125 In practice, however, the post-cooler temperature tends to stabilize around  
 126 193 K during condensation, and thus the post-cooler heater does not power  
 127 on.<sup>2</sup> This reflects the fact that the thermal coupling between the post-cooler  
 128 and the second stage of the cold head has too much thermal resistance. It

---

<sup>2</sup>Nevertheless, the post-cooler was found to be crucial for maintaining a uniform temperature throughout the condenser; see Section 6.

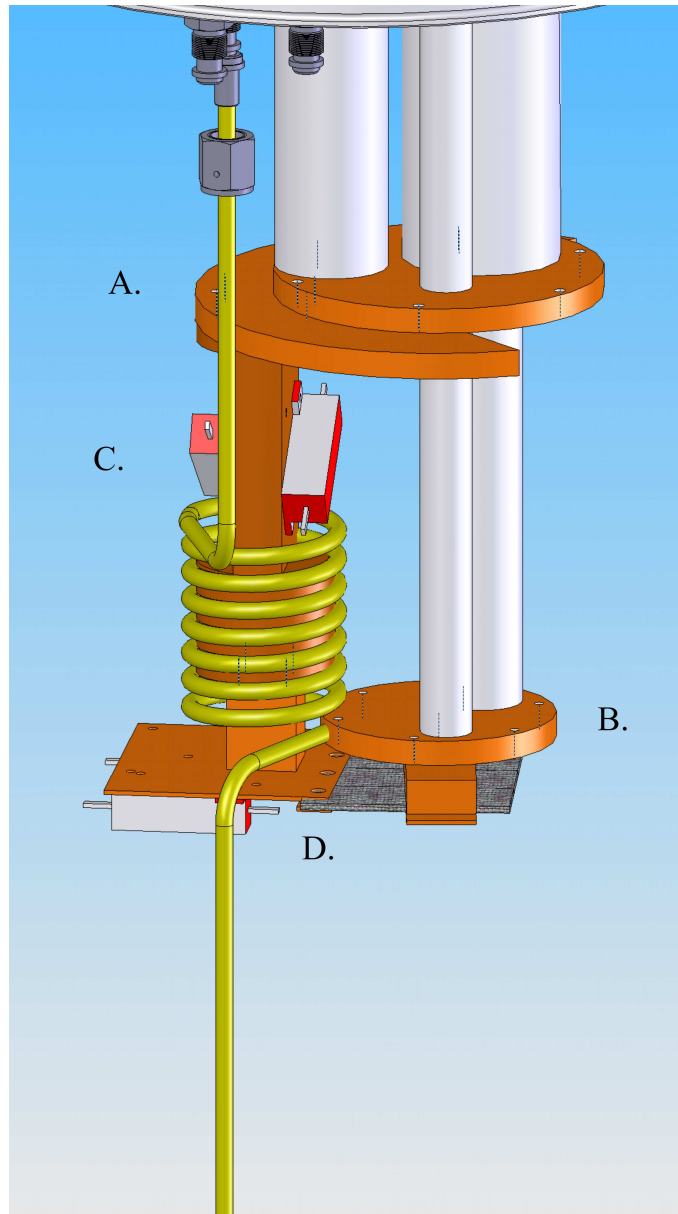


Figure 3: A drawing of the condenser. It is coupled to the cryocooler cold head at both the first stage (A) and second stage (B). The temperature at each coupling is controlled by trim heaters (C) and (D).

129 is likely that this resistance limits the power and efficiency of the condenser,  
130 so we intend to modify this arrangement in the near future. Currently, the  
131 condensation that does occur is sufficient for our purposes.

132 The trim heaters on the pre-cooler allow the condenser to adjust for the  
133 effects of changes in the gas flow rate. At high flow rates, significant heat  
134 is delivered to the condenser by incoming room temperature xenon gas, so  
135 the pre-cooler trim heater reduces its power output to maintain the set point  
136 temperature. In some cases, the incoming heat was found to be sufficient to  
137 warm the condenser to one or two degrees above the PID set point. At zero  
138 flow or low flow rates, the heat delivered to the condenser by the gas is small,  
139 and the pre-cooler trim heater provides compensation to keep the condenser  
140 temperature above the freezing point of xenon. Consequently, one benefit of  
141 the pre-cooler PID feedback loop is that it prevents over-cooling in the event  
142 that a large gas flow suddenly stops.

143 We present some measurements of the cooling power of the second stage of  
144 the PT805 near LXe temperature. (Cryomech, the cryocooler manufacturer,  
145 does not report this data at such high temperatures.) We attached a known  
146 thermal resistance to the second stage of the cold head, applied heat from  
147 the opposite end, and measured the temperature at each end. Once thermal  
148 equilibrium is reached, the cooling power is equal to the power delivered by  
149 the heater, and can also be inferred from the temperature gradient across  
150 the known thermal resistance. These two methods give results in agreement  
151 with each other. We find that the cooling power is 25 W at 65 K, 26.2 W at  
152 75 K, and 32.2 W at 152 K.

## 153 *2.2. Xenon Vessel*

154 The xenon vessel, seen in Figure 4, is a 6" x 6" OD cylinder constructed  
155 from stainless steel, which was chosen for its purity and thermal mass. This  
156 volume is sufficient to hold 10 kg, or about 3 liters, of LXe, as well as a  
157 particle detector. The top and bottom of the vessel are 8" Conflat flanges.  
158 The vessel is suspended from a vertical 3" OD stainless steel tube, 6" in  
159 length, which penetrates its top flange. At its upper end, the 3" tube is  
160 welded to a 4 5/8" Conflat flange and a large stainless steel plate which  
161 provides mechanical support. The plate also acts as part of the cryostat.  
162 The 4 5/8" flange and tube allows direct access to the interior of the xenon  
163 vessel from the laboratory. With a glass viewport attached to this flange,  
164 the LXe can be seen inside the vessel. This access port can also be used to  
165 introduce detectors, materials, or radioactive sources. A smaller, fused-silica

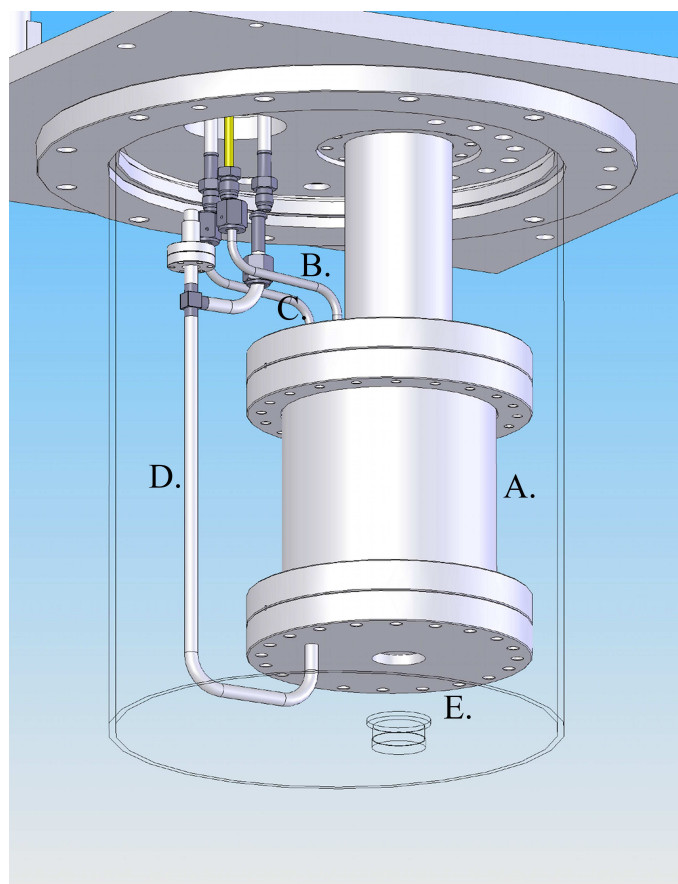


Figure 4: A drawing of the xenon vessel (A) and its associated plumbing. The liquid supply line (B) and the gas return line (C) can be seen at the top of the vessel. The liquid return line (D) connects to the bottom of the vessel and contains the capacitance level meter. The transparent outer cylinder represents the cryostat, and the laser injection windows (E) are visible at the bottom of the vessel and cryostat. The vessel is supported by the large diameter (3") tube which penetrates its top flange and permits access to the vessel interior from the lab.

166 viewport is welded into the bottom of the vessel; it is paired with an identical  
 167 viewport in the cryostat to admit laser light for use in LXe purity tests.

168 The vessel has three plumbing connections for xenon flow: a 1/4" SS  
 169 liquid supply line that enters the vessel at the top, a 3/8" liquid return line  
 170 that drains the vessel from the bottom, and a 1/4" gas return line that exits  
 171 the vessel at the top; see Figures 2 and 4. LXe from the condenser runs down



172 the liquid supply line and collects in the bottom of the vessel and in the liquid  
173 return line. A heater on the liquid return line is used to boil the liquid for  
174 re-circulation or for recovery. The gas return line at the top of the vessel  
175 allows gas to circulate freely through the vessel, especially when the vessel is  
176 filled with liquid. As discussed in Section 3 below, free recirculation of gas is  
177 important for system operation. All three plumbing lines have VCR fittings  
178 so they can be disconnected from the rest of the system, thus allowing the  
179 condenser to be removed for servicing.

### 180 *2.3. Recirculation Pump*

181 An external, custom gas recirculation pump is used to force xenon flow  
182 through the condenser and the xenon vessel. It can achieve controllable xenon  
183 flow rates of up to 10 SLPM. The pump is a bellows-type, made entirely of  
184 stainless steel, except for a teflon sleeve. The pump is driven by a 1/3-HP,  
185 three-phase motor [18]; the motor itself is controlled by an inverter [19],  
186 which allows flow control via adjustment of the repetition rate of the pump.

### 187 *2.4. Level Meter*

188 A capacitive level sensor is integrated into the liquid return line. It has  
189 a co-axial cylindrical geometry formed by suspending an 11.5" x 1/4" OD  
190 stainless steel tube inside the 3/8" OD liquid return line. The inner conductor  
191 is vented to allow the liquid to flow unimpeded. It is wrapped at the top and  
192 bottom with a small amount of Kapton tape to prevent electrical shorting.  
193 The sensor has a capacitance of  $104 \pm 3$  pF in vacuum ( $224 \pm 3$  pF including  
194 cabling ). Because the level meter is in the liquid return line, rather than the  
195 vessel itself, it will only accurately read the liquid level in the vessel when  
196 the pressures are equal. This can be ensured by keeping the gas return valve  
197 in Figure 2 open.

198 Changes in capacitance are measured with a custom circuit. The level  
199 sensor is in series with a resistor, forming a low-pass RC filter. An AC voltage  
200 of amplitude 0.15 V and frequency 8 kHz is input to the filter, and the voltage  
201 across the capacitor is amplified and rectified for computer readout by a  
202 data acquisition board [20]. To maximize sensitivity and dynamic range, the  
203 resistor in the filter was chosen such that 8 kHz is near the knee frequency.

204 The response of the circuit output was calibrated with a set of known  
205 capacitances. A quadratic fit of this data is used to interpolate for future

206 measurements. The capacitance can then be converted to a liquid level mea-  
207 surement based on the known dielectric constant of LXe. Changes in capaci-  
208 tance as small as 1 pF can be measured, corresponding to a height sensitivity  
209 of about 1 mm.

### 210 *2.5. Alarm system*

211 We have constructed an alarm system to notify lab personnel in the event  
212 of a serious system failure, such as a power outage. Three alarm conditions  
213 are considered: loss of electrical power to the cryocooler, loss of electrical  
214 power to the trim heater power supply and/or PID controllers, and an over-  
215 pressure alarm. Each alarm is represented by a switch which closes if the  
216 alarm condition is present. The switch then activates a commercial pager  
217 unit [21], which dials a list of phone numbers until the alarm is acknowl-  
218 edged. The pager unit derives power from a UPS to ensure that it remains  
219 active in the event that electrical power is lost throughout the lab. Test  
220 alarms can be generated with a push button to confirm that the system is  
221 active.

## 222 **3. The Effects of Gas Flow**

223 During normal operation, the recirculation pump is used to force gas  
224 into the condenser, where some fraction of it condenses and travels with  
225 the remaining gas flow down into the vessel. The liquid either cools the  
226 vessel through evaporation or collects in the vessel, depending on the vessel  
227 temperature. The gas usually exits the system through the gas return line,  
228 although it can return through the liquid return line as well if no liquid is  
229 present. Cryo-pumping can also be used as a source of gas flow, but this is  
230 less convenient due to the consumption of liquid nitrogen.

231 We find that forced gas flow significantly increases the cooling power  
232 delivered from the condenser to the vessel. One possible explanation for this  
233 effect is that gas flow is necessary to transport the xenon “dew” from the  
234 cold surfaces of the condenser to the remote storage vessel. Regardless of the  
235 origin of the effect, however, the importance of gas flow for cooling purposes  
236 magnifies the role of the gas return line in the system. During the initial  
237 cool down, it makes no difference if the gas return valve is open or closed,  
238 because the liquid return line provides an alternate return path. But once  
239 liquid has collected, the valve must remain open, or else the liquid will block  
240 the flow of gas, preventing cooling power from being delivered to the vessel.

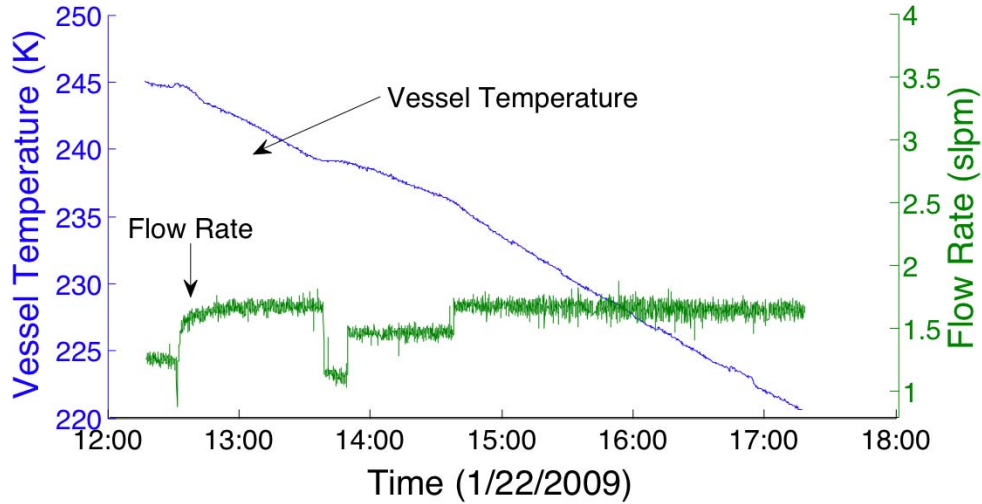


Figure 5: A comparison of the xenon vessel temperature and the gas recirculation rate during a cool-down. Note that the temperature decreases more quickly when the flow rate is increased.

241 *3.1. Gas flow during vessel cooldown*

242 Starting with the system at room temperature and under vacuum, i.e.,  
 243 with no xenon gas present, the vessel temperature drops by only a few de-  
 244 grees when the cryocooler is activated, demonstrating that conductive cool-  
 245 ing through the SS plumbing is negligible. Once xenon gas is introduced,  
 246 the xenon vessel cools through convective heat exchange with the condenser,  
 247 but the vessel temperature levels out at roughly 200 K in the absence of  
 248 forced gas flow. We find that a minimum gas flow rate of about 1.2 SLPM  
 249 is typically necessary to cool the xenon vessel from room temperature to the  
 250 saturation temperature. At higher flow rates, from 2 to 4 SLPM, the vessel  
 251 cools more quickly, as shown in Figure 5. The effects of increasing the flow  
 252 rate even further are unclear. In one instance, increasing the flow rate to

253 5 SLPM decreased the vessel cooling rate, since the large amount of incom-  
254 ing warm gas heated the condenser above LXe temperature. However, in  
255 another instance, flows of 6.8 SLPM led to greater cooling rates in the vessel,  
256 so there may be other effects, such as pressure, that play a role. Cooling  
257 rates are discussed further in Section 5.2.

258 Forced gas flow may only be critical during the initial phase of vessel  
259 cooling. During one cool-down, gas flow and condensation were established,  
260 but the recirculation pump was unexpectedly halted, interrupting flow for  
261  $\sim 1$  hour. Condensation continued during this time. Upon restarting the  
262 pump, the xenon gas pressure rose and would not stabilize until the pump  
263 was shut off again. However, condensation continued and the liquid level  
264 in the vessel rose without the forced flow. It is therefore not clear if the  
265 1.2 SLPM flow rate is an absolute limit for vessel cooling; it may only be  
266 required to initialize condensation.

### 267 *3.2. Gas flow during liquid maintenance*

268 Once the vessel is cold and filled with liquid, it is easy to maintain at  
269 constant temperature and pressure for indefinite periods of time by forcing  
270 gas flow through the condenser and liquid vessel, using the gas return line  
271 as the exhaust. In this arrangement the gas circulates in a closed loop.  
272 This is generally our default configuration during liquid maintenance, and it  
273 allows for continuous purification of the xenon gas with a gas phase purifier.  
274 However, we have also studied the possibility of maintaining the liquid at  
275 constant temperature and pressure with the recirculation pump turned off.  
276 This situation could be important, for example, if the pump were to lose  
277 power unexpectedly. As detailed below, we find that it is possible, but more  
278 difficult, to achieve stability with the recirculation pump turned off.

279 Note that the xenon will thermodynamically recirculate at a small rate (a  
280 few SCCM) if the external gas plumbing allows gas to flow from the system  
281 output to the condenser input. This “thermal recirculation” is driven by the  
282 system heat leak, and it can be augmented by warming the liquid return line  
283 with its heater.

284 We achieve temperature and pressure stability most easily in the absence  
285 of forced gas flow by closing a valve in the external plumbing which prevents  
286 thermal recirculation. Under these conditions, we expect that no gas will  
287 enter the condenser input from the external system. Our flowmeter confirms  
288 this expectation. Nevertheless, we see through the large viewport that a  
289 steady stream of liquid drops falls into the vessel from above. This indicates

290 that gas is counter-flowing up the liquid supply line from the vessel, lique-  
291 fying, and falling back down, creating a closed heat exchange loop. This  
292 behavior is rather sensitive to the condenser temperature. For example, rais-  
293 ing the pre-cooler temperature from 179.5 K to 180.5 K is enough to disturb  
294 the establishment of this heat exchange loop.

295 If we turn off the external recirculation pump without preventing thermal  
296 recirculation, then we usually find that condensing slows or stops, and that  
297 the system temperature and pressure slowly rise. Since the previous tests  
298 show that the system is able to condense the cold counter-flowing gas from  
299 the LXe vessel, this behavior could indicate that the extra heat load from the  
300 room temperature gas is too large. It is possible that this situation could be  
301 remedied by improving the cooling power of the condenser, particularly that  
302 of the post-cooler. Also note that if we encourage thermal recirculation with  
303 the heater, this makes the situation worse, increasing the rate of temperature  
304 and pressure rise.

305 On the other hand, if we establish the internal heat exchange loop by  
306 preventing external thermal recirculation, and we allow the system to run  
307 in this mode for several hours, we find that opening the external valve to  
308 allow thermal recirculation does not disturb the system stability. That is,  
309 the system behaves as if the valve were still closed.

## 310 4. Operation

311 The operation of the system is divided into four stages: preparation,  
312 vessel cool-down and filling, liquid maintenance, and recovery.

### 313 4.1. Preparation

314 To remove impurities in the system, the condenser and vessel are evac-  
315 uated to  $\sim 10^{-7}$  torr using a turbomolecular pump backed by a dry scroll  
316 pump. The cryostat is evacuated to  $\sim 10^{-3}$  torr using a similar configura-  
317 tion, and it is pumped continuously while the cryocooler is running. A xenon  
318 gas supply bottle with a regulator is opened, and the regulator is set to intro-  
319 duce gas at a pressure between 1000 and 1400 torr. The recirculation pump  
320 is turned on and set to a flow rate of at least 1.2 SLPM. The gas circulates  
321 though the condenser, vessel, and external plumbing, and optionally through  
322 a zirconium getter for purification [22]. The pressure regulator on the xenon  
323 gas supply bottle is left open to allow additional gas to enter the system as

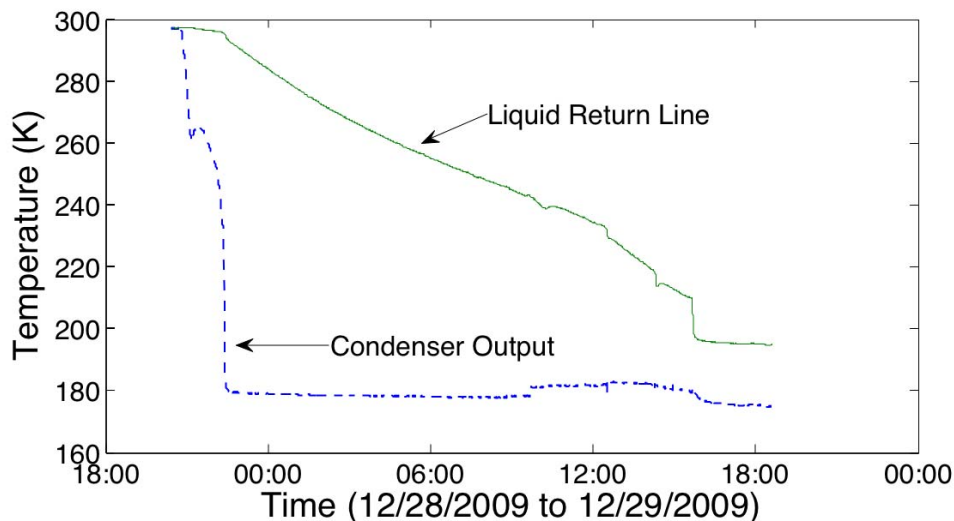


Figure 6: Temperature of the output of the condenser and temperature of the liquid return line of the vessel during a typical cool-down. The condenser output temperature line is kinked at about 260 K when the heaters turn on. Other kinks are due to intentional changes in flow rate.

324 needed. The pre-cooler set point temperature is set to 179.5 K. Once this is  
 325 done, the system is prepared for cool-down.

#### 326 4.2. Vessel cool-down and filling

327 During the vessel cool-down and filling phase, the gas continues to cir-  
 328 culate through the condenser, vessel, and external plumbing. As the gas  
 329 cools and becomes more dense, and later as the gas liquefies, additional gas  
 330 is delivered from the supply bottle to the system as needed to maintain a  
 331 constant pressure. The total amount of xenon in the liquid system can be  
 332 monitored by measuring the remaining pressure in the supply bottle.

333 Figure 6 shows a temperature history of a typical cooldown as recorded

334 by thermocouples at the output of the condenser and at the liquid return  
335 line near the bottom of the vessel. The cryocooler is activated to begin cool-  
336 ing the system, and the condenser temperature immediately drops. When  
337 the temperature reaches 260 K, the trim heaters turn on for the first time,  
338 momentarily warming the condenser. (The heaters are programmed to stay  
339 off when the condenser is above this temperature. This avoids overheating  
340 the cold head). After  $\sim 1$  hour, the temperature of the output of the con-  
341 denser falls sharply, indicating LXe has formed and is pouring down into the  
342 vessel. After this sharp temperature decrease, drops of LXe can be seen via  
343 the viewport falling onto the bottom of the vessel and boiling away, and the  
344 temperature of the LXe return line decreases at about six times the previous  
345 rate. Note that this indicates that most of the cooling power is delivered to  
346 the vessel in the form of liquid xenon, rather than gaseous xenon.

347 The vessel cools this way for some time, typically about 8 hours, but it  
348 can be more or less depending on flow and pressure conditions. When the  
349 vessel temperature is low enough, liquid begins to collect in the bottom of the  
350 vessel. When enough liquid has collected, it overflows a small lip at the drain  
351 of the vessel that leads to the liquid return line. This “splash” brings the  
352 cold liquid into direct contact with the return line, causing the temperature  
353 there to quickly drop, as can be seen in Figure 6 around 16:00 on 12/29. This  
354 drop is correlated with a quick rise in the measured liquid level as the level  
355 meter is filled for the first time.

356 After the “splash” of liquid fills the return line, LXe continues to collect in  
357 the vessel. At this stage, the gas return bypass valve must be open to ensure  
358 that a high rate of gas flow can continue. (We typically keep this valve open  
359 through the entire process, from initial preparation through xenon recovery.)  
360 Condensing 1 kg of xenon takes 2-3 hours, depending on pressure and flow  
361 conditions.

### 362 *4.3. Liquid Maintenance*

363 Liquid can be maintained indefinitely in the vessel at constant tempera-  
364 ture and pressure by keeping the pre-cooler set point temperature at 179.5 K,  
365 and by maintaining a nominal gas flow rate of 1-2 SLPM with the recircula-  
366 tion pump. As described in Section 3.2, the recirculation pump can also be  
367 turned off under certain conditions.

#### 368 4.4. Xenon Recovery

369 The xenon is recovered by cryopumping. The room temperature gas  
370 supply bottle, which is made of aluminum, is placed in a liquid nitrogen  
371 bath, freezing any xenon inside. The supply bottle pressure regulator is fully  
372 opened to allow xenon gas to flow backwards through it. In addition, a  
373 bypass valve which is connected in parallel with the pressure regulator can  
374 be opened to reduce the pumping impedance further, but this is usually not  
375 necessary.

376 Finally the cryocooler is turned off and the LXe is allowed to warm,  
377 raising the vapor pressure. A hand valve between the liquid system and the  
378 cold supply bottle is used to manually regulate the gas flow rate. Care should  
379 be taken to avoid forming xenon ice by pumping too quickly. Recovery can  
380 be made quicker by heating the liquid return line directly with its integrated  
381 heater, or by filling the cryostat insulating vacuum with a dry gas such as  
382 nitrogen or helium to provide a thermal connection to the room temperature  
383 lab.

### 384 5. Measurements

#### 385 5.1. Liquid Level

386 Figure 7 shows a sample plot of the height of LXe in the vessel during a  
387 cool-down. The gas pressure was 1000 torr. The recirculation pump forced  
388 a constant flow rate of 2.8 SLPM from 11:00 to 18:45 PM, after which it  
389 was turned off, and the xenon supply bottle was valved off. The initial rapid  
390 rise from -2.7 cm to 0.6 cm corresponds to the “splash” when liquid first  
391 overflows from the vessel into the return line. After that, a roughly linear  
392 increase in the height of the xenon can be seen as the storage vessel fills with  
393 liquid. A surprising feature is that from 6:45 PM until 9:00 AM the next  
394 day, the liquid level is seen to be rising slowly. This could indicate that the  
395 level sensor circuit has a slowly drifting systematic error.

#### 396 5.2. Cooling Rate of the Vessel

397 The cooling rate of the stainless steel vessel was measured on two separate  
398 trials. During the first run, the flow rate was maintained at about 1.4 SLPM,  
399 with a xenon pressure of 1500 torr. The vessel cooled at 5.6 K/hr. During  
400 another run, the vessel cooled at 6.5 K/hr, with a much higher gas flow rate  
401 of 6.8 SLPM and pressure of 1500 torr. Using a specific heat for stainless steel



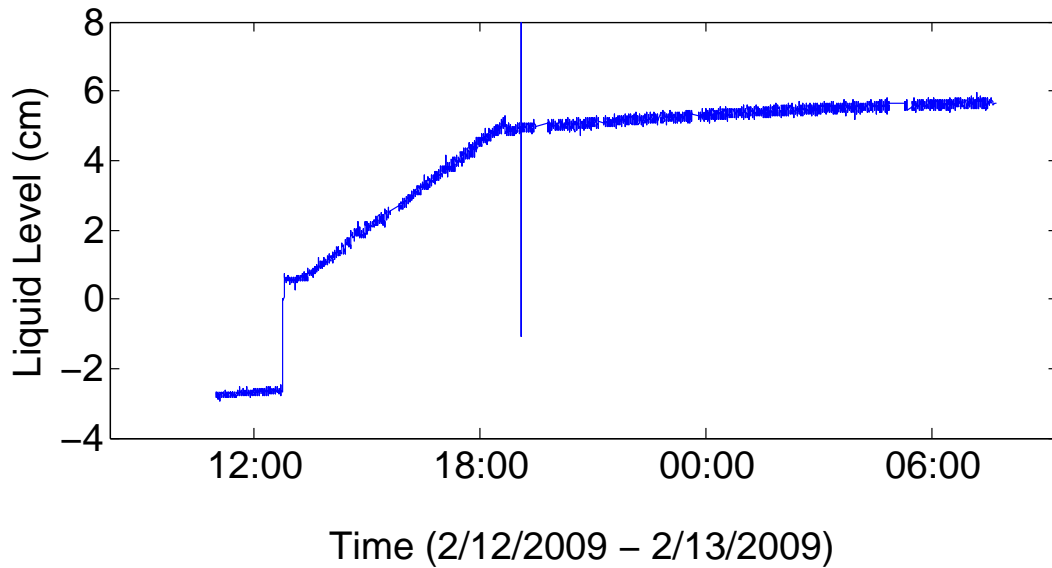


Figure 7: A time history of the LXe level during condensation. The fast rise from negative to positive values around 14:00 is due to the filling of the liquid return line, while the slower rise over the next several hours is the filling of the storage vessel. Note that the spike around 19:00 is spurious.

402 of 500 J/kg/K, and estimating the vessel weight at 20 kg, we find that these  
 403 cooling rates correspond to 15.5 W and 18.9 W of cooling power transferred  
 404 to the vessel.

### 405 5.3. Condensation Rate

406 Two trials are presented in which the condensation rate is measured at  
 407 various pressure and flow parameters. Condensation rate is measured during  
 408 vessel filling in two ways. 1.) The mass of xenon remaining in the gas  
 409 supply bottle is calculated from the known volume of the bottle, the measured  
 410 pressure, and the density of xenon gas at that pressure. (The density of xenon  
 411 gas has a non-linear relationship to the pressure at typical bottle pressures,  
 412 and this dependence must be taken into account.) The rate of decrease of  
 413 source mass was then calculated using central differencing and equated to  
 414 the rate of condensation in the liquid system. 2.) The height of xenon in  
 415 the vessel is measured, and converted into a volume and mass using the  
 416 known cylindrical geometry of the system and the density of LXe. The rate  
 417 of change in liquid height then equates to the condensation rate. These two

418 methods give results in good agreement with each other, although the liquid  
419 height method is complicated by the volume displacement of the irregularly  
420 shaped particle detector in the LXe vessel. In the following, we quote results  
421 based on the bottle pressure method.

422 In the first trial, the xenon gas pressure was set to 1550 torr, and the  
423 gas flow rate was varied. The condensation rates were 0.54 kg/hr for a flow  
424 of 2.65 SLPM and 0.61 kg/hr for a flow 3.89 SLPM. In a second trial, the  
425 gas pressure was set to 1000 torr. The recirculation flow rate was set to  
426 2.8 SLPM for the bulk of the trial. Condensation rates between 0.36 and  
427 0.40 kg/hr were observed.

428 The largest rate of condensation, 0.61 kg/hr, was observed at a high  
429 pressure and large recirculation rate, 3.89 SLPM and 1550 torr. This implies  
430 that 44% of the circulating gas is condensed in a single pass. A condensation  
431 fraction of 58% was also achieved, but only at the cost of a lower condensation  
432 rate of 0.54 kg/hr. Tests at a lower pressure of 1000 torr indicate both lower  
433 condensation rate,  $\sim 0.36$  kg/hr, and a lower efficiency of 36%. Thus, the  
434 condensation fraction depends moderately on pressure and flow.

435 Using the latent heat of vaporization of xenon, 12.64 kJ/mol, and the spe-  
436 cific heat of xenon gas, 20.8 J/mol/K, we can calculate the cooling power im-  
437 plied by our condensing rates. For the largest condensation rate, 0.61 kg/hr,  
438 we find a cooling requirement of 3.5 W and a latent heat removal of 16.4 W,  
439 for a total of 19.9 W. This is similar to the 18 W of cooling power we esti-  
440 mated is delivered to the stainless steel vessel during cooling.

## 441 6. Discussion

442 We find that for reliable operation of our condenser it is essential to  
443 control the temperature at both the top and bottom. Initial tests in which  
444 the condenser was cooled only at the top showed that a large temperature  
445 gradient would appear along its length. In this arrangement, the top of the  
446 condenser must be over-cooled to allow liquefaction to occur in the lower  
447 portion of the condenser. This leaves the condenser prone to ice formation,  
448 particularly if the gas flow rate suddenly decreases, removing a heat source.

449 With a dual control mechanism, both ends of the condenser are cooled,  
450 ensuring that temperature gradients are small. In addition, the two thermo-  
451 dynamic functions of the condenser (cooling warm gas and liquefaction) are  
452 separated spatially in the condenser, allowing for semi-independent tempera-  
453 ture regulation with a pre-cooler and a post-cooler. This gives the condenser

454 the flexibility to adjust to changing conditions, such as a change in the gas  
455 flow rate. In practice, our post-cooler temperature is usually above its set  
456 point temperature during condensation, and therefore its trim heater plays  
457 little role. Nevertheless, the additional cooling provided by the second stage  
458 of the cold head through the post-cooler improves temperature uniformity in  
459 the condenser, leading to more robust operation. To achieve greater control  
460 in the future, we intend to increase the cooling power of the post-cooler by  
461 reducing its thermal resistance.

462 Xenon ice formation is a dangerous problem for external condensers such  
463 as the one described here. Ice can block the flow of liquid and gas to the  
464 xenon vessel. Since the xenon flow is the only cooling mechanism for the  
465 vessel, its interruption can lead to a dangerous rise in system temperature  
466 and pressure. Our design is particularly sensitive to this problem, because  
467 the helical coil of the condenser has a cross-sectional outer diameter of only  
468 1/4". Therefore even a small amount of ice can lead to flow blockage. A  
469 condenser coil made from larger diameter tubing may improve the situation,  
470 or perhaps a condenser with an altogether different geometry may be better.  
471 For our purposes, however, the current design has proved to be adequate  
472 with appropriate safeguards.

473 Gas flow is essential for transferring the cooling power from the condenser  
474 to the xenon vessel, and a xenon gas return line from the top of the vessel  
475 greatly improves the effectiveness. There is a minimum flow rate necessary  
476 to cool the vessel to LXe temperature, and the condensation rate increases,  
477 albeit only slightly, with increasing flow. Cryopumping can serve as the  
478 method for forcing flow, but this is an awkward process that consumes large  
479 amounts of liquid nitrogen, and flow must be interrupted frequently to warm  
480 and cool the supply and recovery bottles. Therefore we find it is very useful to  
481 have a recirculation pump. Using the pump, gas flow can be easily maintained  
482 for days at a time, and with an inverter controller, the flow rate can be dialed  
483 to a desired value for easy testing. The main drawback to the pump is a loss  
484 of purity: our custom pump contains a teflon sleeve which is a source of  
485 outgassing and teflon debris in the system. These problems can be solved  
486 with purifiers and filters, however.

## 487 **7. Conclusion**

488 We have described a system for condensing and storing xenon where the  
489 source of cooling power has been removed from the vicinity of the liquid

490 storage vessel, facilitating the introduction of instruments and materials to  
491 the vessel.

492 Condensation rates as high as 0.61 kg/hr were achieved, after an initial  
493 cool-down period of 8-10 hours. This corresponds to a condensation fraction  
494 of 44% and a cooling power of about 20 W. Changes in condenser design may  
495 be able to improve the condensation fraction and cooling power. Our design  
496 includes a two-stage cooling system for improved temperature uniformity  
497 and control. We find that a nominal gas flow rate is important for delivering  
498 cooling power to the vessel, and that a dedicated gas return line is useful for  
499 maintaining this flow when the vessel is filled with liquid.

## 500 8. Acknowledgments

501 We thank John Carriker for his many contributions to the system de-  
502 scribed in this article. This work was supported by award number 0652690  
503 from the National Science Foundation.

## 504 References

- 505 [1] W. H. Lippincott, *et. al.*, “Scintillation time dependence and pulse  
506 shape discrimination in liquid argon”, *Phys. Rev. C* **78**, 035801 (2008),  
507 arXiv:0801.1531[nucl-ex].
- 508 [2] K. Abe *et. al.* “The XMASS experiment”, *J. Phys.: Conf. Ser.* **120**  
509 042022 (2008)
- 510 [3] G.J. Alner, *et. al.* (ZEPLIN collaboration), “Limits on spin-dependent  
511 WIMP-nucleon cross-sections from the first ZEPLIN-II data”, *Phys.*  
512 *Lett. B* **653**:161-166 (2007), arXiv:0708.1883v1.
- 513 [4] J. Angle, *et. al.* (XENON collaboration), “First Results from the  
514 XENON10 Dark Matter Experiment at the Gran Sasso National Labo-  
515 ratory”, *Phys. Rev. Lett.* **100**:021303 (2008), arXiv:0706.0039v2.
- 516 [5] [www.luxdarkmatter.org](http://www.luxdarkmatter.org).
- 517 [6] P. Benetti, *et. al.* (WARP collaboration), “First results from a Dark  
518 Matter search with liquid Argon at 87 K in the Gran Sasso Under-  
519 ground Laboratory”, *Astropart. Phys.* **28**:495-507 (2008), arXiv:astro-  
520 ph/0701286v2.

- 521 [7] M. Laffranchi and A. Rubbia (ArDM collaboration), “The ArDM  
522 project: a liquid argon TPC for dark matter detection”, arXiv:hep-  
523 ph/0702080v1.
- 524 [8] C. Hall, “Searching for Double Beta Decay with the Enriched Xenon  
525 Observatory”, Proceedings of the 9th Conference on the Intersections of  
526 Particle and Nuclear Physics (CIPANP 2006), AIP Conf. Proc. 870:532-  
527 535, 2006.
- 528 [9] A. Badertscher (ICARUS collaboration), “The ICARUS project: a 3000  
529 t LAr TPC for neutrino physics and a search for nucleon decays” Nucl.  
530 Instr. Meth. A **535**, 129 (2004).
- 531 [10] S. Ritt (MEG collaboration), “Status of the MEG experiment  $\mu \rightarrow e\gamma$ ”  
532 Nucl. Phys. B Proc. Suppl. **162**, 279 (2006).
- 533 [11] M. Danilov, *et. al.* (EXO collaboration), “Detection of very small neu-  
534 trino masses in double beta decay using laser tagging”, Phys. Lett. B  
535 **480**, 12 (2000).
- 536 [12] T. Haruyama, “Progress of xenon liquefaction technology by using a  
537 pulse tube cryocooler”, 2005 IEEE International Conference on Dielec-  
538 tric Liquids **26**, 353 (2005).
- 539 [13] C. Wang, “Helium liquefaction with a 4 K pulse tube cryocooler”, Cryo-  
540 genics, **41**, 491 (2001).
- 541 [14] C. Wang, “Efficient helium recondensing using a 4K pulse tube cry-  
542 ocooler”, Cryogenics **45**, 719 (2006).
- 543 [15] C. Wang and R.G. Scurlock, “Improvement in performance of cryocool-  
544 ers as condensers”, Cryogenics **48**, 169 (2008).
- 545 [16] Model PT805 from Cryomech, Inc., 113 Falso Drive, Syracuse, New  
546 York, 13211, www.cryomech.com.
- 547 [17] Model Cni16D from Omega Engineering, Inc., One Omega Drive, Stam-  
548 ford, Connecticut, 06907-0047, www.omega.com.
- 549 [18] Dayton model 3N694BA.

- 550 [19] Frenic-Mini Type FRNF25C1S-6U from Fuji Electric,  
551 [www.fujielectric.com](http://www.fujielectric.com).
- 552 [20] Model PCI-6259 from National Instruments Corp., 11500 N Mopac Ex-  
553 pwy, Austin, TX 78759-3504, [www.ni.com](http://www.ni.com).
- 554 [21] Model 1512 Sentrifugal from Zetron, Inc., PO Box 97004, Redmond, WA  
555 98073-9704, [www.zetron.com](http://www.zetron.com).
- 556 [22] SAES Monotorr purifier model PS4-MT3, [www.saespuregas.com](http://www.saespuregas.com).