

A LEGO Watt Balance: An apparatus to demonstrate the definition of mass based on the new SI

L.S. Chao, S. Schlamminger, D.B. Newell, and J.R. Pratt

Physical Measurement Laboratory, National Institute of Standards and Technology, Gaithersburg, MD 20899

G. Sineriz, F. Seifert, A. Cao, and D. Haddad

Joint Quantum Institute, University of Maryland, College Park, MD 20742

X. Zhang

Computational Instrumentation Lab, Massachusetts Institute of Technology, Cambridge, MA 02139

(Dated: December 7, 2024)

A global effort to redefine our International System of Units (SI) is underway and the change to the new system is expected to occur in 2018. Within the newly redefined SI, the present base units will still exist but be derived from fixed numerical values of seven reference constants. More specifically, the unit of mass, the kilogram, will be realized through a fixed value of the Planck constant h . For instance, a watt balance can be used to realize the kilogram unit of mass within a few parts in 10^8 . Such a balance has been designed and constructed at the National Institute of Standards and Technology. For educational outreach and to demonstrate the principle, we have constructed a LEGO tabletop watt balance capable of measuring a gram size mass to 1% relative uncertainty. This article presents the design, construction, and performance of the LEGO watt balance and its ability to determine h .

I. INTRODUCTION

The quest for a redefined International System of Units (SI) has been a formidable global undertaking. If the effort concludes as expected, sometime in 2018 our seven base units (meter, kilogram, second, ampere, kelvin, mole, candela) that have formed the foundation of our unit system for over half a century will be redefined via seven reference constants. In terms of mass metrology, the present standard, forged in 1879 and named the International Prototype Kilogram (IPK), is the only mass on Earth defined with zero uncertainty. However, the base unit kilogram will be redefined via a fixed value of the Planck constant h , finally severing its ties to the IPK. Different experimental approaches can be used to realize mass from the fixed value of h . At the National Institute of Standards and Technology (NIST), we have chosen to pursue the watt balance to realize the kilogram in the US after the redefinition.¹

The watt balance, first conceived by Dr. Bryan Kibble in 1975, is a weight-measuring apparatus that forms a virtual comparison between electrical power and mechanical power.² Since power is measured in units of watts, the device has been dubbed a “watt balance”. It is essentially a force transducer that can be calibrated solely in terms of electrical, optical, and frequency measurements. A few watt balances around the world have demonstrated the capability of measuring 1 kg masses with a relative uncertainty of a few parts in 10^8 . Here, under the inspiration of Terry Quinn,³ we describe the construction of a tabletop LEGO⁴ watt balance capable of measuring gram-level masses with a much more modest relative standard uncertainty of 1%. For the instrument described here, the cost of parts totaled about \$650, but, such a device sim-

ilar to it can be built for significantly less. The largest fraction of the cost is in the data acquisition system used to transfer the data to a laptop. A recommended parts-list is provided in Appendix A. We encourage readers to use this manuscript as general guidance for constructing such a device and by no means as a definitive prescription. There are many ways to build a watt balance, and we consider here a concept to highlight general details that are most important for success.

II. BASIC WATT BALANCE THEORY

Although a watt balance might appear functionally similar to an equal arm balance, an equal arm balance is passive, relying on comparing an unknown mass to a calibrated one, while a watt balance is active, relying on compensating the unknown mass with a known force. In this case, the weight of an object is compensated by a precisely adjusted electromagnetic force.

The experiment consists of two modes of operation illustrated in Fig. 1: velocity mode and force mode. Velocity mode is based on Faraday’s Law of Induction. A coil (wire length L) is moved at a vertical speed v through a magnetic field (flux density B) so that a voltage V is induced. The induced voltage is related to the velocity through the flux integral BL :

$$V = BLv. \quad (1)$$

Force mode is based on Lorentz Forces. The gravitational force on a mass m is counteracted by an upward electromagnetic force F generated by the now-current-carrying

coil in a magnetic field:

$$F = BLI = mg \quad (2)$$

where g is the local gravitational acceleration and I is the current in the coil.

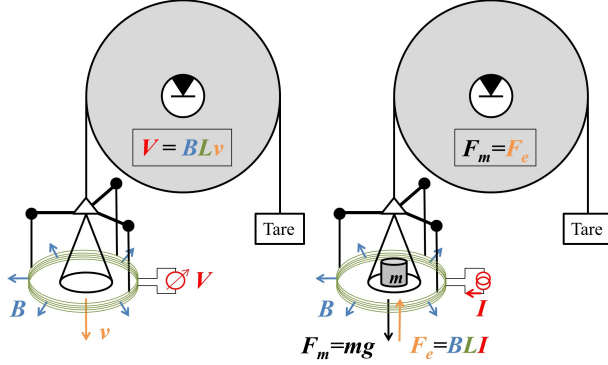


FIG. 1. Left: velocity mode. The coil moves vertically in a radial magnetic field and a voltage V is induced. Right: force mode. The upward electromagnetic force generated by the coil opposes the gravitational force exerted by m .

Canceling out the BL factor common to equations (1) and (2) and rearranging the variables, expressions for electrical and mechanical power are equated and a solution for mass is obtained:

$$VI = mgv \Rightarrow m = \frac{VI}{gv}. \quad (3)$$

The equation above relates mechanical power to electrical power and provides a means to relate mass to electrical quantities. Now, in order to make the connection from mass to the Planck constant through the electrical quantities, it is necessary to understand two quantum physical effects that have revolutionized electrical metrology since the second half of the last century: the Josephson effect and the quantum Hall effect. These two phenomena are what permit the measurement of electrical quantities in terms of the Planck constant to the precision required for the watt balance and redefinition.

The Josephson effect can be observed in a Josephson junction which consists of two superconducting materials separated by thin non-superconducting barrier. At superconducting temperatures, and while irradiating the junction with an electrical field at a microwave frequency f , a bias current is forced through this junction and a voltage of

$$V = \frac{h}{2e} f \equiv K_J^{-1} f \quad (4)$$

will develop across the junction. The quotient $K_J = 2e/h$ is named the Josephson constant in honor of Brian

Josephson who predicted this effect in 1960.⁵ One junction delivers only a small voltage, typically $37 \mu\text{V}$, so, in order to build a practical voltage standard, tens of thousands of these junctions are connected in series on a single chip. At NIST,⁶ a chip with approximately 250,000 junction is used to produce any voltage up to 10 V with an uncertainty of 1 nV.

The quantum Hall effect is a special case of the Hall effect. The Hall effect occurs when a current carrying conductor is immersed in a magnetic field and a Hall Voltage V_H occurs perpendicular to the magnetic flux and the current. While in the classical Hall effect the conductor immersed is a three dimensional object, in the quantum Hall effect, the electrical conduction is confined to two dimensions. In such a system and at sufficiently high magnetic field, the ratio between the Hall Voltage and current, or Hall Resistance R_H , becomes quantized to

$$R_H = \frac{V_H}{I} = \frac{1}{i} \frac{h}{e^2} \equiv \frac{1}{i} R_K \quad (5)$$

where i is an integer. The quotient $R_K = h/e^2$ is named the von Klitzing constant to honor Klaus von Klitzing who discovered this effect first in 1980 (discussion of this effect can be found in reference).⁷ At NIST, the quantum Hall effect is the starting point of resistance dissemination.⁸ Scaling with a cryogenic current comparator allows researchers to measure a 100Ω precision resistor with a relative uncertainty of a few parts in 10^9 .

The watt balance experiment employs these two quantum effects to measure the electrical power. However, instead of directly measuring $P = VI$, the current I is driven through a precisely calibrated resistor R , producing a voltage drop V_R , yielding $P = VV_R/R$. Both voltages are measured by comparing to a Josephson voltage standard, so their values can be expressed in terms of a frequency and the Josephson constant. The resistor is measured by comparing to a quantum Hall resistor, so its value can be expressed in terms of R_K . This can be written as

$$P = VV_R/R = C f_1 f_2 \frac{h}{2e} \frac{h}{2e} \frac{e^2}{h} = \frac{C f_1 f_2}{4} h. \quad (6)$$

Here, C is a known constant that indicates the number of junctions used and the ratio of R to h/e^2 . Combining the above equation with the watt Eq. (3) yields

$$h = \frac{4}{C f_1 f_2} mgv \Rightarrow m = \frac{C f_1 f_2}{4} \frac{h}{gv}. \quad (7)$$

Before the 2018 redefinition of units, the equation on the left is used to measure h from a fixed value of mass. After redefinition, the equation on the right will be used to realize mass from a fixed value of h .

In a classroom setting, quantum electrical standards are typically unavailable. However, it is still possible to

measure the Planck constant due to the way the current unit system is structured. While the SI is used for most measurements, a different system of units has been used for electrical measurements since 1990. For these so-called conventional units, the Josephson and von Klitzing constants were fixed at values adjusted to the best knowledge in 1989.^{9,10} These fixed values are named “conventional Josephson” and “conventional von Klitzing” constants and are abbreviated K_{J-90} and R_{K-90} , respectively. Since 1990, almost all electrical measurements are calibrated in conventional units. By comparing electrical power in conventional units to mechanical power in SI units, h can be determined.

Starting at equation 3, we see that

$$VI = mgv \implies \{VI\}_{90} W_{90} = \{mgv\}_{SI} W_{SI}, \quad (8)$$

where $\{x\}_{90/SI}$ denotes the numerical value of the quantity x in conventional or SI units, respectively. Further, W_{90} and W_{SI} are the units of power (watt) in the conventional and SI system. The above equation can be written as

$$\frac{\{mgv\}_{SI}}{\{VI\}_{90}} = \frac{W_{90}}{W_{SI}} = \frac{h}{h_{90}} \implies h = h_{90} \frac{\{mgv\}_{SI}}{\{VI\}_{90}}, \quad (9)$$

where h_{90} is the conventional Planck constant defined as,

$$h_{90} \equiv \frac{4}{K_{J-90}^2 R_{K-90}} = 6.626\,068\,854 \dots \times 10^{-34} \text{ J s}. \quad (10)$$

Thus, the value of the Planck constant can be determined by multiplying the conventional Planck constant by the ratio of mechanical power in SI units to the electrical power in conventional units.

In practice, the easiest way to compare the quotient of the SI power to the conventional power is by assigning different flux integrals BL to each mode, i.e.,

$$(BL)_V = \frac{V}{v} \text{ and } (BL)_F = \frac{mg}{I} \quad (11)$$

Using these two numbers, the ratio of h/h_{90} is given by

$$\frac{h}{h_{90}} = \frac{(BL)_F}{(BL)_V} = \frac{\{mgv\}_{SI}}{\{VI\}_{90}}. \quad (12)$$

After redefinition, electrical power and mechanical power will be measured in the same units and the schism between units will vanish. Then, referring back to Eq. (3), mass can be realized using a watt balance simply as:

$$m = \frac{VI}{gv} \quad (13)$$

where all quantities are expressed in SI units.

The remaining two variables g and v are measured accurately with an absolute gravimeter and interferometric methods, respectively. However, since this manuscript's

main focus is still a proof-of-principle LEGO watt balance, ultra high precision metrology approaches are unnecessary. Gravity can be estimated by inputting one's geographical coordinates into the web page found in reference.¹¹ Velocity can be determined using a simple optical method we describe in Section V.

Like all great LEGO systems, the LEGO watt balance is versatile, capable of measuring in either mode. It will be a device to measure Planck's constant before redefinition and one to realize mass after redefinition. A capable operator can perform a measurement with a relative uncertainty of 1% with the device described below.

III. LEGO WATT BALANCE MECHANICS

We chose a symmetric design for the LEGO watt balance that conforms to easily recognized notions of an equal-arm beam balance. We reiterate that there are many ways to construct a watt balance. One way is described below. Fig. 2 shows a CAD drawing of our balance. Two weighing pans are suspended from either arm of the balance which pivots about its center. Suspended below each weighing pan is a wire-wound coil immersed in a radial magnetic field.

The magnet system we chose to generate this radial magnetic field consists of a pair of neodymium (N48) ring magnets, one pair per coil. For simplicity, we recommend keeping the system an open-field design, i.e., “yokeless”, meaning no additional ferromagnetic material to guide the magnetic flux direction. The dimensions of the ring magnets were chosen such that they could fit inside the PVC pipe coil former with approximately 0.5 cm clearance all around. A brass threaded rod secured to a non-magnetic base plate (wood, aluminum, etc.) provides the vertical guide for each magnet system, see Fig 6. The magnets are oriented on the brass rod such that they repel each other, and two aluminum nuts on either side of the magnets constrain the repulsive force, also setting their separation distance. This design allows us to adjust the distance between the magnets and the geometrical center of the magnet assembly.

Each coil former was made from a standard 1 inch PVC water pipe with end caps glued to it. Any nonmagnetic, rigid, cylindrical body will suffice in serving as the coil former. The coil was manually wound onto the PVC pipe using a really low speed lathe spindle and each layer of wire was potted with spray glue. We chose to use AWG¹²-35 wire with 400 windings. In our system, a current of 2.7 mA generated about 0.1 N of force. The total resistance of the wire was 450 Ohms.

The coils can be constructed without a lathe by either hand winding or by using a battery powered drill. Using a lathe to turn down the PVC pipe is an optional step which we chose to execute because it allowed a deeper groove for more windings. Increasing the number of windings on the coil increases the vertical electromagnetic force generated, hence increasing the BL factor.

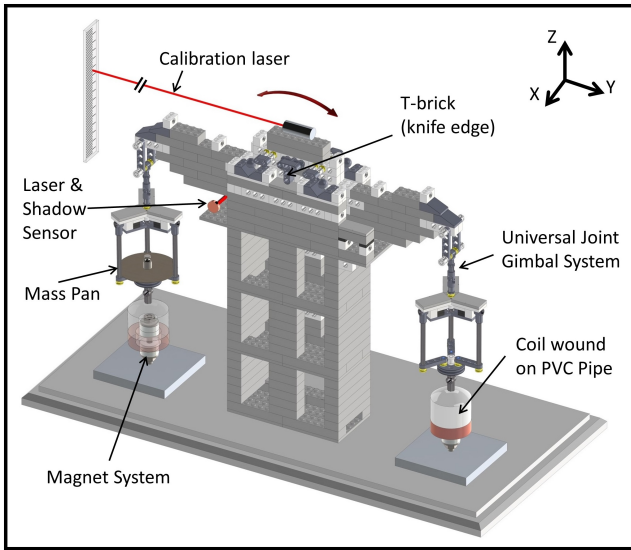


FIG. 2. CAD model of the LEGO watt balance. The balance pivots about the T-block at the center. Two PVC endcaps with copper windings hang from universal joints off either side of the balance beam. Coil A is on the left and Coil B is on the right. A 10 gram mass sits on the Coil A mass pan and each coil is concentric to its own magnet system. Two lasers are used to calibrate and measure the linear velocity of each coil.

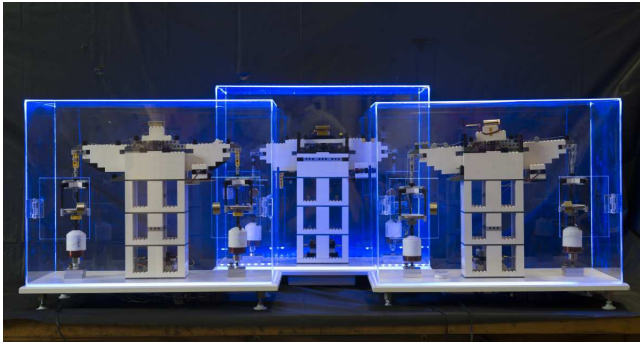


FIG. 3. Image of three similar versions of the LEGO watt balance. The acrylic cases are backlit with blue LEDs and serve the purpose of blocking out disturbances from air currents. Two hinged doors on the front panel allow for small masses to be placed and removed from the mass pans. All the electronics are mounted below the wooden base board. Four adjustable feet are used for leveling the balance.

A small hole was drilled into each end cap top where a LEGO cross axle was attached vertically, allowing each to hang rigidly beneath their corresponding mass pan, see Fig. 6. The mass pan was suspended from three rigid rods linking to a LEGO universal joint (part no. 61903). This dual-gimbal system hangs from a set of two freely pivoting axles parallel to the central pivot (part no. 4208204) connecting to the balance arm. The central pivot (T-brick part no. 4211713) has a "knife edge" radius of approximately 3.1 mm and rests on a smooth

surface.

The whole balance measures approximately 43 cm x 36 cm x 10 cm and has a mass of 4 kg, including the wooden base board. Fig. 3 shows a photograph of three balances that we have built.

IV. ELECTRONICS AND DATA ACQUISITION

We employ two USB devices, a U6 from Labjack and a 1002.0 from Phidget, to connect the LEGO watt balance to a laptop. The U6 is used to measure the position of the balance beam, the induced voltage, and current in each coil. We connect a sixth input to a LEGO handheld controller (potentiometer) that allows students to manually tare the balance, making a popular exhibit at science fairs and demonstrations. The 1002.0 is a four channel analog output that is capable of producing a voltage between -10 V and +10 V. Each channel can source up to 20 mA. One channel is used for each coil. One channel is connected to a double throw, double pole relay. This relay allows the analog output to disconnect from either coil. One coil serves as the sine-driven actuator while the induced voltage can be measured in the other. The relay toggles between the two coils allowing the operator to select which one is the driver. The last output channel is used to remove the bias voltage in the photodiode, later explained in Section V. This allows the use of a smaller gain setting on the analog-to-digital converter that reads the photodiode.

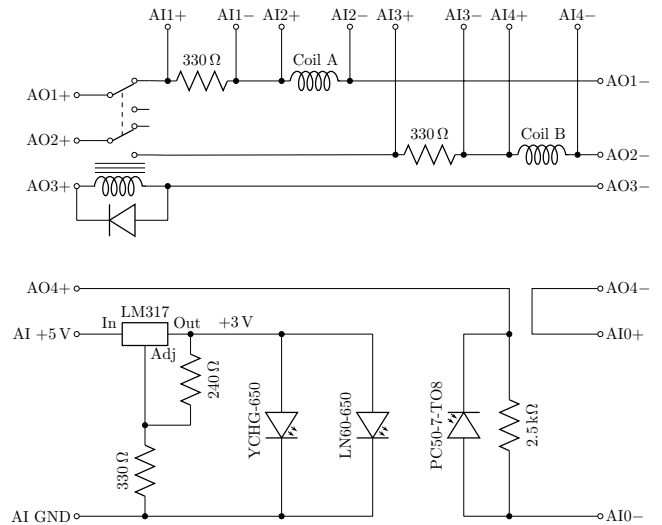


FIG. 4. Circuit diagram for the LEGO watt balance. The top diagram connects one of the two coils to the analog output via a double pole relay. The bottom diagram shows the power supply for the two laser pointers and the photodiode.

We designed the circuit to keep the part count low (seven resistors, one relay and one voltage regulator) to allow for easy construction. Figure 4 shows the circuit

diagram. The top circuit shows the circuitry that is used to measure the induced voltage and current in each coil.

The circuit on the bottom left provides the 3V for the two laser diodes (see Section V for functions of the optical system). The diagram on the right reads the position in the following way: The photo current produced in the photodiode is proportional to the balance position. The photo current flows through R5, the 2.5 k Ω resistor. The voltage drop across R5 is added to the analog output voltage produced by AO4 and the sum is measured. By setting AO4 negative, 0V can be obtained when the balance is at the nominal weighing position.

A custom executable program has been designed to control the LEGO watt balance. If interested in obtaining the free file, please contact the authors. A screenshot of the main interface is shown in Figure 5.

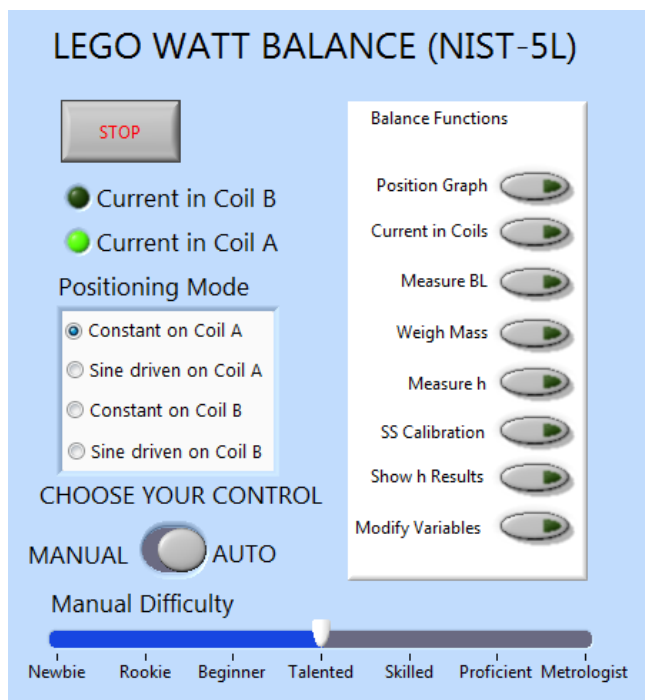


FIG. 5. The front panel of the LEGO watt balance control system. It allows the operator to calibrate the system, weigh small masses, and measure Planck’s constant. Weighing mode (or force mode) can be done either automatically or manually.

V. MEASUREMENT

Like a real watt balance, the LEGO watt balance must undergo a series of alignments and calibrations prior to starting the experiment, detailed in the following 4-step procedure. It is important to calibrate and sense the balance’s angular position. Again, there are many ways to achieve this, but here the angular position of the balance was monitored using a shadow sensor. The system consists of a laser pointer and a photodiode near

the lower edge of one arm of the balance. When the balance moves, it gradually obstructs the optical path of the laser, thereby changing the intensity of light hitting the photo detector. A second laser pointer mounted on top of the balance serves as an optical lever for calibrating the shadow sensor, as we will describe shortly.

Once these prerequisites have been achieved, a complete determination of mass or the value of the Planck constant can begin using a common A-B-A measurement technique. This repetition method is used to cancel the time dependent drift associated with measurements. For instance, one can interleave velocity mode, then force mode, then velocity mode again. Ideally, these measurements are done such that the instrument undergoes as little change as possible, or by performing the measurements in quick succession, neither moving nor tinkering with the balance in between measurements. Once properly aligned and calibrated, a full determination of h through measuring $(BL)_V$ in velocity mode and $(BL)_F$ in force mode is possible. For reference, our experienced operators could perform the following alignment, calibration, and measurement procedure in about 10 minutes.

A. Calibrating the Shadow Sensor

If a shadow sensor and optical lever are indeed chosen for position sensing, a four step process is advised to prepare the balance for calibration.

1. Place the LEGO watt balance on a flat, level surface located a few meters from a wall or vertical structure.
2. Shine the laser pointer mounted on top of the balance to a wall a few meters away as in Figure 6. Ideally, a ruler or grid paper is taped to the wall where the laser spot is located. Measure the distance d from the pivot point of the balance to the wall.
3. Align the balance beam to its support tower. Ensure the balance is not rotated around the y and z axes (the coordinate system is shown in Fig. 2). Looking from the top, the clearances between the beam and the support tower have to be evenly spaced on each side. Our version of the balance has several auxiliary parts, i.e., balls that engage in V grooves and a swivel bracket that aid in the balance alignment. However, it is also possible to perform alignment without these parts. Also, it is good practice to check if the balance is fairly leveled when absent of masses.
4. Concentrically align each magnet system to its corresponding coil. Each magnet system is mounted on X-Y adjustable plates that may slide around until concentricity is reached. Each plate can be clamped down after. It is very important to ensure the coils are not touching the magnets.

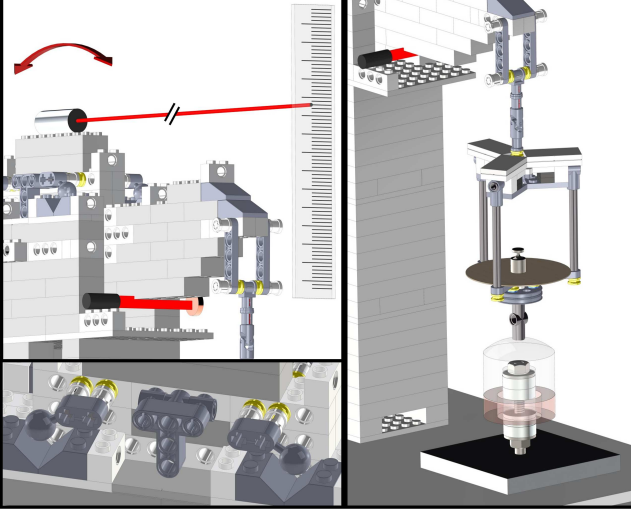


FIG. 6. Top left: The calibration laser projects onto a ruler a few meters away. The shadow sensor detects angular motions of the balance and outputs an oscillatory voltage signature. Right: A transparent view of the PVC coil assembly show its concentricity to the magnet assembly. A stainless steel 10 g mass sits centered on the mass pan and the gimbal system above the mass pan is shown. Bottom left: LEGO T-block serving as the central pivot with balls and V-blocks for kinematic realignment. An identical set exists on the opposing face of the balance.

After these four alignment steps, the instrument is ready for calibration. The balance was servoed to a few different angles which caused the shadow sensor to detect differing light intensities and convert them into voltages V_i . For each voltage, the position of the light spot on the ruler was measured, x_i . In addition to these points, we noted the position of the light when the balance was horizontal, x_0 . The balance angle was then determined as $\theta_i = (x_i - x_0)/d$ and the coil height was calculated by multiplying the balance angle by the effective radius, or $z_i = r_{\text{eff}}\theta_i$. The effective radius was found by measuring the distance of the knife edge to the mass pan universal joint. For the balance described here, $r_{\text{eff}} = 175$ mm.

The optical sensing described here was contrived to drive the measurement uncertainty down to reach our 1% goal. If this goal is not required, easier methods can be used, i.e., directly measuring the coil height for differing servoed positions.

Within reasonable range, the voltage produced by the shadow sensor is a linear function of the coil height. Hence the coil height can be obtained as $z(V) = b(V - V_o)$. A best-fit line to the data (z_i, V_i) yields b and V_o . Fig. 7 shows an example of such a calibration.

B. Velocity Mode Measurement

As stated before, velocity mode measurement ($(BL)_V$) is the key for characterizing the electromagnetic properties of the balance, the first measurement step toward

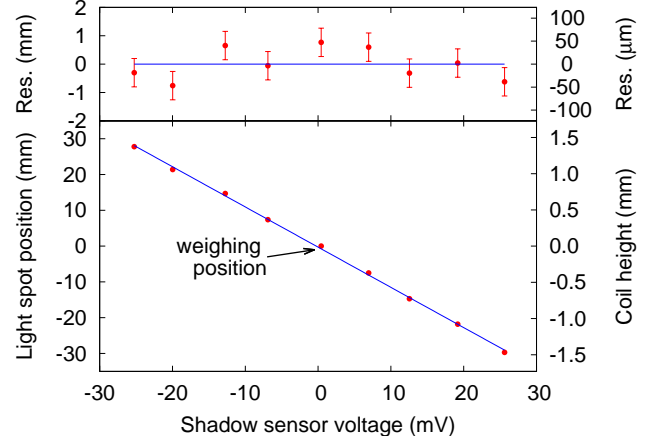


FIG. 7. Calibration of the shadow sensor. The balance is servoed to 9 different shadow sensor voltages. For each voltage, the position of the light spot on the wall, in this case 3489 mm away, is measured. The relationship between the position of the light spot on the wall and the shadow sensor voltage is almost linear. The solid line indicates the best fit to the data. The upper graph shows the residuals between the fit line and the measured data points. We attributed an uncertainty of 0.5 mm (represented by the error bars) to the position measurement of the light spot. Judging from the residuals, that seems reasonable.

obtaining an h or mass value. Our chosen method was to use the information from our calibrated optical displacement sensor and simply take its time derivative to calculate velocity.

If one wants to perform a watt balance experiment using Coil A, then Coil B will be used to drive the balance in a sinusoidal motion, see Fig. 2. Again, there are many ways to actuate velocity mode. We chose a symmetric design such that either arm could be the driver, but other ideas such as a LEGO miniature piston engine have also been tried.³ Because we arbitrarily chose Coil A as the measurement coil and Coil B as the driver, we will continue this nomenclature for consistency and clarity. Using the language of control theory, Coil B was the input driven with an open-loop sinusoidal voltage, and the output balance position was detected by the shadow sensor.

A sinusoidal driving signal with an amplitude of 1 mm and a period of 1.5 s seemed to be a good starting point for our balance. We sampled the Labjack analog input device at a rate of $\Delta = 1$ ms and obtained values for the induced voltage on Coil A, $V(i\Delta)$, and the shadow sensor voltage $V_{SS}(i\Delta)$, where i is the sample number. The coil position $z(i\Delta)$ was then extracted from the shadow sensor voltage. The sampled data were filtered and the coil velocity was obtained by taking the numerical derivative:

$$v(i\Delta) = \frac{z((i+1)\Delta) - z((i-1)\Delta)}{2\Delta} \quad (14)$$

For the pairs of voltages and velocities, a best-fit line

was calculated whose slope was $(BL)_V$. For simplicity, we assumed that $(BL)_V$ did not vary significantly along the coil's trajectory. Since the coil only moved 2 mm, this seems like a reasonable assumption.

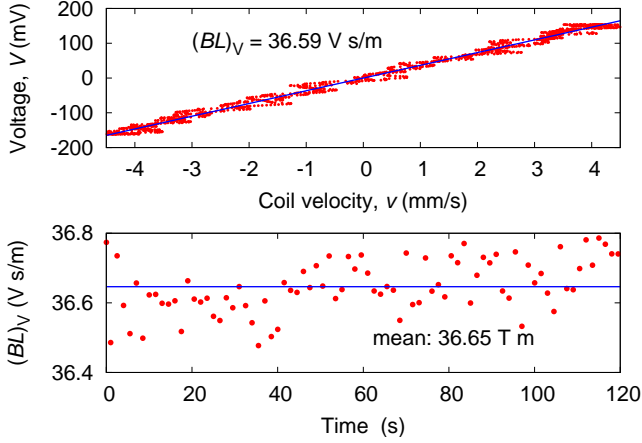


FIG. 8. The top graph shows the measured voltages and velocities of the coil for one period, i.e., 1.5 s. The slope of the solid line which is the best linear fit to the data gives the measured flux density, BL . The bottom graph shows the result of 80 such determinations. The relative standard deviation of the data is 0.2%. For a possible explanation of the small drift see the main text.

The top graph in Fig. 8 shows the measured values of the induced voltage versus the calculated coil velocity. The data shown were taken during one period of a sinusoidal motion of the coil. The slope¹³ is 36.59 V s/m. The bottom graph shows 80 of these measurements for a total of 120 s. A value of the flux integral is determined every motion period. From this data we obtained a mean value

$$(BL)_V = 36.65 \text{ V s/m}. \quad (15)$$

The relative standard deviation of the data was 0.21%. A small relative drift of $3 \times 10^{-5} \text{ s}^{-1}$ was apparent in the data. This drift can be explained by small temperature changes of the magnet. The remanent field of Nd-FeB magnet changes relatively by about $-1 \times 10^{-3} \text{ K}^{-1}$. Hence, a temperature change of 0.01 K/s would explain the observed drift in the $(BL)_V$. Here, we ignore the observed drift of the $(BL)_V$ and assign the mean value.

C. Force Mode Measurement

Coil A is used in the force mode to apply an electromagnetic force to one arm of the balance. The force is easily created by running a current through the coil, but keep in mind that the magnitude must be controlled somehow to hold the balance in its null position after masses are added or removed from the mass pans. The most direct way to control the current is to simply watch

the balance and manually adjust the magnitude of the current until balance is restored. This option is available using the LEGO potentiometer. Simply connect the potentiometer to the coil in series with a battery to form a closed circuit. The projected laser spot on the wall, used to calibrate the shadow sensor, can be used as a target for restoring the balance by manually adjusting the potentiometer.

For users more familiar with control theory and application, the manual feedback can be automated to achieve more consistent results. For instance, the output position of the balance can be detected by the shadow sensor and employed as the control variable for an analog or digital controller. The measured position is then continuously compared by the controller to a desired position, or set-point (typically a null position), and the error between the two used to continuously update the current input to Coil A. In our system, a digital feedback control software tool operates on the data acquisition and control laptop. The controller generates currents that are proportional to the measured error and the integral and derivative of this value with respect to time. Such a scheme is referred to as PID control, where the acronym stands for proportional, integral, and derivative control.

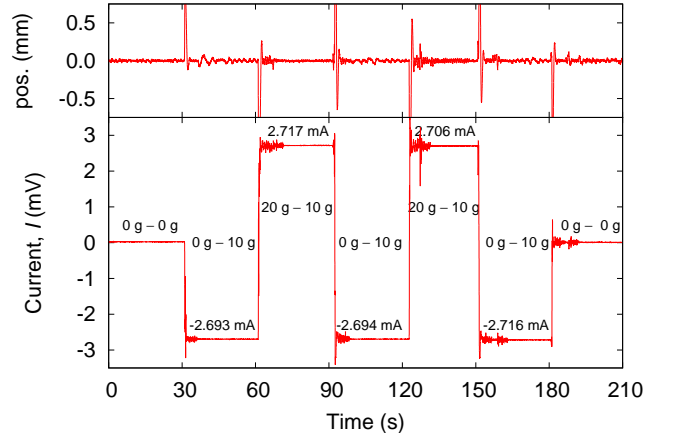


FIG. 9. Force mode in the time domain. The lower graph shows the current required to maintain the balance at a nominal position for seven different load states. The load states are abbreviated by differences. The minuend denotes the mass on the mass pan above Coil A and the subtrahend the mass on the mass pan above Coil B. The mass difference multiplied by the local gravitational acceleration is equal to the force produced by the coil. The software PID controller that is used to servo the balance employs two different gain settings. The change in noise in the measured current occurs when the gain is switched. The top graph shows the position of the coil as a proxy of the balance angle. Adding and removing a mass leads to a spike in position up to 2 mm. The servo quickly reestablishes the nominal weighing position.

Fig. 9 shows the measurement sequence in the force mode. In this example, the measurement was performed in seven steps, each lasting 30 s. The steps were:

1. Both balance pans are empty and the current required to hold the balance at its weighing position is small.
2. A tare mass $m_T = 10\text{ g}$ is added to the pan above Coil B. The exact mass is irrelevant as it will drop out in the final equation. The current $I_1 = -2.693\text{ mA}$ is necessary to maintain the balance position. The current is given by

$$I_1(BL)_F = -m_T g \quad (16)$$

3. The calibrated mass, here $m = 20.2\text{ g}$ is added to the pan above Coil A. This time a positive current, $I_2 = 2.717\text{ mA}$ is required to servo the balance. The equation

$$I_2(BL)_F - mg = -m_T g \quad (17)$$

describes this weighing. Subtracting Eq. 17 from Eq. 16 is sufficient to get an estimate of $(BL)_F$,

$$mg = (I_2 - I_1)(BL)_F \implies (BL)_F = \frac{mg}{I_2 - I_1}. \quad (18)$$

However, to cancel out drift and to get an idea how big the drift is, it is always a good idea to perform a couple more weighings.

4. Another weighing with the Coil A calibrated mass removed determines

$$I_3(BL)_F = -m_T g \quad (19)$$

5. A second weighing with the Coil A calibrated mass added to the pan yields

$$I_4(BL)_F - mg = -m_T g \quad (20)$$

6. A third weighing with the Coil A calibrated mass removed yields

$$I_5(BL)_F = -m_T g \quad (21)$$

7. We finally remove both masses and check if the balance is back at the nominal position and if the current to servo the balance with no weights on either pan remained stable.

Using the above observations, the following number can be calculated:¹⁴

$$I = -\frac{1}{3}(I_1 + I_3 + I_5) + \frac{1}{2}(I_2 + I_4) = \frac{mg}{(BL)_F}. \quad (22)$$

In order to obtain the flux integral from the force mode, one needs the local acceleration g . Your local gravitational acceleration can be obtained by a website provided by the National Oceanic and Atmospheric Administration.¹¹ For our geographical coordinates at NIST

Gaithersburg (Latt: 39.1261°N , Long: 77.2211°W , Elevation: 124.304 m), the website yielded $g = 9.80103\text{ m/s}^2$ with a relative uncertainty of 2×10^{-6} . This uncertainty was well below what we needed for a 1 % level measurement.

With the above numbers and $m = 20.2\text{ g}$, we obtain

$$(BL)_F = \frac{0.0202\text{ kg} \cdot 9.80103\text{ m/s}^2}{0.0054125\text{ A}} = 36.58\text{ N/A} \quad (23)$$

D. A value for h

To obtain a value for h , the ratio of the flux integrals obtained in force mode and velocity mode was multiplied with h_{90} as described in Eq. 12. Here, we obtained

$$h = h_{90} \frac{36.58\text{ N/A}}{36.65\text{ Vs/m}} = 0.998 h_{90} = 6.61 \times 10^{-34}\text{ Js}. \quad (24)$$

Determining a value for the Planck value is half the work, the other half is to determine the measurement's uncertainty. We believe that a measurement of the Planck constant with a relative uncertainty of 1 % or less is possible with this LEGO watt balance.

For example, a source of uncertainty that is easy to understand comes from the distance measurement in the shadow sensor calibration. If the distance from the laser diode to the wall is measured with a measuring tape to be 3000 mm with an uncertainty of 3 mm , the relative uncertainty associated with this would be 0.1% . Also, if the laser spot oscillates by $\pm 30\text{ mm}$ and the vertical wall ruler can be read with an uncertainty of 0.5 mm , then the relative uncertainty of this measurement is 0.83% ; clearly the dominant source of uncertainty in this measurement. A good metrologist will identify the largest sources of uncertainty and will try to improve these.

Large contributors to bias and measurement uncertainty are offset forces produced in the large-radius knife edge, parasitic motions of the coil during velocity mode, and horizontal forces in force mode which arise from misalignments. An uncertainty analysis is beyond the scope of this article and we leave this exercise to the interested reader. Several inspiring articles can be found in the literature that provide details on how to assess these uncertainties, e.g.,^{15–17}.

E. Summary

In 2013 and 2014, we built three LEGO watt balances. These balances were demonstrated and received with enthusiastic responses in science fairs, classrooms, and to visitors coming to NIST. It prompted us to write this article to promote building these devices for STEM education. What will building and operating such a device accomplish?

In the new SI, the kilogram will be defined via fixing the Planck constant,[?] a way that is not apparent to most

people. The LEGO watt balance makes it easy to demonstrate how electrical power relates to mechanical power. One can show how it is possible to generate a mechanical force, whose value is precisely given by electrical measurements. Unfortunately, it still requires some abstraction to explain how electrical power is related to the Planck constant via the Josephson effect and the quantum Hall effect. But once that bridge is crossed, the relationship between mass and h can be made clear. From there, it is fascinating to ponder the implications of mass redefinition: the Planck constant, a natural constant found in quantum mechanics, can be directly used to define mass on a macroscopic scale! At any scale for that matter. With this LEGO balance, one can determine the absolute mass of small objects without needing any comparison or traceability to reference mass standards. It is also interesting to consider how an apparatus assembled from plastic bricks can measure h in a classroom or living room setting with an uncertainty of only 1 %.

Besides these top level concepts, we found that the LEGO watt balance provides ample teachable moments. For example, our balance controls included a “Manual Feedback” option where an operator could rotate a potentiometer to try and null the balance. Most people who find it hard to control the balance are amazed to see how seemingly effortless and precise a PID controller achieves the task. This provides a nice segue into control theory.

Closer to home, questions at the heart of metrology such as “How does one measure something? Are all measurements comparisons? What is accuracy and precision? What is the error in this measurement?” arise while constructing such an experiment. Answering these questions is an opportunity to teach the audience about the impor-

tance of measurements for society. Getting the audience interested in metrology is the intent of the LEGO watt balance. This goal has been achieved every single time we demonstrated our glorified toy.

Finally, the LEGO watt balance combines three important things: Science, Technology, and Fun.

Appendix A: List of parts

The table below shows the components necessary to rebuild our LEGO watt balance. Each LEGO engineer or scientist is encouraged to explore other building materials and components to create an even more optimized and personalized instrument. The prices are accurate as of early 2014 and do not include shipping and handling. Though we spent over \$600 on this project, it can be built for significantly less. For example, one may choose to implement the National Instruments NI-USB-6001 (\$189), which replaces both the Labjack DAQ and the Phidget Analog 4 Output, driving the total cost down by \$200. The internet addresses and vendors are suggestions and are not endorsed.

In addition to the parts listed below, a wooden base, wire to connect the electrical circuits, and wires to wind the coil are required and can be purchased from a variety of vendors.

Construction and commissioning a LEGO watt balance can be a lot of fun. Such a project is suitable for small teams at school or a private person at home. We hope to encourage many enthusiasts to build a LEGO watt balance and have fun toying with the science of measurement, metrology.

Part Name	Part No.	Quantity	Total Price (\$)
Custom LEGO Watt Balance Software		1	Free
Order from http://shop.lego.com/en-US/Pick-A-Brick-ByTheme			
Brick 2x4	300101	65	19.50
Brick 2x8	6033776	73	36.50
Brick 1x2 with cross hole	4233487	12	4.20
T-Beam 3x3 w/hole O4.8	4552347	2	0.60
Technic Brick 1x2 O4.9	370026	16	2.40
Technic Brick 1x4 O4.9	4211441	48	12.00
Technic Brick 1x6 O4.9	389426	2	0.80
Technic Brick 1x8 O4.9	4211442	2	1.00
Technic Ang. Beam 3x5 90 Deg.	4211713	2	0.60
Plate 8x8	4210802	9	9.90
Plate 1x2	4211398	15	1.50
Plate 1x4	4211445	10	1.50
Plate 2x3	4211396	6	1.20
Cross Axle 2M W. Groove	4109810	8	0.80
Cross Axle 3M	4211815	6	0.14
Cross Axle 5M	4211639	6	1.20
Cross Axle 8M	370726	8	1.60
Bush for Cross Axle	4211622	14	2.10
1/2 Bush for Cross Axle	4211573	28	2.80
Double Bush 3M O4.9	4560175	6	1.20
Roof Tile 2x2/45 deg	303926	4	0.80

Part Name	Part No.	Quantity	Total Price (\$)
Roof Tile 2x2/45 deg Inv.	366026	2	0.40
Roof Tile 2x3/25 deg	4211106	2	0.40
Roof Tile 2X3/25 deg Inv.	374726	4	0.80
Connector Peg W. Friction 3M	4514553	8	2.00
Connector Peg/Cross Axle	4666579	6	0.60
Catch w. Cross Hole	4107081	8	1.60
Flat Tile 2x4	4560178	4	1.20
Hinge 1x2 Lower Part	383101	6	1.50
Hinge 1x2 Upper Part	6011456	6	1.50
Double Conical Wheel Z12 1M	4177431	2	0.60
Angle Element, 180 Degrees [2]	4107783	2	0.40
Order from http://atrbricks.brickowl.com/			
Technic Beam 1 x 4 x 0.5 with Boss	2825 / 32006	6	0.30
Technic Beam 2 Beam w. Angled Ball Joint	50923 / 59141	2	0.13
Beam 5 x 0.5	32017	4	0.76
Wedge Belt Wheel	2786 / 4185	4	1.00
Gear with 8 Teeth (Narrow)	3647	2	0.20
Order from http://thatspecialbrick.brickowl.com/			
Universal Joint	61903	2	0.94
Gear with 16 teeth	94925	2	0.34
Bevel Gear with 12 teeth	6589	2	0.26
Order from http://www.labjack.com/u6			
Multifunction DAQ with USB - 16 Bit	U6	1	299.00
Order from http://www.phidgets.com			
PhidgetAnalog 4 Output	1002_0	1	90.00
Order from http://www.apinex.com			
Focus. Line Red Laser Module <1mW	YCHG-650	1	15.00
Line Laser Module (650nm) - <1mW	LN60-650	1	15.00
Order from http://www.mouser.com			
Photodiode 7.98mm Dia Area	718-PC50-7-TO8	1	61.63
USB Cables USB A - MINI-B	538-88732-8702	1	1.65
USB A-TO-B Shielded 2.09m	538-88732-9200	1	3.58
Low signal Relay	769-TXS2-4.5V	1	4.58
Resistors 240ohms	291-240-RC	1	0.10
Resistors 330ohms	291-330-RC	4	0.40
Resistors 1500ohms	291-1.5k-RC	1	0.10
Linear Voltage Regulators	511-LM317T	1	0.72
Order from http://www.magnet4less.com			
N48 grade - 3/4 (OD) x 1/4 (ID) x 3/4 in. ring magnet NR011-1		4	15.96
Order from http://www.mcmaster.com			
Brass Threaded Rod - 1/4"-20 Thread, 1' length	98812A039	1	2.65
White PVC Pipe Fitting	4880K53	2	1.00
White PVC Unthreaded Pipe	48925K93	1	5.27
Total			633.77

¹ F. Seifert, A. Panna, S. Li, B. Han, L. Chao, A. Cao, D. Haddad, H. Choi, L. Haley, and S. Schlamminger, "Construction, Measurement, Shimming, and Performance

of the NIST-4 Magnet System," *Instrum. Meas., IEEE Transactions on*, vol. 63, no. 12, pp. 3027–3038, 2014.

² B. Kibble, "A Measurement of the Gyromagnetic Ratio of

- the Proton by the Strong Field Method,” *Atomic Masses and Fundamental Constants vol. 5*, ed J H Sanders and A H Wapstra (New York: Plenum), pp. 545–551, 1976.
- ³ T. Quinn, L. Quinn, and R. Davis, “A simple watt balance for the absolute determination of mass,” *Physics Education*, vol. 48, no. 5, p. 601–606, 2013.
 - ⁴ Certain commercial equipment, instruments, or materials are identified in this paper in order to specify the experimental procedure adequately. Such identification is not intended to imply recommendation or endorsement by the National Institute of Standards and Technology, nor is it intended to imply that the materials or equipment identified are necessarily the best available for the purpose.
 - ⁵ J. Clarke, “The Josephson Effect and e/h ,” *American Journal of Physics*, vol. 38, no. 9, pp. 1071–1095, 1970.
 - ⁶ Y. Tang, V. N. Ojha, S. Schlamminger, A. Rüfenacht, C. Burroughs, P. Dresselhaus, and S. Benz, “A 10V programmable Josephson voltage standard and its applications for voltage metrology,” *Metrologia*, vol. 49, no. 6, pp. 635–643, 2012.
 - ⁷ J. Eisenstein, “The quantum Hall effect,” *American Journal of Physics*, vol. 61, no. 2, pp. 179–183, 1993.
 - ⁸ R. E. Elmquist, M. E. Cage, Y.-h. Tang, A.-M. Jeffery, J. R. J. Kinard, R. F. Dziuba, N. M. Dziuba, and E. R. Williams, “The Ampere and Electrical Standards,” *J. Res. Natl. Inst. Stand. Technol.*, vol. 106, pp. 65–103, 2001.
 - ⁹ T. Quinn, “News from the BIPM,” *Metrologia*, vol. 26, no. 1, pp. 69–74, 1989.
 - ¹⁰ B. Taylor and T. Witt, “New international electrical reference standards based on the Josephson and quantum Hall effects,” *Metrologia*, vol. 26, no. 1, pp. 47–62, 1989.
 - ¹¹ “Surface Gravity Prediction by NGS Software Requests.” https://www.ngs.noaa.gov/cgi-bin/grav_pdx.prl.
 - ¹² American Wire Gauge.
 - ¹³ We deliberately use V s/m as units for $(BL)_V$, to emphasize that V is a conventional unit, while m and s are SI units. With the unit of T m, we would have lost this distinction. Analogous, we use the unit N/A for the quantity $(BL)_F$.
 - ¹⁴ H. E. Swanson and S. Schlamminger, “Removal of zero-point drift from AB data and the statistical cost,” *Measurement Science and Technology*, vol. 21, no. 11, p. 115104–115110, 2010.
 - ¹⁵ R. L. Steiner, E. R. Williams, D. B. Newell, and R. Liu, “Towards an electronic kilogram: an improved measurement of the Planck constant and electron mass,” *Metrologia*, vol. 42, no. 5, p. 431–441, 2005.
 - ¹⁶ I. Robinson, “Alignment of the NPL Mark II watt balance,” *Meas. Sci. Technol.*, vol. 23, no. 12, p. 124012–124029, 2012.
 - ¹⁷ B. Kibble and I. Robinson, “Principles of a new generation of simplified and accurate watt balances,” *Metrologia*, vol. 51, no. 2, p. S132–S139, 2014.