

Fourier plane optical microscopy and spectroscopy

Adarsh B. Vasista¹, Deepak K. Sharma¹ and G V Pavan Kumar^{1,2,*}

¹ Department of Physics, Indian Institute of Science Education and Research, Pune-411008, India

² Centre for Energy Science, Indian Institute of Science Education and Research, Pune-411008, India

*E-mail: pavan@iiserpune.ac.in

Abstract

Intensity, wavevector, phase, and polarization are the most important parameters of any light beam. Understanding the wavevector distribution has emerged as a very important problem in recent days, especially at nanoscale. It provides unique information about the light-matter interaction. Back focal plane or Fourier plane imaging and spectroscopy techniques help to measure wavevector distribution not only from single molecules and single nanostructures but also from metasurfaces and metamaterials. This review provides a birds-eye view on the technique of back focal imaging and spectroscopy, different methodologies used in developing the technique and applications including angular emission patterns of fluorescence and Raman signals from molecules, elastic scattering etc. We first discuss on the information one can obtain at the back focal plane of the objective lens according to both imaging and spectroscopy viewpoints and then discuss the possible configurations utilized to project back focal plane of the objective lens onto the imaging camera or to the spectroscope. We also discuss the possible sources of error in such measurements and possible ways to overcome it and then elucidate the possible applications.

Keywords: Back focal plane, Spectroscopy, Single molecule, Light scattering

Submitted to Wiley-VCH Encyclopedia of Applied Physics (<http://dx.doi.org/10.1002/3527600434>)

Table of Contents

1. Introduction.....	2
1.1. What is the back-focal plane of an objective lens?.....	3
1.2. What information can we get using back focal plane of an objective lens?.....	3
1.1.1. <i>Imaging view point</i>	3
1.1.2. <i>Spectroscopy view point</i>	4

2. Imaging the back focal plane.	4
2.1. How to image the BFP using an optical microscope?	5
2.2. What are the possible sources of error?	6
2.2.1. Placement of the relay lenses:	6
2.2.2. Chromatic and spherical aberrations:	6
2.2.3. Choice of objective lens:	6
2.2.4. Placement of associated optics:	6
3. Spectroscopy in back focal plane	7
4. Applications	7
4.1. Single molecule fluorescence imaging	7
4.2. Raman scattering	8
4.3. Plasmonic scattering/ Optical antenna	9
4.4. Secondary emission via optical antenna	11
4.5. Nonlinear scattering	14
5. Conclusions	14
References	15
Suggested further reading	18

1. Introduction

Fourier plane imaging , in other words , back focal plane (BFP) imaging and spectroscopy has emerged as an important tool, of late, in nanophotonics. It provides an alternative set of information which is not available in conventional real plane imaging. In a BFP image, one usually has information encoded in angular co-ordinates in contrast to the real plane image, where the resolution is in spatial co-ordinates. Measuring the angular emission pattern has been utilized in variety of applications including measuring radiation patterns from single quantum emitters (1-5), nanoantennas(6-12), understanding coupling mechanisms in cavities (13), in cathodoluminescence(14, 15), secondary emission(16-21) and much more.

BFP imaging has been utilized in conoscopy (22-24) since early days. But the application of this technique to microscopy is a recent development, which needs attention. The advancements in micro/nanofabrication methodologies have made handling light at subwavelength scales possible. This calls for sophisticated microscopy techniques to understand and harness light at

subwavelength scales. BFP imaging has emerged as an important tool along with scanning near field optical microscopy(25) , super-resolution microscopy(26), dark field microscopy(27), etc. In addition to this, BFP spectroscopy provides wavelength resolution to the angular information of the light emanated from the sample. This technique coupled to a microscope helps to measure energy-momentum dispersion relation of the outcoupled light from structures down to subwavelength limit (13, 28-31).

Here we explain, in detail, the concept of BFP imaging and spectroscopy, different methodologies and techniques involved in BFP measurement with a microscope, and its applications in nanophotonics, optical communication and nonlinear optics.

1.1. What is the back-focal plane of an objective lens?

In a typical microscope, the light emanating from the sample is collected using an objective lens and then routed to an imaging camera through external optics. Now, the light emerging from the sample plane can be expressed as weighted sum of different plane wave components, $u^i(x,y)$, emanating at different angles, θ^i 's. When this composite light passes through the objective lens, the different plane wave components get separated. The objective lens transforms the plane waves to paraboloidal waves which are focused at a point in the back focal plane. Any plane wave arriving at an angle θ_x, θ_y gets focussed onto $(\theta_x f, \theta_y f)$, where f is the focal length of the objective lens. The objective lens maps each direction (θ_x, θ_y) to a single point $(\theta_x f, \theta_y f)$ onto the back focal plane and hence gives information about the angular distribution of the emanated light from the sample plane.

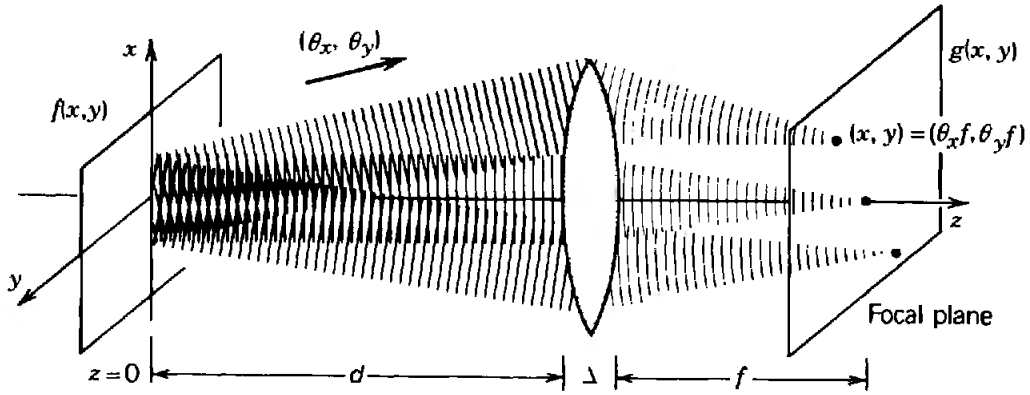


Fig1: Conceptual schematic of Fourier transform using a lens. Reproduced with permission from(32)

1.2. What information can we get using back focal plane of an objective lens?

1.1.1. Imaging view point

The most important information which can be availed is the angular spectrum of the emanated light. In other words, far field emission pattern of the emission from a nano-object. As the angular distribution of light is intimately related to its momentum ($\hbar\mathbf{k}$), the BFP image can be viewed as inplane momentum (k_x, k_y) distribution of the outcoupled light. This enables to measure angle resolved

scattering/emission from individual/array of nanostructures(33-39), perform single molecule orientation imaging and measure far field molecular emission patterns(40-43), designing optical/opto-mechanical sensors(44, 45) etc.

1.1.2. Spectroscopy view point

As the back focal plane has resolution in terms of inplane momentum vectors, one can select specific wavevectors and disperse them to obtain wavelength information. In other words, one can experimentally measure energy-momentum dispersion relation. BFP spectroscopy is an important tool to distinguish competing molecular processes at nanoscale(13), perform angle resolved absorption and emission spectroscopy(12, 46-48) , measuring radiation patterns of Raman modes from a molecular system(49, 50) etc.

2. Imaging the back focal plane.

As discussed above, the back focal plane of the objective lens contains information about the emission wavevectors from the sample plane. Imaging the BFP of the objective lens is tricky as physically the BFP of the objective lens lies near-inside the metallic case. Hence in a typical optical microscope, relay optics is utilized to project the BFP of the objective lens to the imaging camera.

2.1. How to image the BFP using an optical microscope?

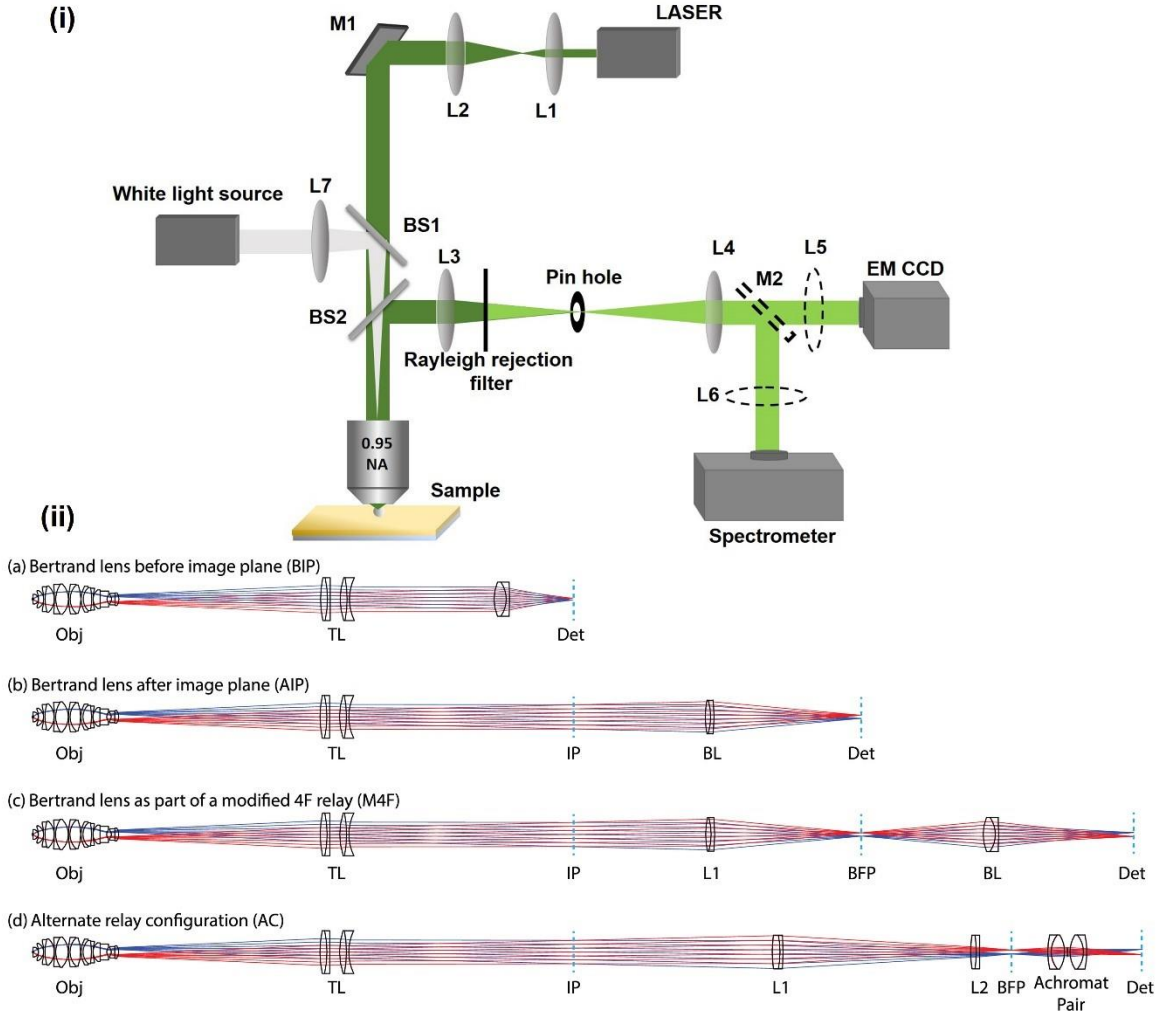


Fig. 2: (i) Schematic of imaging back focal plane using an optical microscope (ii) (a) – (d) Different relay configurations to image the back focal plane of the objective lens. Reproduced with permission from (51)(Optical Society of America).

Figure 2 (i) shows a typical setup used to perform back focal plane imaging and spectroscopy. This is a prototype setup and can be expanded and utilized with modifications. Here, the incoming laser beam is focused onto the sample plane using a high numerical aperture objective lens and the back scattered light is collected using the same objective lens (upright configuration). The collected light is then passed through lenses L3, L4 and L5 and projected to the CCD. The beam which is emanated from the infinity corrected objective lens is focused using lens L3, which creates a conjugate image plane. Free optics like spatial filters can be placed at the conjugate image plane to filter out certain region of the sample. This configuration is very important as it can be extended for confocal microscopy, where a confocal pinhole is placed at the conjugate image plane in the place of spatial filter. Spectral filters like long pass / notch filters can be placed in the optical path to reject the incident laser, for secondary emission related measurements. This conjugate image plane now acts as the object plane for the lens L4, which projects the Fourier transform of the image plane onto its focal plane, where an imaging camera will be placed. Thus the relay

of lenses L3 and L4 projects the BFP of the objective lens onto the imaging camera outside the microscope. Combination of lenses L4 and L5 will create real plane image at the imaging camera.

The configuration explained above is one among many such possible configurations to project the BFP of the objective lens onto the camera as seen in fig 2 (ii). Kurvits et.al.,(51) lists different possible configuration to project BFP onto the camera.

2.2. What are the possible sources of error?

As different optical elements like lenses, mirrors and beam-splitters are placed in the optical path to image the BFP of the objective lens, there are different possibilities of generation of errors. Abberations (chromatic and spherical), vignetting, alignment issues, non-uniform intensity pattern of the BFP etc are possible errors in the imaging of BFP. We try to enumerate some possibilities. The list is not exhaustive.

2.2.1. Placement of the relay lenses:

The tube lens essentially performs inverse Fourier transform and slight misplacement of the tube lens can actually make the phase information of the BFP lost (32). Since we usually measure intensity, where the phase information is integrated out, this will not be visualised. But the placement can have significance when BFP measurements are coupled with phase sensitive measurements. Also the choice of the focal length of the lenses play a crucial role in accurately projecting the BFP to the camera. Choice of very large (small) focal length of the lenses can create problem in relative placement of lens and magnification, thus affecting the formation of the conjugate BFP at the camera.

2.2.2. Chromatic and spherical aberrations:

If the collection optical path is slightly off the optical axis of the lenses, then spherical aberrations can occur. This will affect the formation of BFP at the camera. To avoid this, proper care has to be taken during alignment to ensure that the collected light from the objective lens passes through the optical axis of the relay lenses. In the case of polychromatic light (white light imaging, imaging broadband secondary emission) chromatic aberrations can occur, which can be avoided by the use of achromatic lenses.

2.2.3. Choice of objective lens:

Selection of objective lenses are also crucial in minimising aberrations in the imaging of the BFP. Kurvits et al.,(51) provides a detailed description for the choice of objective lenses for high numerical aperture imaging.

2.2.4. Placement of associated optics:

Placement of bi-refrigrant optics like polarizers can introduce errors in the measurement of BFP. Polarizers, in general, placed anywhere at the optical path can shift the position of BFP on the camera, which creates difficulty in measurements like BFP polarimetry. Hence polarizers need to be placed at the conjugate BFP. This conjugate BFP should be projected onto the camera using another set of relay lenses.

3. Spectroscopy in back focal plane

In the previous sections, we discussed the importance and technique of performing BFP imaging. However when the collected light from the objective lens has specific dispersion relation, then it would be necessary to have wavelength resolution for the wavevector distribution. Spectroscopy in the BFP essentially provides the extra dimension needed for measuring dispersion relation. Here we project the BFP of the objective lens onto the spectrometer after filtering certain part of the BFP. Filtering can be achieved either by introducing a slit at the conjugate BFP(28) as shown in fig.3 or projecting the BFP onto the slit of the spectrometer and utilizing the slit as the filter(13, 30).

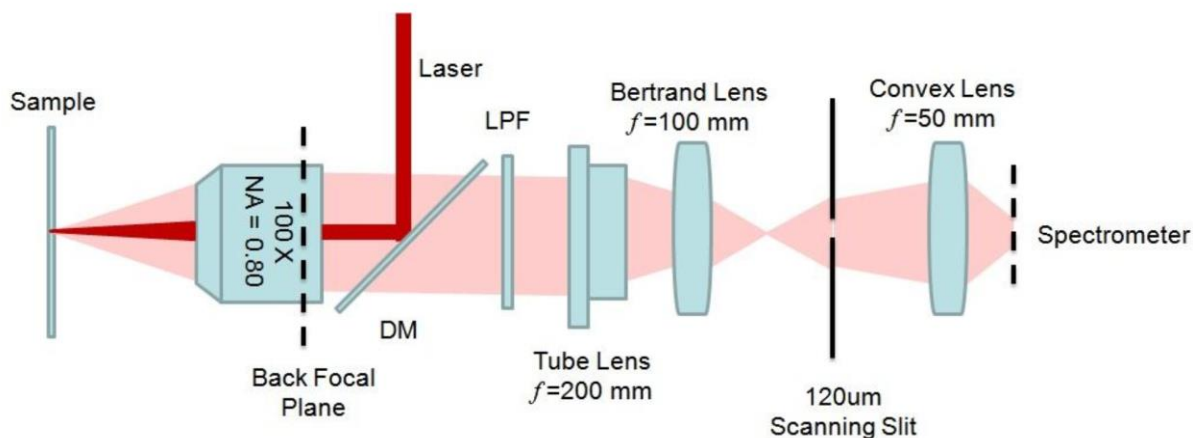


Fig.3: Schematic representation of back focal plane spectrometry. Reproduced with permission from (28)

The information in a typical BFP image contains convolution of all possible wavelengths with angular resolution. If a process contains multiple modes of emission, like SERS, which can possibly occupy different wavevectors can only visualized by performing BFP spectroscopy (18, 52). Also such technique can be expanded to understand competing processes, Raman and fluorescence, from an ensemble of molecules(13) as these processes happen to emit at the same wavelength and cannot be separated in the real plane (in continuous wave limit). This provides an added advantage over real plane spectroscopy technique. The possible sources of error enumerated in the last section equally applies to the BFP spectroscopy also. In addition, the resolution of the spectrum now will also depend on the width of the slit used to filter the wavevectors in the BFP. Hence there will be a trade-off between the throughput and the resolution, which depends on the size of the slit, as in the case of real plane spectroscopy.

4. Applications

4.1. Single molecule fluorescence imaging

Back focal plane imaging is utilized extensively in molecular emission imaging and spectroscopy as intensity pattern in the back focal plane image is intricately related to the molecular orientation. Lieb et.al.,(3) has shown that back focal plane imaging can be utilized to probe single molecule orientations. By understanding the radiation patterns of molecular dipole single molecular orientations can be probed, as shown in figure 5.

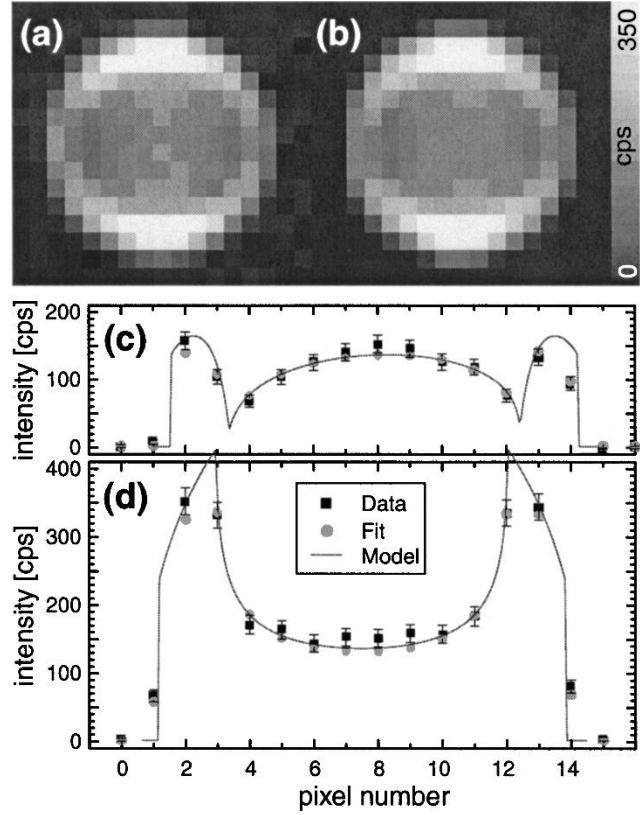


Fig.5: (a) Measured single molecule fluorescence emission pattern captured by imaging the back focal plane of the objective lens. (b) Fitted emission pattern of the molecular emission. (c) and (d) are intensity cross sections along the horizontal and vertical line through the center of the pattern respectively. Reproduced with permission from (3)(Optical Society of America)

4.2. Raman scattering

Raman scattering provides unique information about the molecule. Angular distribution of the Raman scattering signals carry vital value, as one can determine molecular bonding orientation, interaction strength etc. Especially in anisotropic 2D materials like Graphene and Phosphorene, BFP images of different Raman modes is of crucial value. Budde et al.(18), have probed the angular radiation patterns from Graphene and demonstrated the important differences between the G and 2D bands of the material. BFP imaging coupled with output polarization resolution provides information about the polarizability of the material in addition to the angular radiation pattern of the modes of emission as shown in figure 6.

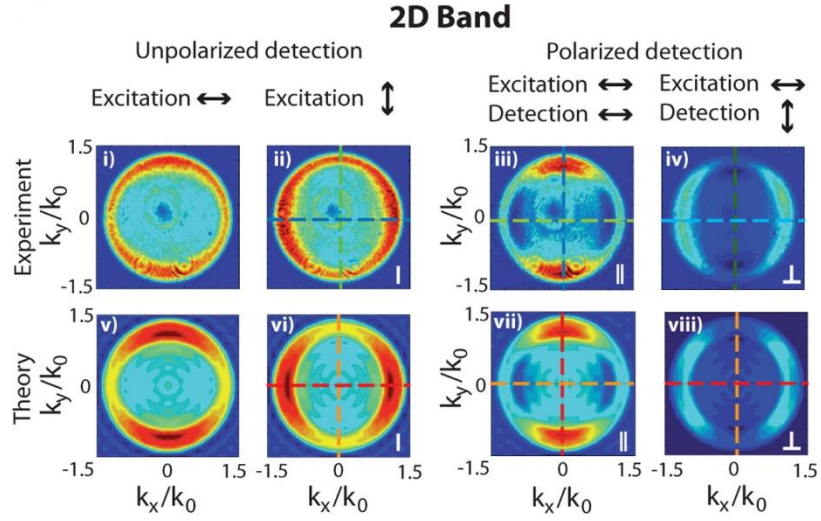


Fig.6: Back focal plane images captured by spectrally filtering 2D band of graphene, showing strongly polarized emission. Reproduced with permission from (18)

4.3. Plasmonic scattering/ Optical antenna

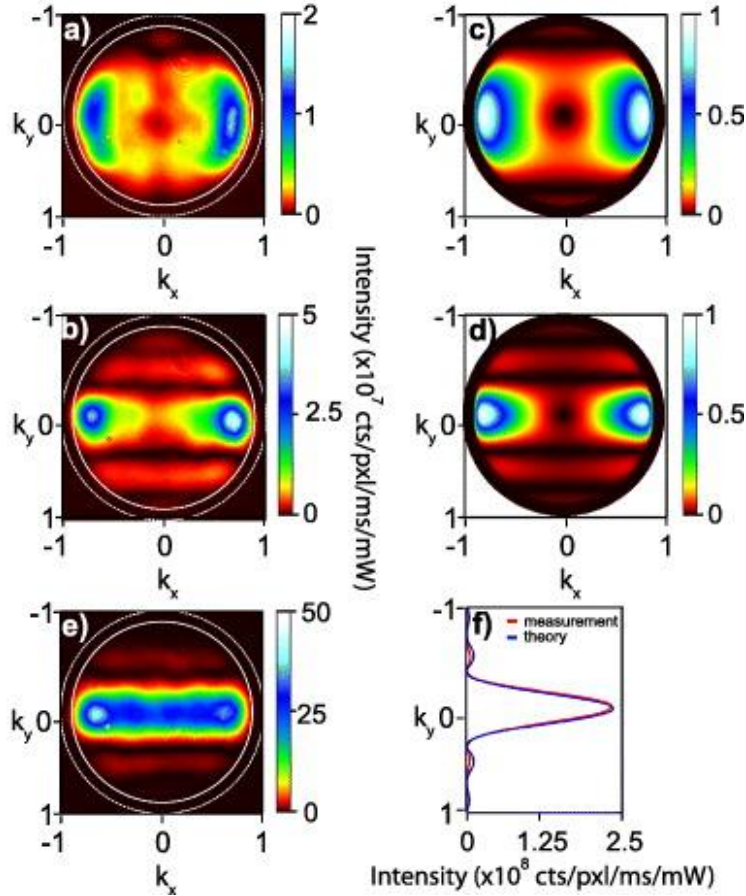


Fig.7: Back focal plane images captured from single plasmonic scatterers using evanescent field excitation. (a) and (b) are measured back focal plane images from 1 μm long and 2 μm long Au nanorods placed over glass substrate respectively . The width of the rod was 50 nm and excitation was p-polarized (c) and (d) are calculated BFP patterns. (e) measured BFP image for a 2 μm long, 50 nm wide Au nanorod on glass substrate with s-polarized excitation. Reproduced with permission from (6)

Farfield scattering of light from nano objects helps in understanding various optical phenomena at subwavelength scale which can further help in accommodating nanostructures in an integrated nano photonic device. This is essential as one of the important task in photonic circuits is to direct or guide light in particular direction and hence knowing the wave vector distribution is necessary. In the past, back focal plane microscopy has been utilized to understand elastic scattering from various nanostructures (39, 53-59). Sersic et.al (6) has shown how farfield scattering profile from a nanorod changes as a function of polarization of light and length of the nanorod, as shown in figure 7. It is clear from back focal plane images that farfield pattern changes drastically with the change in either geometry of nanorod or polarization of light.

On the other hand, Shegai et.al., (8) probed wave guiding properties of Ag nanowire and analyzed the outcoupled light using back focal plane microscopy. Back focal plane image in figure 8(c) reveals that outcoupled light from distal end of wire is unidirectional which is not evident by just looking at real plane propagation image of wire given in figure 8(b). Spread in the back focal plane image tells about the divergence of out coupled light which is an important information about a photonic device.

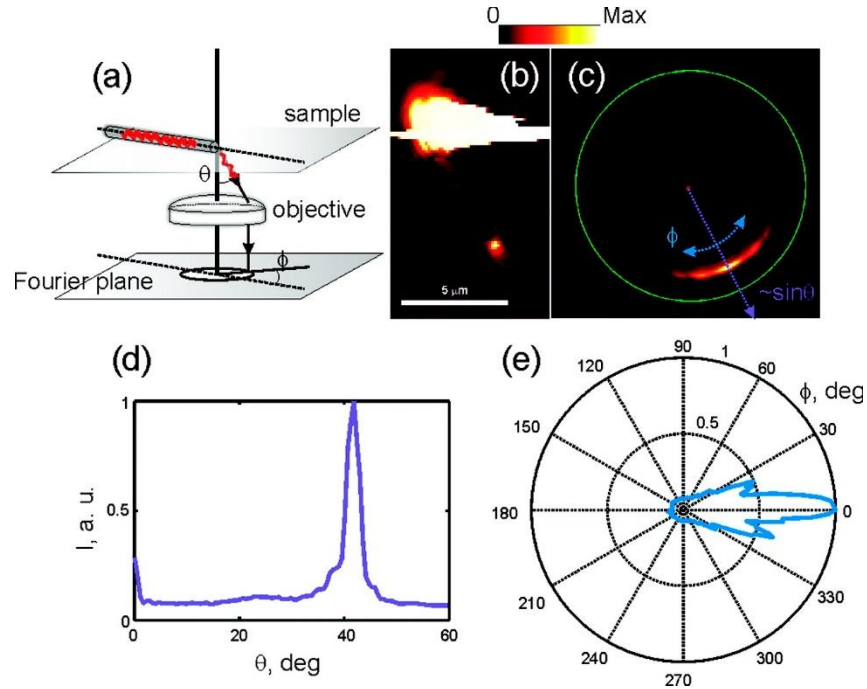


Fig.8: Unidirectional emission from single Ag nanowire. (a) Schematic of the emission configuration. (b) Real plane image showing the distal end of the wire. (c) Back focal plane image captured from distal end

of an individual Ag nanowire. (d) and (e) are intensity profiles plotted along the radius of the BFP image and circumference at the maximum emission angle. Reproduced with permission from (8)

4.4. Secondary emission via optical antenna

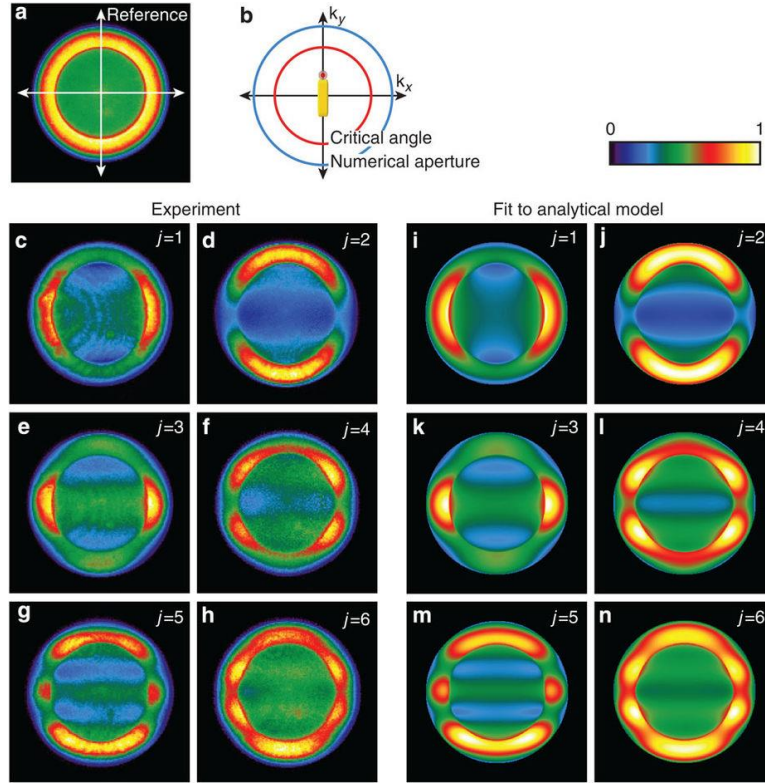


Fig.10 : Back focal plane fluorescence images captured from single quantum dot coupled to plasmonic nanorod. (a) Reference structure (b) Schematic of an angular pattern highlighting the critical angle of the glass–air interface and the numerical aperture of the objective (maximum collection angle). (c–h) Experimental patterns of the first six resonant antenna modes. (i–n) The patterns are reproduced with the analytical resonator model, including a rotationally symmetric background. Reproduced with permission from (1)

We discussed about probing molecular emission process using back focal plane imaging. This technique has been extended quite extensively to study molecular emission coupled to nanostructures. BFP imaging technique is of critical importance in probing antenna effects in plasmon coupled molecular emission. Recent reports (60-63) show the importance of BFP imaging in studying plasmon coupled molecular fluorescence. BFP imaging not only provides information about the emission properties of the molecule but also that of the plasmonic nanostructure, as studied by (1).

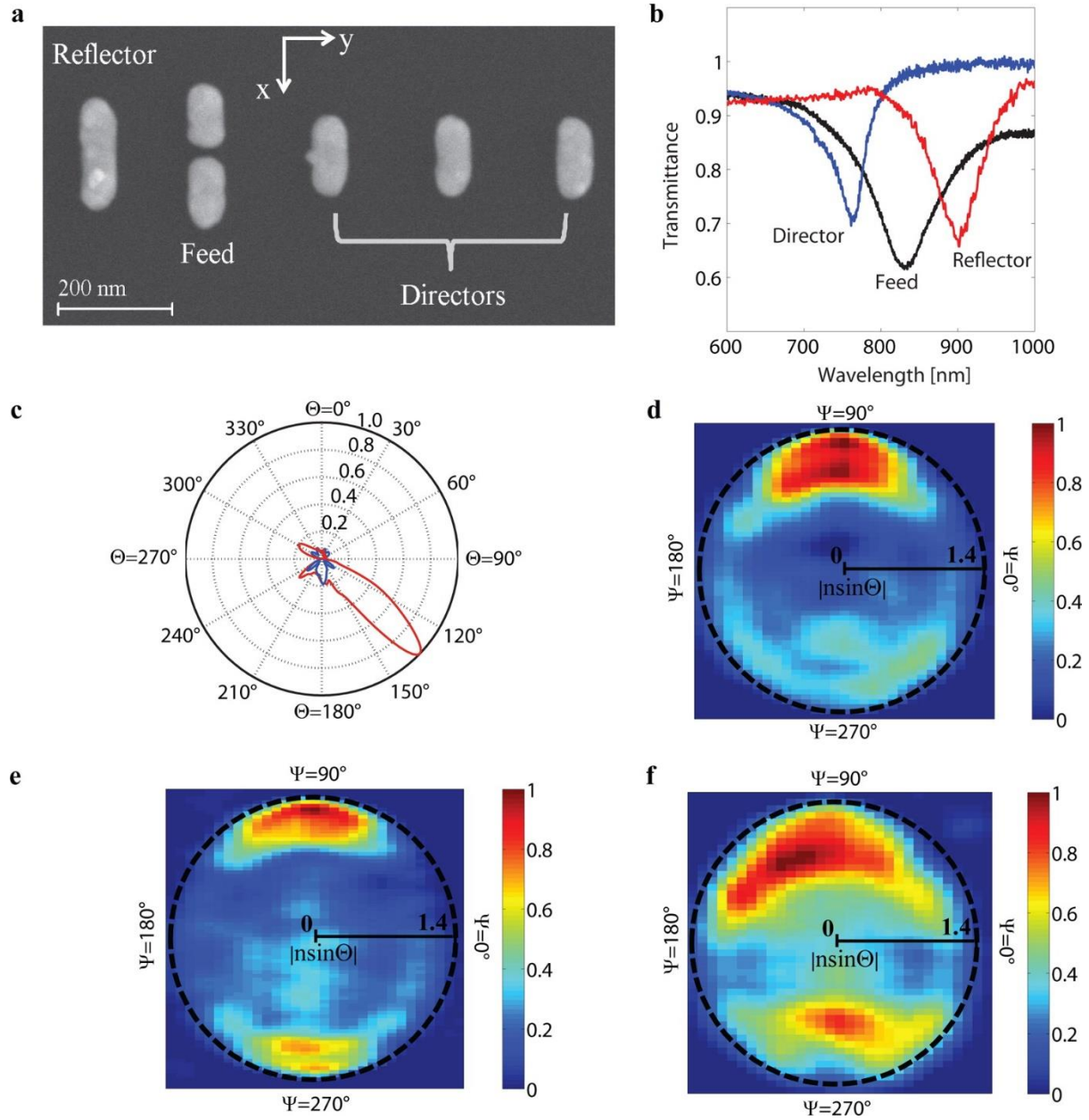


Fig.11: SERS emission pattern from Yagi Uda antenna. (a) SEM image of the Yagi Uda antenna. (b) Measured transmittance of arrays of directors, feed antennas, and reflectors. (c) FEM simulation of emission pattern resulting from an electric dipole (free space wavelength $\lambda = 857$ nm) being placed in center of gap in feed element of YU antenna; red: emission pattern in yz -plane and blue: emission pattern in xz -plane. (d) Emission pattern of thiophenol 1074 cm^{-1} Raman line retrieved from E-k measurements. Measured emission patterns of (e) 415 and (f) 1586 cm^{-1} Raman lines. Reproduced with permission from (28)

In the similar manner molecular vibrational transitions, especially Raman, gets affected when the molecule is coupled to a nanostructure. Wavevector distribution of SERS, in such cases, can alter drastically. BFP imaging and spectroscopy has been utilized to understand SERS with particular emphasis to wavevector distribution, as shown in figure 11(28). This has been further extrapolated to study coupling strengths and molecular orientations, which has applications in quantum electrodynamics.

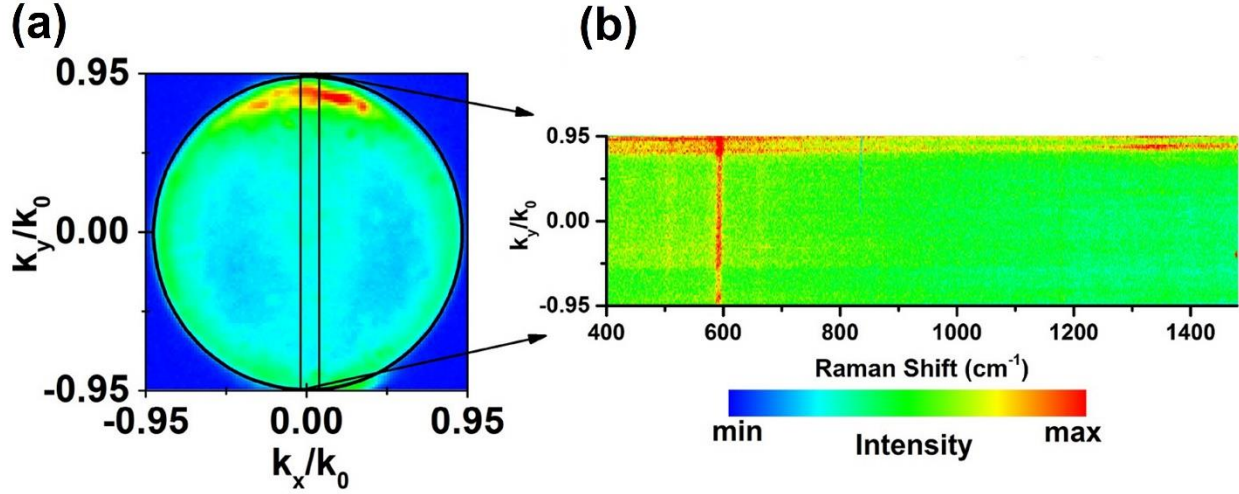


Fig.4: Imaging remote excitation SERS radiation pattern by filtering narrow range of wavevectors from the BFP. (a) Back focal plane image captured from the distal end of a film coupled plasmonic nanowire cavity, showing unidirectional emission pattern. (b) Energy-momentum spectrum captured by filtering a narrow range of the back focal plane image near $k_x/k_0=0$. Reproduced with permission from (13)

Another important aspect of BFP imaging is to understand the coupling of molecules to nanocavities. Nanocavities have specific wavevector signatures which can alter the emission wavevectors of the molecules coupled to such cavities. This in turn provides a unique way to measure the coupling strength and efficiency of molecules. Certain nanocavities, eg. Elongated cavities like film coupled nanowire cavity, have very complex polarization signatures. This makes polarization sensitive emission process like SERS couple differently to the system when compared to molecular fluorescence, as shown by Vasista et.al.(13).

4.5. Nonlinear scattering

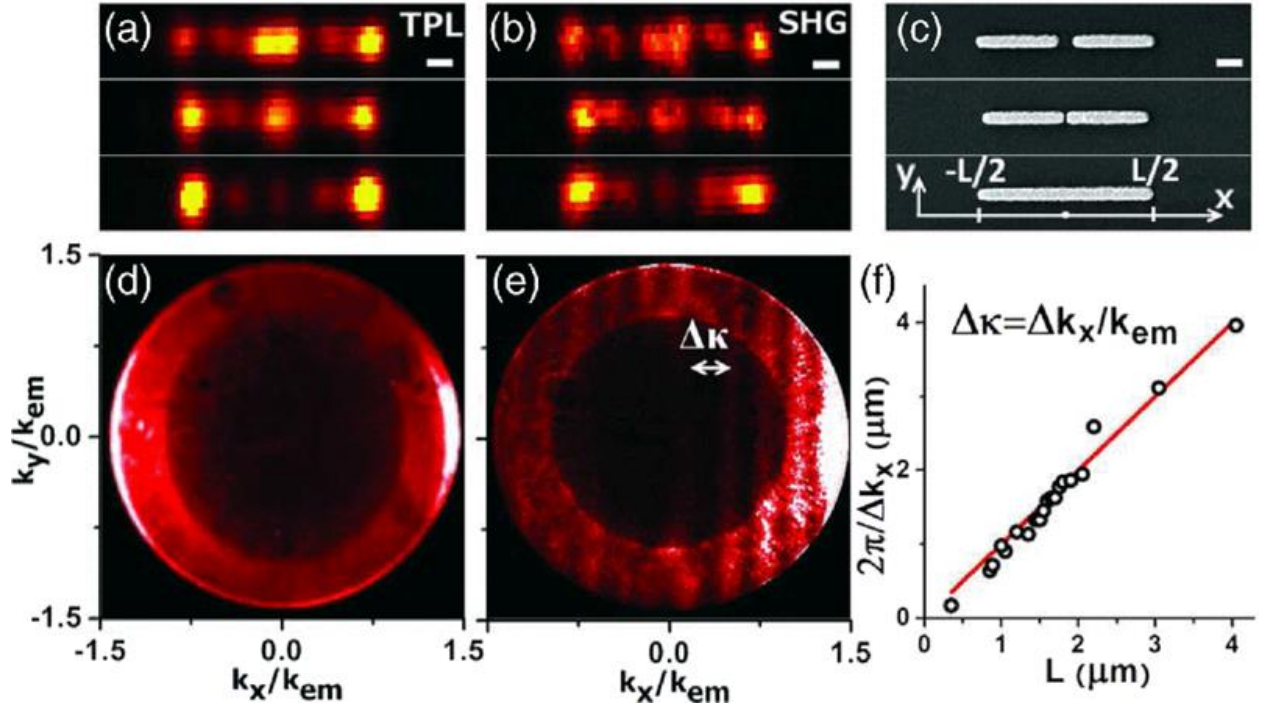


Fig.12: Back focal plane images of Nonlinear signals from Au nanoantenna. a)–(c) TPL and SHG confocal maps and SEM micrographs of three gold rod antennas. (d) and (e) are angular distributions of TPL and SHG signal from a single nanowire [bottom antenna in (a)–(c)]. The laser is focused at the left extremity. (f) Dependence of $2\pi/\Delta k_x$ on L (black \circ) and $2\pi/\Delta k_x=L$ (line). The axes are normalized in units of $k_{em}=nm2\pi/\lambda_{em}$, where λ_{em} is the emission wavelength of the nonlinear process in vacuum and $nm=1.55$ is the refractive index of glass and oil. The interfringe distance Δk is denoted with the white double-headed arrow. Reproduced with permission from (64)

Depending on the generation nonlinear signals can be incoherent like two photon luminescence (TPL) or coherent like second harmonic generation (SHG), third harmonic generation etc. Both coherent and incoherent nonlinear signals have unique spectral properties and hence are distinguishable spectrally. But just observing the real plane images (shown in Fig. 12 (a) and (b)) one cannot tell about the far field emission patterns. Back focal plane images shown in figure 12 (d) for TPL and figure 12 (e) for SHG shows that two signals have different far field patterns. Formation of fringes in back focal plane in figure 12 (e) again confirms the coherent nature of the SHG(64).

5. Conclusions

In summary, we discussed the technique of back focal plane imaging and spectroscopy in detail. We have also pondered upon different applications of the technique with specific emphasis on recent developments. Back focal plane imaging and spectroscopy forms a crucial factor in understanding emission processes like Raman spectroscopy(18), molecular fluorescence (65, 66), bio-imaging(67), non linear scattering(68), to understand strong coupling(69) etc along with elastic scattering by nano/micro

structures(70). This knowledge can be further extrapolated to quantum optics and quantum electrodynamics to understand coupling mechanisms between atom and cavity etc.

References

1. A. G. Curto *et al.*, Multipolar radiation of quantum emitters with nanowire optical antennas. *Nat. Commun.* **4**, 1750 (2013).
2. A. G. Curto *et al.*, Unidirectional emission of a quantum dot coupled to a nanoantenna. *Science* **329**, 930-933 (2010).
3. M. A. Lieb, J. M. Zavislan, L. Novotny, Single-molecule orientations determined by direct emission pattern imaging. *J.Opt.Soc.Am. B* **21**, 1210-1215 (2004).
4. S. Kühn, G. Mori, M. Agio, V. Sandoghdar, Modification of single molecule fluorescence close to a nanostructure: radiation pattern, spontaneous emission and quenching. *Mol. Phys.* **106**, 893-908 (2008).
5. D. Wang, W. Zhu, M. D. Best, J. P. Camden, K. B. Crozier, Directional Raman scattering from single molecules in the feed gaps of optical antennas. *Nano Lett.* **13**, 2194-2198 (2013).
6. I. Sersic, C. Tuambilangana, A. F. Koenderink, Fourier microscopy of single plasmonic scatterers. *New J. Phys.* **13**, 083019 (2011).
7. L. Novotny, N. Van Hulst, Antennas for light. *Nat. Photonics* **5**, 83 (2011).
8. T. Shegai *et al.*, Unidirectional broadband light emission from supported plasmonic nanowires. *Nano Lett.* **11**, 706-711 (2011).
9. Z. Li *et al.*, Directional light emission from propagating surface plasmons of silver nanowires. *Nano Lett.* **9**, 4383-4386 (2009).
10. P. Biagioni, J.-S. Huang, B. Hecht, Nanoantennas for visible and infrared radiation. *Rep. Prog. Phys.* **75**, 024402 (2012).
11. G. Grzela *et al.*, Nanowire antenna emission. *Nano Lett.* **12**, 5481-5486 (2012).
12. G. Grzela *et al.*, Nanowire antenna absorption probed with time-reversed fourier microscopy. *Nano Lett.* **14**, 3227-3234 (2014).
13. A. Vasista *et al.*, Differential Wavevector Distribution of Surface-Enhanced Raman Scattering and Fluorescence in a Film-Coupled Plasmonic Nanowire Cavity. *Nano Lett.* **18**, 650-655 (2018).
14. T. Coenen, E. J. R. Vesseur, A. Polman, A. F. Koenderink, Directional emission from plasmonic Yagi-Uda antennas probed by angle-resolved cathodoluminescence spectroscopy. *Nano Lett.* **11**, 3779-3784 (2011).
15. T. Coenen, E. J. R. Vesseur, A. Polman, Angle-resolved cathodoluminescence spectroscopy. *Appl. Phys. Lett.* **99**, 143103 (2011).
16. J. A. Schuller *et al.*, Orientation of luminescent excitons in layered nanomaterials. *Nat. Nanotechnol.* **8**, 271 (2013).
17. J. J. Baumberg *et al.*, Angle-resolved surface-enhanced Raman scattering on metallic nanostructured plasmonic crystals. *Nano Lett.* **5**, 2262-2267 (2005).
18. H. Budde *et al.*, Raman radiation patterns of graphene. *ACS Nano* **10**, 1756-1763 (2015).
19. B. Fazio *et al.*, Coherent backscattering of Raman light. *Nat. Photonics* **11**, 170-176 (2017).
20. T. Shegai, B. r. Brian, V. D. Miljkovic, M. Kall, Angular distribution of surface-enhanced Raman scattering from individual Au nanoparticle aggregates. *ACS Nano* **5**, 2036-2041 (2011).
21. D. Piatkowski *et al.*, Silver nanowires as receiving-radiating nanoantennas in plasmon-enhanced up-conversion processes. *Nanoscale* **7**, 1479-1484 (2015).
22. S. Shatalin, R. Juškaitis, J. Tan, T. Wilson, Reflection conoscopy and micro-ellipsometry of isotropic thin film structures. *J. Microsc.* **179**, 241-252 (1995).

23. F. E. Wright, Transmission of light through transparent inactive crystal plates, with special reference to observations in convergent polarized light. *Am. J. Sci. Series 4 Vol. 31*, 157-211 (1911).
24. H. Kubota, S. Inoué, Diffraction images in the polarizing microscope. *J. Opt. Soc. Am.* **49**, 191-198 (1959).
25. P. Muehlschlegel, H.-J. Eisler, O. J. Martin, B. Hecht, D. Pohl, Resonant optical antennas. *Science* **308**, 1607-1609 (2005).
26. E. Betzig, J. Trautman, T. Harris, J. Weiner, R. Kostelak, Breaking the diffraction barrier: optical microscopy on a nanometric scale. *Science* **251**, 1468-1470 (1991).
27. M. Villiger, C. Pache, T. Lasser, Dark-field optical coherence microscopy. *Opt. Lett.* **35**, 3489-3491 (2010).
28. W. Zhu, D. Wang, K. B. Crozier, Direct observation of beamed Raman scattering. *Nano Lett.* **12**, 6235-6243 (2012).
29. Y. N. Polivanov, Raman scattering of light by polaritons. *Phys. Usp.* **21**, 805-831 (1978).
30. R. Wagner, L. Heerklotz, N. Kortenbruck, F. Cichos, Back focal plane imaging spectroscopy of photonic crystals. *Appl. Phys. Lett.* **101**, 081904 (2012).
31. A. Brzezinski *et al.*, Enhanced emission from fcc fluorescent photonic crystals. *Phys. Rev. B* **77**, 233106 (2008).
32. B. E. Saleh, M. C. Teich, B. E. Saleh, *Fundamentals of photonics*. (Wiley New York, 1991), vol. 22.
33. C. Regan, A. Krishnan, R. Lopez-Boada, L. Grave de Peralta, A. Bernussi, Direct observation of photonic Fermi surfaces by plasmon tomography. *Appl. Phys. Lett.* **98**, 151113 (2011).
34. L. Langguth, D. Punj, J. Wenger, A. F. Koenderink, Plasmonic band structure controls single-molecule fluorescence. *ACS Nano* **7**, 8840-8848 (2013).
35. C. I. Osorio, A. Mohtashami, A. F. Koenderink, K-space polarimetry of bullseye plasmon antennas. *Sci. Rep.* **5**, 9966 (2015).
36. V. D. Miljković, T. Shegai, P. Johansson, M. Käll, Simulating light scattering from supported plasmonic nanowires. *Opt. Express* **20**, 10816-10826 (2012).
37. A. Powell, N. Hjerrild, A. Watt, H. Assender, J. Smith, Directional plasmonic scattering from metal nanoparticles in thin-film environments. *Appl. Phys. Lett.* **104**, 081110 (2014).
38. A. B. Vasista, R. P. Tripathi, S. K. Chaubey, S. Tiwari, G. P. Kumar, Doughnut-shaped emission from vertical organic nanowire coupled to thin plasmonic film. *Opt. Lett.* **43**, 923-926 (2018).
39. D. K. Sharma, V. Kumar, A. B. Vasista, S. K. Chaubey, G. P. Kumar, Spin-Hall effect in the scattering of structured light from plasmonic nanowire. *Opt. Lett.* **43**, 2474-2477 (2018).
40. D. G. Zhang, X. Yuan, A. Bouhelier, Direct image of surface-plasmon-coupled emission by leakage radiation microscopy. *Appl. Opt.* **49**, 875-879 (2010).
41. N. Hartmann, D. Piatkowski, R. Ciesielski, S. Mackowski, A. Hartschuh, Radiation channels close to a plasmonic nanowire visualized by back focal plane imaging. *ACS Nano* **7**, 10257-10262 (2013).
42. A. Dasgupta, D. Singh, R. P. Tripathi, G. P. Kumar, Directional Fluorescence Emission Mediated by Chemically-Prepared Plasmonic Nanowire Junctions. *J. Phys. Chem. C* **120**, 17692-17698 (2016).
43. A. Drezet, C. Genet, Imaging surface plasmons: from leaky waves to far-field radiation. *Phys. Rev. Lett.* **110**, 213901 (2013).
44. J. Grandidier *et al.*, Dielectric-loaded surface plasmon polariton waveguides: figures of merit and mode characterization by image and Fourier plane leakage microscopy. *Phys. Rev. B* **78**, 245419 (2008).
45. V. Garbin *et al.*, Mie scattering distinguishes the topological charge of an optical vortex: a homage to Gustav Mie. *New J. Phys.* **11**, 013046 (2009).
46. G. Lozano, T. Barten, G. Grzela, J. G. Rivas, Directional absorption by phased arrays of plasmonic nanoantennae probed with time-reversed Fourier microscopy. *New J. Phys.* **16**, 013040 (2014).

47. C. M. Dodson, J. A. Kurvits, D. Li, R. Zia, Wide-angle energy-momentum spectroscopy. *Opt. Lett.* **39**, 3927-3930 (2014).
48. T. H. Taminiau, S. Karaveli, N. F. Van Hulst, R. Zia, Quantifying the magnetic nature of light emission. *Nat. Commun.* **3**, 979 (2012).
49. Y. Chu, W. Zhu, D. Wang, K. B. Crozier, Beamed Raman: directional excitation and emission enhancement in a plasmonic crystal double resonance SERS substrate. *Opt. Express* **19**, 20054-20068 (2011).
50. D. Wang, W. Zhu, Y. Chu, K. B. Crozier, High directivity optical antenna substrates for surface enhanced Raman scattering. *Adv. Mater.* **24**, 4376-4380 (2012).
51. J. A. Kurvits, M. Jiang, R. Zia, Comparative analysis of imaging configurations and objectives for Fourier microscopy. *J. Opt. Soc. Am. A* **32**, 2082-2092 (2015).
52. N. Coca-López *et al.*, Remote excitation and detection of surface-enhanced Raman scattering from graphene. *Nanoscale*, (2018).
53. D. Singh, A. Dasgupta, V. Aswathy, R. P. Tripathi, G. P. Kumar, Directional out-coupling of light from a plasmonic nanowire-nanoparticle junction. *Opt. Lett.* **40**, 1006-1009 (2015).
54. M. Song *et al.*, Selective excitation of surface plasmon modes propagating in Ag nanowires. *Opt. Express* **25**, 9138-9149 (2017).
55. K. Hassan *et al.*, Momentum-space spectroscopy for advanced analysis of dielectric-loaded surface plasmon polariton coupled and bent waveguides. *Phys. Rev. B* **87**, 195428 (2013).
56. C. Huang *et al.*, Gain, detuning, and radiation patterns of nanoparticle optical antennas. *Phys. Rev. B* **78**, 155407 (2008).
57. J. Berthelot *et al.*, Determinant role of the edges in defining surface plasmon propagation in stripe waveguides and tapered concentrators. *J. Opt. Soc. Am. B* **29**, 226-231 (2012).
58. D. Vercruysse *et al.*, Unidirectional side scattering of light by a single-element nanoantenna. *Nano Lett.* **13**, 3843-3849 (2013).
59. H. Yang, M. Qiu, Q. Li, Identification and control of multiple leaky plasmon modes in silver nanowires. *Laser Photonics Rev.* **10**, 278-286 (2016).
60. R. P. Tripathi, A. B. Vasista, R. Chikkaraddy, G. Kumar, Radiative Channeling of Nanowire Frenkel Exciton Polaritons through Surface Plasmons. *Adv. Opt. Mater.* **5**, (2017).
61. T. Taminiau, F. Stefani, F. B. Segerink, N. Van Hulst, Optical antennas direct single-molecule emission. *Nat. Photonics* **2**, 234 (2008).
62. H. Aouani *et al.*, Plasmonic antennas for directional sorting of fluorescence emission. *Nano Lett.* **11**, 2400-2406 (2011).
63. S.-H. Cao, W.-P. Cai, Q. Liu, Y.-Q. Li, Surface plasmon-coupled emission: what can directional fluorescence bring to the analytical sciences? *Annu. Rev. Anal. Chem.* **5**, 317-336 (2012).
64. S. Viarbitskaya, O. Demichel, B. Cluzel, G. C. des Francs, A. Bouhelier, Delocalization of nonlinear optical responses in plasmonic nanoantennas. *Phys. Rev. Lett.* **115**, 197401 (2015).
65. D. Zhang *et al.*, Back focal plane imaging of directional emission from dye molecules coupled to one-dimensional photonic crystals. *Nanotechnology* **25**, 145202 (2014).
66. A. L. Mattheyses, D. Axelrod, Fluorescence emission patterns near glass and metal-coated surfaces investigated with back focal plane imaging. *J. Biomed. Opt.* **10**, 054007 (2005).
67. M. W. Allersma, F. Gittes, R. J. Stewart, C. F. Schmidt, Two-dimensional tracking of ncd motility by back focal plane interferometry. *Biophys. J* **74**, 1074-1085 (1998).
68. D. Wolf, T. Schumacher, M. Lippitz, Shaping the nonlinear near field. *Nat. Commun.* **7**, 10361 (2016).
69. G. Beane, B. S. Brown, P. Johns, T. Devkota, G. V. Hartland, Strong Exciton-Plasmon Coupling in Silver Nanowire Nanocavities. *J. Phys. Chem. Lett.* **9**, 1676-1681 (2018).

70. P. Johns, G. Beane, K. Yu, G. V. Hartland, Dynamics of surface plasmon polaritons in metal nanowires. *J. Phys. Chem. C* **121**, 5445-5459 (2017).

Suggested further reading

Saleh, Bahaa EA, Malvin Carl Teich, and Bahaa E. Saleh. *Fundamentals of photonics*. Vol. 22. New York: Wiley, 1991.

Jenkins, Francis A., and Harvey E. White. *Fundamentals of optics*. Tata McGraw-Hill Education, 1937.

Goodman, Joseph W. *Introduction to Fourier optics*. Roberts and Company Publishers, 2005.

Novotny, Lukas, and Niek Van Hulst. "Antennas for light." *Nature photonics* 5, no. 2 (2011): 83.

Novotny, Lukas, and Bert Hecht. *Principles of nano-optics*. Cambridge university press, 2012.

Hecht, Eugene. *Optics*. Pearson Education, 2016.