

Sample Motion for Structured Illumination Fluorescence Microscopy

RUIMING CAO^{1,4,*}, GUANGHAN MENG^{2,3}, AND LAURA WALLER²

¹Department of Bioengineering, University of California, Berkeley, Berkeley CA, 94720

²Department of Electrical Engineering and Computer Sciences, University of California, Berkeley, Berkeley CA, 94720

³School of Optometry and Vision Science, University of California, Berkeley, Berkeley CA, 94720

⁴Chan Zuckerberg Imaging Institute, Redwood City CA, 94605

*rcao@berkeley.edu

Compiled February 7, 2025

Structured illumination microscopy (SIM) uses a set of images captured with different patterned illumination to computationally reconstruct resolution beyond the diffraction limit. Here, we propose an alternative approach using a static speckle illumination pattern and relying on inherent sample motion to encode the super-resolved information in multiple raw images, for the case of fluorescence microscopy. From a set of sequentially captured raw images, we jointly estimate the sample motion and the super-resolved image. We demonstrate the feasibility of the proposed method both in simulation and in experiment.

1. INTRODUCTION

Structured Illumination Microscopy (SIM) recovers a super-resolved image from a set of raw images captured with different illumination patterns (either sinusoidal [1] or speckle [2, 3]). The structured illumination shifts the sample's higher-frequency content to lower spatial frequencies that are within the diffraction-limited bandwidth of the optical system. By capturing a set of raw images with different illumination patterns and solving a computational inverse problem, SIM can reconstruct up to $2\times$ super-resolution. The varying illumination patterns serve as an encoding mechanism to ensure that diverse information covering all spatial frequencies is captured in the raw images. For sinusoidal SIM, ~ 9 raw images are required, with sinusoidal illumination patterns having different orientations and phase shifts. For speckle SIM, dozens of raw images are needed for a well-posed inverse problem.

Updating the illumination pattern dynamically for each exposure causes SIM systems to be complicated and expensive. Furthermore, acquisition times can be long, especially if the optomechanical components are not as fast as the sensor framerate [4–6]. This is particularly limiting for imaging live samples that moves in-between different raw images, causing the reconstruction to suffer from motion artifacts. Our previous study developed a computational framework to estimate sample motion and eliminate motion artifacts [7].

In this paper, we explore an alternative approach, called Speckle Flow SIM, that uses a single structured illumination pattern during the acquisition of a sequence of raw images (Fig. 1). Instead of varying the illumination over time, Speckle Flow SIM relies on sample motion to encode diverse spatial frequency information at different exposure times. Consider, for example,

the case of a sample moving with translational motion. As the sample moves across the field-of-view, each part of the sample sees a different speckle illumination pattern at each timepoint. By collecting a series of images as the sample moves, the super-resolution reconstruction problem becomes equivalent to that of having a laterally-shifting speckle illumination, which is straightforward to solve [3]. However, in our case the motion may be unknown and arbitrary.

In order to accommodate unknown motion of the sample, we solve an inverse problem simultaneously for a sequence of super-resolved images, one for each timepoint in the captured raw data. This inverse problem is under-constrained because the spatial bandwidth of each raw image covers only a part of the bandwidth needed for the super-resolved reconstruction. Nevertheless, it can be solved if the motion is smooth or low-rank in some way. To do so, we turn to a neural space-time model (NSTM) [8] for modeling the spatiotemporal relationship of the moving sample. The NSTM assumes that the super-resolved image at one timepoint can be obtained by applying a motion transformation to the super-resolved image at a different timepoint. In essence, this means that result can be represented by a single super-resolved reference image along with motion maps corresponding to each timepoint. The NSTM then simultaneously estimates the sample motion and the super-resolved image. While our Speckle Flow SIM idea has been previously demonstrated for coherent imaging [8], this paper focuses on its application in fluorescence imaging.

We choose to use a speckle illumination pattern because it is non-directional and therefore suitable for achieving isotropic resolution gain. Sinusoidal illumination, in contrast, can only modulate high-frequency information along one orientation into

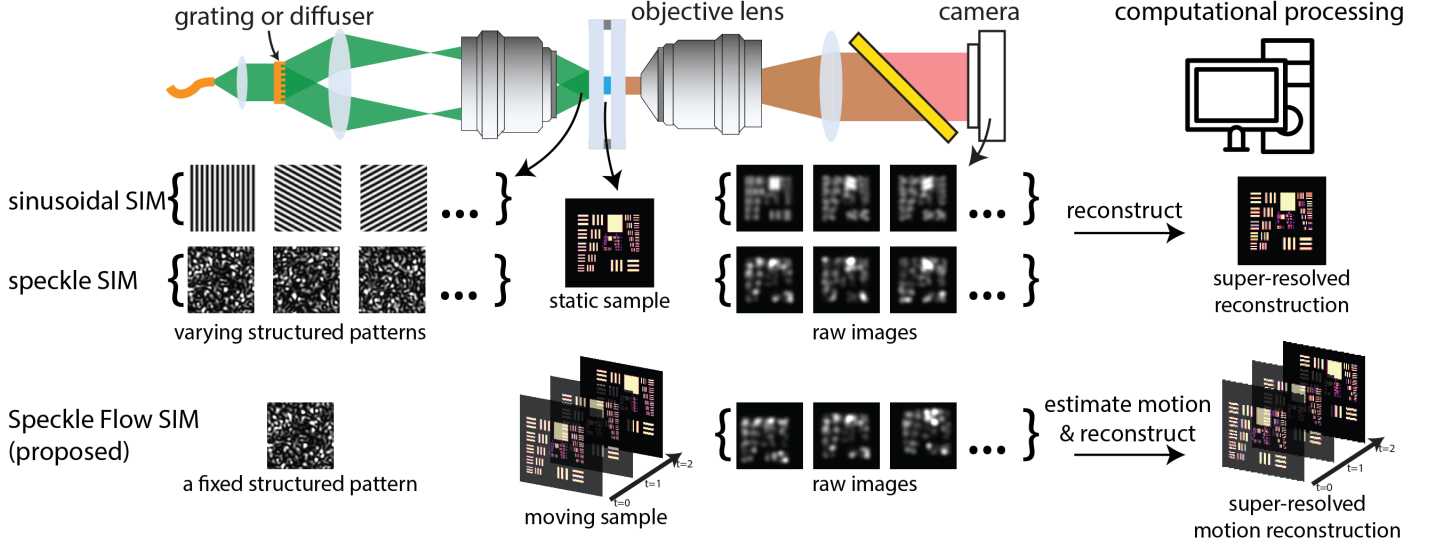


Fig. 1. Existing SIM methods (e.g. sinusoidal SIM, speckle SIM) vary the structured excitation patterns dynamically to encode super-resolution information into multiple raw images. We propose to implement SIM using a fixed speckle illumination pattern and a moving sample, such that each part of the sample is illuminated by multiple different speckle patterns over time. Our algorithm then uses the captured images to solve for both the object’s motion and the super-resolved reconstruction.

the system bandwidth. Speckle patterns are also simple to generate by propagating a coherent beam through a random phase element, such as a diffuser. Unlike previous speckle SIM approaches that either rely on statistical properties [2] or require estimating the pattern during reconstruction [3, 9], our method’s use of a static speckle pattern allows for pre-calibrating the actual pattern before the imaging experiment. Therefore, the reconstruction algorithm for Speckle Flow SIM can operate with known illumination, requiring significantly fewer raw images (~ 10) than previous speckle SIM (~ 50).

2. METHODS

A. SIM forward model

SIM methods use patterned illumination, I , to excite a fluorescent sample, S . The emitted light is convolved with the diffraction-limited system’s point-spread function (PSF) before being captured by the camera. The raw data, D , is expressed as

$$\begin{aligned} D(\mathbf{r}) &= (I(\mathbf{r}) \cdot S(\mathbf{r})) * \text{PSF}(\mathbf{r}) \\ &= \mathcal{F}^{-1} [(\tilde{I}(\mathbf{u}) * \tilde{S}(\mathbf{u})) \cdot \text{OTF}(\mathbf{u})], \end{aligned} \quad (1)$$

where $*$ denotes convolution, \cdot denotes Hadamard product (element-wise matrix product), \mathcal{F} , $\tilde{\cdot}$ denote 2D Fourier transform, and \mathbf{r} and \mathbf{u} denote 2D coordinates in real and frequency space, respectively. The optical transfer function (OTF) is the Fourier transform of the PSF; it acts as a low-pass filter here that zeros out high-frequency information beyond the diffraction limit.

If the illumination pattern contains any non-zero frequency content, then object frequencies beyond the diffraction-limited bandwidth will be encoded into D due to the convolution between \tilde{I} and \tilde{S} . Both sinusoidal and speckle patterns meet this condition; for example, a sinusoidal illumination with a spatial frequency corresponding to the edge of the system passband results in capture of spatial frequencies from the sample up to $2 \times$ the diffraction limit. Using a sinusoidal pattern in a single orientation can only modulate frequencies along that orientation,

so multiple orientations (and phase shifts) are needed to fully encode the super-resolution content. Speckle Flow SIM achieves isotropic modulation by using a speckle pattern, which provides more uniform coverage across all orientations. In both cases, the super-resolution information can be recovered computationally from multiple raw images, as illustrated in Fig. 1. Speckle Flow SIM pre-calibrates the static speckle pattern, so that its reconstruction is more data-efficient compared with speckle SIM methods whose patterns are not fully known.

B. Sample motion enabled SIM

Speckle SIM encodes a broad set of spatial frequency information from the sample in each raw image, but still requires multiple raw images in order to disentangle spatial frequencies from the convolution in Eq. 1. In our setup with a static speckle illumination, we collect a sequence of raw images as the sample moves; thus, additional super-resolution information is encoded into the raw data at different timepoints. For example, in the simple case of a sample undergoing known translational motion with a relative displacement $\delta\mathbf{r}$, the raw image can be expressed as:

$$\begin{aligned} D(\mathbf{r}) &= (I(\mathbf{r}) \cdot S(\mathbf{r} + \delta\mathbf{r})) * \text{PSF}(\mathbf{r}) \\ &= \mathcal{F}^{-1} [(\tilde{I}(\mathbf{u}) * (\tilde{S}(\mathbf{u}) e^{-2\pi i \delta\mathbf{r} \cdot \mathbf{u}})) \cdot \text{OTF}(\mathbf{u})]. \end{aligned} \quad (2)$$

The motion term, $e^{-2\pi i \delta\mathbf{r} \cdot \mathbf{u}}$, introduces additional frequency components as the displacement $\delta\mathbf{r}$ changes, enriching the encoded super-resolution information. As more raw images are acquired at varying sample displacements, the super-resolved reconstruction becomes well-conditioned.

C. Joint reconstruction of motion and sample

Our neural space-time model (NSTM) [7] has been shown capable of estimating unknown arbitrary deformable motion, as long as it is relatively smooth. NSTM models a dynamic scene using two multi-layer perceptrons: a motion network, $f(\cdot, \theta_{\text{motion}})$, and a scene network, $f(\cdot, \theta_{\text{scene}})$. At each

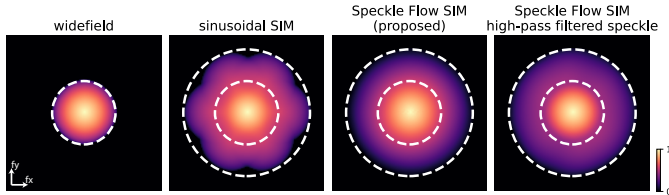


Fig. 2. Comparison of the transfer functions for widefield, sinusoidal SIM and Speckle Flow SIM (with and without high-pass filtering). The inner dashed white circle indicates the system’s diffraction-limited bandwidth, and the outer dash white circle indicates the theoretical $2\times$ super-resolution bandwidth for SIM. Gamma corrections of 0.3 are applied to all OTFs for visibility.

timepoint t , the moving sample is represented by $S(\mathbf{r}, t) = \text{NSTM}(\mathbf{r}, t; \theta_{\text{motion}}, \theta_{\text{scene}})$. With a set of raw images, D_{t_1}, \dots, D_{t_n} taken at t_1, \dots, t_n , the reconstruction becomes

$$\min_{\theta_{\text{motion}}, \theta_{\text{scene}}} \sum_t \|D_t(\mathbf{r}, t) - (I(\mathbf{r}) \cdot \text{NSTM}(\mathbf{r}, t; \theta_{\text{motion}}, \theta_{\text{scene}})) * \text{PSF}(\mathbf{r})\|^2. \quad (3)$$

We use gradient descent to jointly optimize the weights of the scene and motion networks.

3. SIMULATION RESULTS

A. Encoding information with motion

To demonstrate the super-resolving capability of Speckle Flow SIM, we present two simulation studies: first, assessing the ability to encode super-resolution information during the forward process, and second, evaluating the performance of the inverse reconstruction process in decoding super-resolved images. To assess encoding capabilities, we compare the forward process transfer function of our method with other SIM methods. As SIM methods require more than one raw image, the transfer function for each method measures how much signal can be obtained for each possible input frequency as an impulse. Mathematically, the transfer function, $TF(f_x, f_y) = \sum_t \|D_t(\delta(f_x, f_y)) - D_t(\mathbf{0})\|$, where δ is an impulse at a particular spatial frequency and $D_t(\mathbf{0})$ is the background raw image. Figure 2 shows the simulated OTF for traditional widefield microscopy, sinusoidal SIM, and Speckle Flow SIM. Sinusoidal SIM uses nine raw images: three orientations and three phase shifts of the sinusoid per orientation. Speckle Flow SIM also uses nine raw images with the sample shifting laterally between images. Figure 2 demonstrates that Speckle Flow SIM indeed encodes information beyond the diffraction limit, though its encoding power at the highest spatial frequencies is weaker than sinusoidal SIM because the power of a speckle pattern is more concentrated in the lower frequency region compared with a sinusoidal pattern. This limitation may be mitigated by applying high-pass filtering to the speckle pattern, where convenient (Fig. 2).

B. Decoding analysis

Next, we compared the reconstruction capabilities of Speckle Flow SIM with sinusoidal SIM. Both methods used the same gradient descent algorithm and SIM forward model for super-resolution reconstruction from raw images, assuming the sample

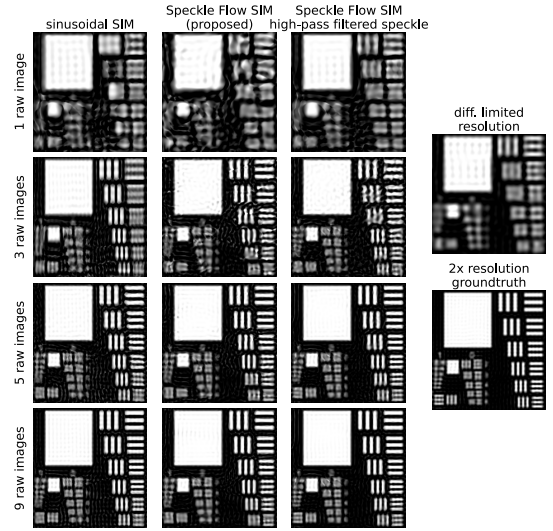


Fig. 3. Reconstruction of a USAF-1951 resolution target using different SIM methods in simulation. Each row compares the reconstruction with a different number of captured raw images. The sample is shifted laterally between raw images. All reconstructions here are obtained using gradient descent optimization with identical parameters and known motion parameters. For sinusoidal SIM, the 3 raw images are taken at the first sinusoidal orientation with three phase shifts, and the 5 raw images contains the 3 images at the first orientation and 2 images at the second orientation.

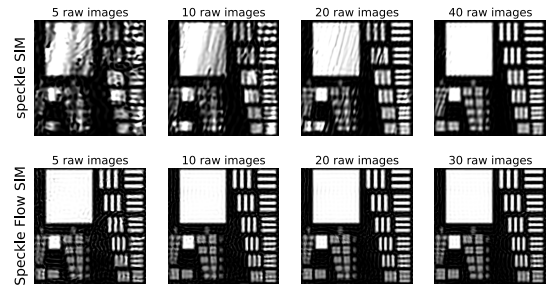


Fig. 4. Reconstruction comparison between speckle SIM and Speckle Flow SIM using different number of raw images in simulation.

motion is known or has been estimated accurately. Figure 3 shows the reconstruction with different numbers of raw images for each method. Due to the directional nature of sinusoidal patterns, the reconstruction quality of sinusoidal SIM suffers from anisotropic resolution when fewer raw images are captured. As shown in Fig. 3, with three raw images, sinusoidal SIM fails to resolve any high-frequency vertical features as the sinusoidal patterns only modulate in the horizontal direction. On the other hand, speckle patterns have isotropic frequency coverage, and thus are robust with a reduced number of raw images.

Similarly, assuming known motion, we compare the data efficiency between speckle SIM, which jointly resolves for the speckle pattern during the reconstruction, and Speckle Flow SIM that knows the exact speckle pattern from the pre-calibration. As shown in the simulation in Fig. 4, Speckle Flow SIM requires much less raw images than speckle SIM, to obtain a high quality reconstruction, because of the a priori knowledge of the speckle. Speckle Flow SIM using 10 raw images achieves better recon-

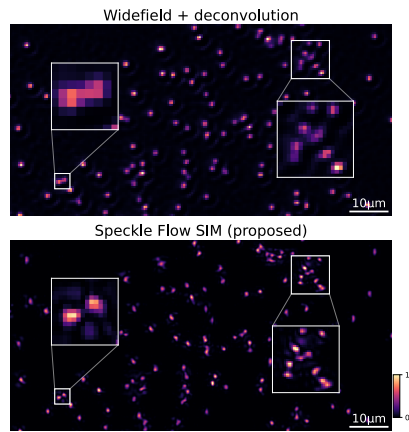


Fig. 5. Experimental result of Speckle Flow SIM on a fluorescence-labeled microbead sample. The left shows the image with widefield excitation and image deconvolution using the theoretical system PSF. The right shows the reconstruction of Speckle Flow SIM. During the acquisition, the sample moves continuously on a motorized lateral stage. The full video is in Supplement Video 1.

struction than speckle SIM using 40 images.

4. EXPERIMENTAL RESULTS

We built a transmission mode setup to experimentally demonstrate the concept of Speckle Flow SIM. A 488nm laser (Coherent OBIS LS) serves as the light source. As in Fig. 1, the beam passes through a beam expander (Thorlabs GBE05-A), and then the speckle pattern is generated by light interference after two layers of Scotch tape. This speckle pattern excites the sample and the emitted light is magnified by a 4- f system consisting of a $10\times$ 0.25NA objective (Nikon Plan) and a tube lens (250mm focal length). The excitation light is filtered by a dichroic mirror and an emission band-pass filter (Edmund Optics 67-079 and 86-984), and a sCMOS (PCO Edge 5.5) collects the fluorescent signal. To calibrate the system, the excitation pattern is pre-captured by the same camera without the sample and the emission band-pass filter in place.

A $0.71\mu\text{m}$ -diameter fluorescence-tagged microbead sample (Thermo Scientific) is placed on a motorized lateral stage (Thorlabs PT1-Z8). During acquisition, a sequence of raw images is captured as the lateral stage continuously moves the sample. Although the sample motion is controlled by the stage, we assume that the stage motion is unknown to the reconstruction algorithm. NSTM is used to jointly solve for the motion and the scene by optimizing the objective defined in Eq. 3. We run gradient descent optimization for 10,000 epochs with an initial learning rate of 0.0001 for both the motion and scene networks, and the learning rate has a decay of 0.8 every 1,000 epochs.

Figure 5 presents the reconstruction result from the Speckle Flow SIM experimental system and compares it with the widefield image using PSF deconvolution. Speckle Flow SIM successfully resolves individual microbeads that the widefield image fails to capture, demonstrating that the contrast introduced by sample motion enables super-resolution imaging. With the joint motion and scene reconstruction from NSTM, Speckle Flow SIM recovers the dynamic super-resolved images for a moving sample with unknown motion. The full dynamic super-resolution reconstruction is in Supplemented Video 1.

5. DISCUSSION

Sample motion between captured images has traditionally been viewed as a problem for super-resolution microscopy, causing motion blurring artifacts. Here, we have utilized it as a novel contrast mechanism for SIM. Rather than relying on time-varying excitation patterns, we rely on sample motion to encode diverse information into the measurements even with a fixed excitation pattern. This greatly simplifies the implementation of SIM systems. With a fixed pattern, optomechanical components or spatial light modulators are no longer needed to update the excitation patterns, consequently reducing the equipment cost as well as eliminating the time delay for the pattern update. Furthermore, a fixed speckle pattern is also easier to pre-calibrate compared to previous speckle SIM methods [2], which enhances the data efficiency of the reconstruction.

A limitation of our work is that the motion contrast mechanism may not encode sufficient information when some parts of a sample stay static. Then, the raw data will encode identical information over multiple exposures, resulting in an ill-posed inverse problem. This problem can be resolved by introducing external global motion to the sample, like shifting the sample laterally as in Fig. 5.

The proposed motion contrast mechanism is enabled by recent advances in dynamic imaging reconstruction, allowing us to jointly estimate both the motion and the scene. This joint estimation relies on a temporal redundancy assumption, which means that the scene at any given timepoint is highly correlated with the scene at the previous timepoint. However, we want to note that this assumption does not always hold true as discussed in [7], potentially causing the failure of motion estimation. This limitation restricts the use of Speckle Flow SIM to samples exhibiting smooth motion, such as cells in microfluidic channels, rather than samples with abrupt, on/off motion like neuron firing. In future research, other types of motion assumptions and modelings can be explored to alleviate this motion constraint.

In addition, fast sample motion may also blur individual raw images. Since NSTM assumes the sample remains motionless during each exposure, any movement within a single exposure cannot be accounted for. This can lead to information loss during encoding and cause model mismatch in motion estimation. Therefore, it is essential to carefully select the exposure settings to minimize blur in the raw images, which is a common requirement in other microscopy methods.

In summary, we introduced Speckle Flow SIM, which employs a static speckle illumination pattern to encode super-resolution information in fluorescence imaging. We demonstrated its feasibility through simulation analysis and proof-of-concept experiments. In our experiments, Speckle Flow SIM achieved super-resolution reconstruction on a continuously moving fluorescently-labeled sample. This motion contrast mechanism creates new opportunities for innovative optical designs for SIM systems.

REFERENCES

1. M. G. Gustafsson, "Surpassing the lateral resolution limit by a factor of two using structured illumination microscopy," *J. microscopy* **198**, 82–87 (2000).
2. E. Mudry, K. Belkebir, J. Girard, *et al.*, "Structured illumination microscopy using unknown speckle patterns," *Nat. Photonics* **6**, 312–315 (2012).
3. L.-H. Yeh, S. Chowdhury, and L. Waller, "Computational structured illu-

- mination for high-content fluorescence and phase microscopy,” *Biomed. optics express* **10**, 1978–1998 (2019).
4. R. Fiolka, L. Shao, E. H. Rego, *et al.*, “Time-lapse two-color 3d imaging of live cells with doubled resolution using structured illumination,” *Proc. National Acad. Sci.* **109**, 5311–5315 (2012).
 5. R. Förster, H.-W. Lu-Walther, A. Jost, *et al.*, “Simple structured illumination microscope setup with high acquisition speed by using a spatial light modulator,” *Opt. Express* **22**, 20663–20677 (2014).
 6. D. Dan, M. Lei, B. Yao, *et al.*, “Dmd-based led-illumination super-resolution and optical sectioning microscopy,” *Sci. reports* **3**, 1116 (2013).
 7. R. Cao, N. S. Divekar, J. K. Nuñez, *et al.*, “Neural space–time model for dynamic multi-shot imaging,” *Nat. Methods* pp. 1–6 (2024).
 8. R. Cao, F. L. Liu, L.-H. Yeh, and L. Waller, “Dynamic structured illumination microscopy with a neural space-time model,” in *2022 IEEE International Conference on Computational Photography (ICCP)*, (IEEE, 2022), pp. 1–12.
 9. J. Min, J. Jang, D. Keum, *et al.*, “Fluorescent microscopy beyond diffraction limits using speckle illumination and joint support recovery,” *Sci. reports* **3**, 2075 (2013).