

Fiducial Inference for Random-Effects Calibration Models: Advancing Reliable Quantification in Environmental Analytical Chemistry

Soumya Sahu

Division of Epidemiology and Biostatistics, University of Illinois Chicago

Thomas Mathew

Department of Mathematics and Statistics, University of Maryland Baltimore County

Robert Gibbons

Department of Public Health Sciences, The University of Chicago

Dulal K. Bhaumik

Division of Epidemiology and Biostatistics, University of Illinois Chicago

Department of Psychiatry, University of Illinois Chicago

Abstract

This article addresses calibration challenges in analytical chemistry by employing a random-effects calibration curve model and its generalizations to capture variability in analyte concentrations. The model is motivated by specific issues in analytical chemistry, where measurement errors remain constant at low concentrations but increase proportionally as concentrations rise. To account for this, the model permits the parameters of the calibration curve, which relate instrument responses to true concentrations, to vary across different laboratories, thereby reflecting real-world variability in measurement processes. Traditional large-sample interval estimation methods are inadequate for small samples, leading to the use of an alternative approach, namely the fiducial approach. The calibration curve that accurately captures the heteroscedastic nature of the data, results in more reliable estimates across diverse laboratory conditions. It turns out that the fiducial approach, when used to construct a confidence interval for an unknown concentration, produces a slightly wider width while achieving the desired coverage probability. Applications considered include the determination of the presence of an analyte and the interval estimation of an unknown true analyte concentration. The proposed method is demonstrated for both simulated and real interlaboratory data, including examples involving copper and cadmium in distilled water.

Keywords: Interlaboratory analysis, Fiducial quantity, Hierarchical designs, Confidence width, Rocke-Lorenzato model.

1 Introduction

A critically important problem in interlaboratory calibration is the accurate interval estimation of the unknown true concentration of an analyte in the context of quantitative analysis in analytical chemistry. The problem is especially challenging when dealing with heterocedastic measurement errors at different concentration levels of the analyte. In practice, analytical measurements are routinely used for various purposes, often without considering the associated level of uncertainty. For instance, in environmental monitoring programs, the concentration of low-level environmental pollutants in a single sample is directly compared with regulatory standards without accounting for the uncertainty in those measurements (see Gibbons, 1994, for reviews of statistical issues in environmental monitoring applications). The uncertainty can be quantified by providing confidence bounds for the true concentration level, given measured concentrations obtained under a hierarchical structure, such as concentrations measured across different laboratories. It is common in environmental monitoring to use multiple laboratories for routine monitoring or to cross-validate an anomalous result. This nested structure, where the laboratory is considered a random factor, complicates the computation of confidence bounds for calibration problem.

The practical importance of the problem is two-fold. First, in the standard situation, a particular laboratory provides an estimate of an unknown concentration, in rare cases, however, the laboratory may compute a corresponding confidence interval for its use (see Rocke and Lorenzato (1995)). The laboratory could benefit from borrowing strength from other laboratories using the same analytical procedure. Both the point estimate and confidence interval would be more realistic if they accounted for interlaboratory variation, thereby reducing the potential bias of the laboratory. Second, it is common in routine analytical practice for a sample to be “split” into a series of subsamples sent to different

laboratories (e.g., a state laboratory, a commercial laboratory, and a U.S. Environmental Protection Agency [EPA] laboratory). The methodology presented here provides an estimate of the true concentration and a corresponding confidence interval based on data from multiple laboratories, borrowing strength from all laboratories involved. The main goal of this article is to construct a confidence interval for a true concentration by the fiducial approach using measurements made by a single laboratory or split samples by multiple laboratories. The proposed confidence interval is constructed using the highest empirical density, it provides the target coverage probability. Confidence intervals constructed using fiducial quantities have coverage probabilities close to the nominal level even with small sample sizes, where traditional large-sample frequentist methods fail.

The Rocke and Lorenzato model has become a widely adopted model in analytical chemistry Rocke and Lorenzato (1995). This model incorporates both additive and multiplicative error components, and the model effectively accommodates heteroscedasticity in measurement errors across the calibration range. Specifically, it accounts for a proportional increase in error with increasing analyte concentration, and a nearly constant error at low analyte levels. This dual-component approach has demonstrated robust fitting across a broad range of analyte concentrations, and the model has garnered significant attention in the literature, leading to various adaptations and extensions; Gibbons and Bhaumik (2001), Bhaumik and Gibbons (2005), and Zhao et al. (2021).

While originally formulated for single-laboratory data, adaptations of the Rocke-Lorenzato model have extended its utility to scenarios involving data from multiple laboratories. Gibbons and Bhaumik (2001) generalized the model to the interlaboratory scenario by using a weighted random-effects regression model to accommodate non-constant variability over the observed concentration range. This extension introduces complexities such as het-

eroscedasticity across laboratories, where variances differ across labs. Calibration, a primary application of the model, involves estimating unknown analyte concentrations based on the corresponding instrument responses, such as spectroscopic measurements. Further, it builds an interlaboratory confidence region for a true concentration, given a series of new measured concentrations from multiple laboratories. Laboratory data used for estimating an unknown true concentration may be from a single laboratory, or from a subset of the laboratories used in estimating the calibration function parameters. For example, all qualified laboratories in a state could be used for estimating the calibration function and associated variances. Borrowing strength from these estimates, we can obtain both point and interval estimates of a true unknown concentration.

Previous work in this area for linear models includes studies by Osborne (1991); Eberhardt and Mee (1994); Johnson and Krishnamoorthy (1996); Cameron et al. (2022). In addition, data analysis using the Rocke-Lorenzato model and its derivatives typically employs likelihood-based methods, maximizing the integrated likelihood to estimate parameters and their asymptotic variances. Bhaumik and Gibbons (2005) obtained moment-based estimators of the model parameters, studied their properties, and used them to construct confidence intervals using data from multiple laboratories. However, conventional large-sample confidence intervals derived from these methods often exhibit inadequate coverage probabilities, unless sample sizes are large. In response to this limitation, Zhao et al. (2021) proposed a modified likelihood approach using a parametric bootstrap based on a standardized signed log-likelihood ratio test to improve confidence interval accuracy. Their method simplifies likelihood functions for heteroscedastic multi-laboratory scenarios and extends to models with varying calibration parameters, achieving improved coverage probabilities.

In the present work, we adopt an alternative approach based on the concept of a fiducial

quantity, which has seen renewed interest in recent years; Hannig et al. (2016) and Bhaumik et al. (2021). Earlier, Weerahandi (1993), Xu and Wang (2008) and Bhaumik et al. (2013) developed and successfully applied the fiducial approach, referring to a generalized pivotal quantity.

The paper is organized as follows: Section 2 presents two motivating applications, followed by a non-linear mixed-effects model for analyzing inter-laboratory calibration data in Section 3. Section 4 introduces the fiducial approach, including fiducial quantities, confidence intervals, and point estimates. Section 5 compares our fiducial approach with other approaches, including the bootstrap, Wald-type large sample method through simulations, while Section 6 revisits the applications. The paper concludes with discussion in Section 7.

2 Motivating Examples

Our first example is an interlaboratory study to analyze cadmium concentrations in water. The example is based on a blind study carried out by the Ford Motor Company, in collaboration with Michigan State Drinking Water Certification laboratories. The study used blind samples prepared by an independent source, randomized and submitted weekly over five weeks. Cadmium was analyzed using inductively coupled plasma atomic emissions spectroscopy (EPA method 200.7). The dataset included five replicates at concentrations of 0, 20, and 100 $\mu\text{g}/L$ across five laboratories.

Our second example is based on an interlaboratory study, also conducted by the Ford Motor Company, in order to analyze copper concentration in water (John Phillips of the Ford Motor Company, personal communication). These data were generated as part of a blind interlaboratory study among the laboratories that hold Michigan State Drinking Water Certifications for the parameters tested. Samples were prepared by an independent

source, randomized, and submitted on a weekly basis over a five-week period. Copper was analyzed by Inductively Coupled Plasma Atomic Emissions Spectroscopy (ICP/AES) using Environmental Protection Agency method 200.7. The dataset consisted of five replicates at each of five concentrations (0 , 2, 10, 50, and 200 $\mu\text{g/L}$) in each of the seven laboratories.

The goals of the above projects were (a) to estimate the calibration curve for each laboratory, (b) estimate the true unknown analyte concentration, and construct corresponding $100(1 - \alpha)\%$ confidence intervals based on data from multiple laboratories, (c) illustrate the use of the confidence interval in environmental monitoring applications where measurements are routinely compared to regulatory thresholds (RT) without incorporating uncertainty in their estimated concentration.

3 Choice of Model

Increasing measurement variation with increasing analyte concentration is commonly observed in analytical data. Denoting the true concentration of an analyte by x , the traditional simple linear calibration model specified as $y = \alpha + \beta x + e$, with the standard normality assumption on the errors, is inappropriate, as it fails to explain increasing measurement variation with increasing analyte concentration. An alternative model can be a log-linear model-for example, $y = xe^\eta$, where η is a normal variable with mean 0 and standard deviation σ_η . However, the main problem with this model is that it fails to explain near-constant measurement variation for a low true concentration level x . Rocke and Lorenzato (1995) proposed a model that better explains the calibration curve, it has two types of errors: an error applicable to low concentration levels and another error for larger concentrations. Accordingly, Gibbons and Bhaumik (2001) consider the following model for the interlaboratory calibration data denoted by y_{ijk} , the k th measurement at the j th true

concentration level x_j in the i th laboratory, $i = 1, 2, \dots, q, j = 1, 2, \dots, r, k = 1, 2, \dots, N_{ij}$:

$$y_{ijk} = \alpha_i + \beta_i x_j \exp(\eta_{ijk}) + \epsilon_{ijk}, \quad (1)$$

where x_j 's are the true concentration values, $j = 1, \dots, r$. We assume that $\eta_{ijk} \stackrel{i.i.d}{\sim} N(0, \sigma_\eta^2)$, $\epsilon_{ijk} \stackrel{i.i.d}{\sim} N(0, \sigma_\epsilon^2)$ and η_{ijk} and ϵ_{ijk} are independent. Also, the α_i and β_i are the calibration parameters for the i th laboratory. We further assume that observations corresponding to zero concentration (i.e. $x = 0$) are available which is often the case in such experiments Bhaumik and Gibbons (2005). The ϵ_{ijk} 's are the additive errors that are present primarily to explain the variability at low-level concentrations. On the other hand, the η_{ijk} 's represent proportional error at higher true concentrations, and we assume that the distributions of the proportional errors remain the same for all laboratories. A model similar to (1) for a single laboratory was proposed by Rocke and Lorenzato (1995) and illustrated using a gas-chromatography/mass-spectrometry example.

Note that under the model (1) corresponding to only one laboratory (i.e., $\beta_i = \beta$ and $\alpha_i = \alpha$), the data near 0 (i.e., $x \simeq 0$) can be used to estimate σ_ϵ^2 , and data for larger concentrations can be used to estimate σ_η^2 . Model (1) also indicates that the variability for observations at larger concentrations is mainly explained by $e^{\eta_{ijk}}$, which follows a lognormal distribution, and observations at low concentrations (i.e., $x \simeq 0$) are normally distributed, which is consistent with our experience Bhaumik and Gibbons (2005).

Parameter Estimation: Multiple estimation methods have been proposed in the literature to estimate the parameters of the model (1). For a single interlaboratory model, Rocke and Lorenzato (1995) used maximum likelihood and applied the technique to toluene by gas-chromatography/mass-spectrometry, and cadmium by atomic absorption spectroscopy. For multiple laboratories, Gibbons and Bhaumik (2001) estimated model parameters by iteratively reweighted maximum marginal likelihood estimation. Performance of the method

was examined via a simulation study and then it was applied to a typical interlaboratory calibration example for copper in distilled water. Employing the concept of data separation technique, i.e. separately considering the role of observations with near zero concentration and that of all other observations corresponding to higher concentrations, Bhaumik and Gibbons (2005) used the moment estimation method to estimate all model parameters. The main advantage of this method is that it gives an explicit expression for each parameter estimate, including the variance components. Inferential procedures for the calibration parameters, and confidence intervals for the true unknown concentrations were developed using large sample properties, and implemented via Wald statistics. Observing some limitations of the standard likelihood based methodology in small sample size situations Zhao et al. (2021), proposed a modified version of the signed log-likelihood ratio test statistic. The modified signed log-likelihood ratio test statistic provides a satisfactory approach for developing confidence intervals for the model parameters and the true unknown concentrations. The goal of this paper is to implement the fiducial approach to achieve accurate results for (i) determining if an analyte is present in a new sample, and (ii) deriving confidence bounds for a true concentration given new measurement(s).

4 Fiducial Approach

The notion of fiducial inference is based on a fiducial quantity, introduced by Fisher in the 1930s. Unlike Bayesian inference, it provides probabilistic statements about the unknown parameters without relying on prior distributions. The fiducial quantity is a function of random variables and corresponding observed values such that (i) the distribution of the fiducial quantity given the data does not depend on any unknown parameters, and (ii) the observed value of the fiducial quantity obtained by replacing each random variable by the

corresponding observed data is free of the nuisance parameters, and is often equal to the parameter of interest. Applications based on fiducial quantities are available in numerous articles. In particular, in the context of lognormal distributions, such applications are given in Krishnamoorthy and Mathew (2003), Bebu and Mathew (2008), and Bhaumik et al. (2021).

Remark: Interest in the fiducial approach was revived by the seminal works of Tsui and Weerahandi (1989) and Weerahandi (1993); a fiducial quantity was referred to as a generalized pivotal quantity in the latter work. Several later authors have developed and successfully applied the fiducial approach, referring to the fiducial quantity as a generalized pivotal quantity; see, for example, the articles by Krishnamoorthy and Mathew (2003), Bebu and Mathew (2008), Xu and Wang (2008), Bhaumik et al. (2013) and the book by Krishnamoorthy and Mathew (2009). The fiducial quantities we have considered in our work are in the spirit of the generalized pivotal quantities considered by these authors. We have referred to them as fiducial quantities. However, what we have developed are approximate fiducial quantities since their implementation relies on the consistency of the maximum likelihood estimate of the model parameter σ_η . In what follows, we derive fiducial quantities for each model parameter and for an unknown true concentration.

4.1 Fiducial Quantities for Model Parameters

To obtain fiducial quantities, we rewrite the model as follows:

$$y_{ijk} = \alpha_i + \beta_i x_j \exp(\sigma_\eta z_{ijk}^\eta) + \sigma_\epsilon z_{ijk}^\epsilon, \quad (2)$$

where $z_{ijk}^\eta \stackrel{i.i.d}{\sim} N(0, 1)$, $z_{ijk}^\epsilon \stackrel{i.i.d}{\sim} N(0, 1)$ and $z_{ijk}^\eta, z_{ijk}^\epsilon$ are independent. Note that (2) has observed data on the left-hand side and a function of model parameters along with random quantities on the right-hand side.

4.1.1 Initial Fiducial Quantities for α_i and σ_ϵ

First, we obtain the fiducial quantities of α_i , $i = 1, \dots, q$, and σ_ϵ based on the data corresponding to $x_j = 0$. In the next subsections, we update these fiducial quantities based on the complete data (corresponding to all concentrations x_j) to obtain updated fiducial quantities of α_i , $i = 1, \dots, q$, and σ_ϵ . We assume $x_1 = 0$ and for $j = 1$ the model (2) can

be written as $y_{i1k} = \alpha_i + \sigma_\epsilon z_{i1k}^\epsilon$. By taking an average over all replications from laboratory i we obtain, $\frac{1}{N_{i1}} \sum_{k=1}^{N_{i1}} y_{i1k} = \alpha_i + \sigma_\epsilon \frac{1}{N_{i1}} \sum_{k=1}^{N_{i1}} z_{i1k}^\epsilon$. Thus, we obtain, $\sum_{k=1}^{N_{i1}} \left(y_{i1k} - \frac{1}{N_{i1}} \sum_{k=1}^{N_{i1}} y_{i1k} \right)^2 = \sigma_\epsilon^2 \sum_{k=1}^{N_{i1}} \left(z_{i1k}^\epsilon - \frac{1}{N_{i1}} \sum_{k=1}^{N_{i1}} z_{i1k}^\epsilon \right)^2$. By summing both sides over all laboratories, so as to utilize in-

formation across all labs, we get $\sum_{i=1}^q \sum_{k=1}^{N_{i1}} \left(y_{i1k} - \frac{1}{N_{i1}} \sum_{k=1}^{N_{i1}} y_{i1k} \right)^2 = \sigma_\epsilon^2 \sum_{i=1}^q \sum_{k=1}^{N_{i1}} \left(z_{i1k}^\epsilon - \frac{1}{N_{i1}} \sum_{k=1}^{N_{i1}} z_{i1k}^\epsilon \right)^2$.

Thus, an initial fiducial quantity for σ_ϵ , say $\tilde{\sigma}_\epsilon^{init}$, can be obtained by solving the above

equation for σ_ϵ : $\tilde{\sigma}_\epsilon^{init} = \sqrt{\frac{\sum_{i=1}^q \sum_{k=1}^{N_{i1}} \left(y_{i1k} - \frac{1}{N_{i1}} \sum_{k=1}^{N_{i1}} y_{i1k} \right)^2}{\sum_{i=1}^q \sum_{k=1}^{N_{i1}} \left(z_{i1k}^\epsilon - \frac{1}{N_{i1}} \sum_{k=1}^{N_{i1}} z_{i1k}^\epsilon \right)^2}}$. In the expression for $\tilde{\sigma}_\epsilon^{init}$ given above, the y_{i1k} s are to be treated as fixed, and the random quantities are the random variable z_{i1k}^ϵ .

An initial fiducial quantity for α_i can be obtained as $\tilde{\alpha}_i^{init} = \frac{1}{N_{i1}} \sum_{k=1}^{N_{i1}} y_{i1k} - \tilde{\sigma}_\epsilon^{init} \times \frac{1}{N_{i1}} \sum_{k=1}^{N_{i1}} z_{i1k}^\epsilon$.

Remark. In the expression for $\tilde{\sigma}_\epsilon^{init}$ given above, the quantity $\sum_{i=1}^q \sum_{k=1}^{N_{i1}} \left(z_{i1k}^\epsilon - \frac{1}{N_{i1}} \sum_{k=1}^{N_{i1}} z_{i1k}^\epsilon \right)^2$ clearly has a chisquare distribution. However, while implementing the fiducial methodology, we will be generating the independent standard normal quantities z_{i1k}^ϵ . The reason is that these quantities also appear in later expressions; see for example equation (3).

4.1.2 Fiducial Quantities for β_i

We construct a total of q structural equations, one for each β_i , $i = 1, \dots, q$. These equations are linear in each β_i , $i = 1, \dots, q$; thus, solutions are obtained in closed forms.

We use (2) to construct structural equations for each β_i , $i = 1, \dots, q$, by averaging (2) over all replications corresponding to different concentrations, $\frac{1}{r} \sum_{j=1}^r \frac{1}{N_{ij}} \sum_{k=1}^{N_{ij}} y_{ijk} = \alpha_i +$

$\beta_i \frac{1}{r} \sum_{j=1}^r \frac{1}{N_{ij}} \sum_{k=1}^{N_{ij}} x_j \exp(\sigma_\eta z_{ijk}^\eta) + \sigma_\epsilon \frac{1}{r} \sum_{j=1}^r \frac{1}{N_{ij}} \sum_{k=1}^{N_{ij}} z_{ijk}^\epsilon$, $i = 1, \dots, q$. Parameters α_i , $i = 1, \dots, q$, and σ_ϵ^2 in these structural equations are replaced by their initial fiducial quantities, $\tilde{\alpha}_i^{init}$, $i = 1, \dots, q$, and $\tilde{\sigma}_\epsilon^{init}$ respectively. We replace σ_η^2 with its consistent point estimate $\hat{\sigma}_\eta^2$ which can be obtained by the maximum likelihood or the method of moments proposed by Bhaumik et al. (2004). Using $\hat{\sigma}_\eta^2$, the fiducial quantities that we exhibit for the β_i , $i = 1, \dots, q$, are only approximate fiducial quantities; however, we shall henceforth refer to the approximate fiducial quantities simply as fiducial quantities. Thus, the fiducial pivots for β_i , $i = 1, \dots, q$, can be obtained as

$$\tilde{\beta}_i = \frac{\sum_{j=1}^r \sum_{k=1}^{N_{ij}} y_{ijk} - r N_{ij} \tilde{\alpha}_i^{init} - \tilde{\sigma}_\epsilon^{init} \sum_{j=1}^r \sum_{k=1}^{N_{ij}} z_{ijk}^\epsilon}{\sum_{j=1}^r \sum_{k=1}^{N_{ij}} x_j \exp(\hat{\sigma}_\eta z_{ijk}^\eta)}, \quad i = 1, \dots, q. \quad (3)$$

4.1.3 Updated Fiducial Quantities for α_i

We can re-write the model (2) as follows for $\{y_{ijk} : x_j \neq 0, i = 1, \dots, q; j = 1, \dots, r; k = 1, \dots, N_{ij}\}$; since we assume $x_1 = 0$, this set is the same as $\{y_{ijk} : i = 1, \dots, q; j = 2, \dots, r; k = 1, \dots, N_{ij}\}$.

$$\frac{y_{ijk}}{x_j} = \frac{\alpha_i}{x_j} + \beta_i \exp(\sigma_\eta z_{ijk}^\eta) + \sigma_\epsilon \frac{z_{ijk}^\epsilon}{x_j}, \quad (4)$$

The method for constructing equations that lead to fiducial quantities for α_i is similar to that for β_i (for $i = 1, \dots, q$), as discussed in the previous subsection. In the following, we construct a total of q equations, one for each α_i , by summing equation (4) over all replications at different non-zero concentrations for each lab, $\sum_{j=2}^r \sum_{k=1}^{N_{ij}} \frac{y_{ijk}}{x_j} = \alpha_i N_{ij} \sum_{j=2}^r \frac{1}{x_j} + \beta_i \sum_{j=2}^r \sum_{k=1}^{N_{ij}} \exp(\sigma_\eta z_{ijk}^\eta) + \sigma_\epsilon \sum_{j=2}^r \sum_{k=1}^{N_{ij}} \frac{z_{ijk}^\epsilon}{x_j}$, $i = 1, \dots, q$. Again, these are linear in each α_i , $i = 1, \dots, q$. We now replace β_i , $i = 1, \dots, q$, and σ_ϵ^2 by $\tilde{\beta}_i$, $i = 1, \dots, q$, and $\tilde{\sigma}_\epsilon^{init}$ respectively. We replace σ_η^2 with its consistent point estimate $\hat{\sigma}_\eta^2$. Thus, the updated fiducial pivots for α_i ,

$i = 1, \dots, q$, can be obtained as

$$\tilde{\alpha}_i = \frac{\sum_{j=2}^r \sum_{k=1}^{N_{ij}} \frac{y_{ijk}}{x_j} - \tilde{\beta}_i \sum_{j=2}^r \sum_{k=1}^{N_{ij}} \exp(\hat{\sigma}_\eta z_{ijk}^\eta) - \tilde{\sigma}_\epsilon^{init} \sum_{j=2}^r \sum_{k=1}^{N_{ij}} \frac{z_{ijk}^\epsilon}{x_j}}{N_{ij} \sum_{j=2}^r \frac{1}{x_j}}, \quad i = 1, \dots, q. \quad (5)$$

4.1.4 Fiducial Quantity for σ_η^2 and Updated Fiducial Quantity σ_ϵ^2

In this section, we construct two structural equations and solve them jointly for σ_η^2 and σ_ϵ^2 to find their respective fiducial quantities. As these two parameters represent variances of random effects, appropriate sums of squares based on the observed data $\{y_{ijk} : i = 1, \dots, q; j = 1, \dots, r; k = 1, \dots, N_{ij}\}$, can be used to construct structural equations. From the model (1), the variance of y_{ijk} can be expressed as $\beta_i^2 x_j^2 (e^{\sigma_\eta^2} - 1) e^{\sigma_\eta^2} + \sigma_\epsilon^2$. Observe that for high values of x_j , the variability in the data inflates significantly. This inflates the variability in the solutions of the structural equations, leading to instability in root-finding algorithms and resulting in wide confidence intervals. To avoid this issue, we shall use the model (4) as the structural equation for non-zero concentrations and (2) for the zero concentration. We treat $\{\frac{y_{ijk}}{x_j} : i = 1, \dots, q; j = 2, \dots, r; k = 1, \dots, N_{ij}\}$ as observed quantities. Two structural equations are constructed based on these observed quantities—(i) using combined within-lab variabilities for all labs (structural equation (6)) and (ii) using combined within-concentration variabilities for all concentrations (structural equation (7)). Due to this reason, while constructing structural equation (7), we also utilized within concentration variability of the observations corresponding to zero concentrations, which are denoted by $\{y_{i1k} : i = 1, \dots, q; k = 1, \dots, N_{i1}\}$. The structural equations for σ_η^2 and σ_ϵ^2 can be expressed as follows:

$$\sum_{i=1}^q \sum_{j=2}^r \sum_{k=1}^{N_{ij}} \left(\frac{y_{ijk}}{x_j} - \frac{1}{r-1} \sum_{j=2}^r \frac{1}{N_{ij}} \sum_{k=1}^{N_{ij}} \frac{y_{ijk}}{x_j} \right)^2 = \sum_{i=1}^q \sum_{j=2}^r \sum_{k=1}^{N_{ij}} \left(\frac{\tilde{\alpha}_i}{x_j} + \tilde{\beta}_i \exp(\sigma_\eta z_{ijk}^\eta) + \sigma_\epsilon \frac{z_{ijk}^\epsilon}{x_j} - \frac{1}{r-1} \sum_{j=2}^r \frac{1}{N_{ij}} \sum_{k=1}^{N_{ij}} \left(\frac{\tilde{\alpha}_i}{x_j} + \tilde{\beta}_i \exp(\sigma_\eta z_{ijk}^\eta) + \sigma_\epsilon \frac{z_{ijk}^\epsilon}{x_j} \right) \right)^2, \quad (6)$$

$$\begin{aligned}
& \sum_{i=1}^q \sum_{k=1}^{N_{i1}} \left(y_{i1k} - \frac{1}{q} \sum_{i=1}^q \frac{1}{N_{ij}} \sum_{k=1}^{N_{i1}} y_{i1k} \right)^2 + \sum_{j=2}^r \sum_{i=1}^q \sum_{k=1}^{N_{ij}} \left(\frac{y_{ijk}}{x_j} - \frac{1}{q} \sum_{i=1}^q \frac{1}{N_{ij}} \sum_{k=1}^{N_{ij}} \frac{y_{ijk}}{x_j} \right)^2 = \\
& \sum_{i=1}^q \sum_{k=1}^{N_{i1}} \left(\tilde{\alpha}_i + \sigma_\epsilon z_{i1k}^\epsilon - \frac{1}{q} \sum_{i=1}^q \frac{1}{N_{i1}} \sum_{k=1}^{N_{i1}} (\tilde{\alpha}_i + \sigma_\epsilon z_{i1k}^\epsilon) \right)^2 + \sum_{j=2}^r \sum_{i=1}^q \sum_{k=1}^{N_{ij}} \\
& \left(\frac{\tilde{\alpha}_i}{x_j} + \tilde{\beta}_i \exp(\sigma_\eta z_{ijk}^\eta) + \sigma_\epsilon \frac{z_{ijk}^\epsilon}{x_j} - \frac{1}{q} \sum_{i=1}^q \frac{1}{N_{ij}} \sum_{k=1}^{N_{ij}} \left(\frac{\tilde{\alpha}_i}{x_j} + \tilde{\beta}_i \exp(\sigma_\eta z_{ijk}^\eta) + \sigma_\epsilon \frac{z_{ijk}^\epsilon}{x_j} \right) \right)^2. \quad (7)
\end{aligned}$$

In the above structural equations, $\alpha_i, \beta_i, i = 1, \dots, q$, are replaced by their corresponding fiducial quantities $\tilde{\alpha}_i, \tilde{\beta}_i, i = 1, \dots, q$, respectively. These structural equations are nonlinear in σ_η^2 and σ_ϵ^2 , and can be solved numerically to obtain a fiducial quantity for σ_η^2 , and an updated fiducial quantity for σ_ϵ^2 ; denoted by $\tilde{\sigma}_\eta^2$ and $\tilde{\sigma}_\epsilon^2$, respectively.

4.2 Fiducial Quantity for an Unknown Concentration

Let $y_{ijk}^*, i = 1, \dots, q^*, j = 1, \dots, r^*, k = 1, \dots, N_{ij}^*$, denote observations of new samples from laboratory i for an unknown concentration $\mathcal{X}_j$ at replication k . We assume that the laboratories where these samples are obtained are among the q laboratories used to estimate the model parameters in Model (1). In the previous section, we obtained the fiducial quantities of all the model parameters. If we replace all the model parameters with their corresponding fiducial quantities in the structural equation (2), it becomes a function of the unknown concentrations. We rewrite (2) for the new samples as

$$y_{ijk}^* = \alpha_i + \beta_i \mathcal{X}_j \exp(\sigma_\eta z_{ijk}^{\eta*}) + \sigma_\epsilon z_{ijk}^{\epsilon*} \quad (8)$$

where, $z_{ijk}^{\eta*} \stackrel{i.i.d}{\sim} N(0, 1)$, $z_{ijk}^{\epsilon*} \stackrel{i.i.d}{\sim} N(0, 1)$ and $z_{ijk}^{\eta*}, z_{ijk}^{\epsilon*}$ are independent. A structural equation for an unknown concentration $\mathcal{X}_j, j = 1, \dots, r^*$, can be obtained as follows by summing (8) over all replications and all laboratories, $\sum_{i=1}^{q^*} \sum_{k=1}^{N_{ij}^*} y_{ijk}^* = \sum_{i=1}^{q^*} \sum_{k=1}^{N_{ij}^*} \alpha_i + \mathcal{X}_j \sum_{i=1}^{q^*} \sum_{k=1}^{N_{ij}^*} \beta_i \exp(\sigma_\eta z_{ijk}^{\eta*}) + \sum_{i=1}^{q^*} \sum_{k=1}^{N_{ij}^*} \sigma_\epsilon z_{ijk}^{\epsilon*}$. We solve this for $\mathcal{X}_j$ while replacing $\alpha_i, \beta_i, i =$

$1, \dots, q^*$, σ_η and σ_ϵ by their corresponding fiducial quantities $\tilde{\alpha}_i$, $\tilde{\beta}_i$, $i = 1, \dots, q$, $\tilde{\sigma}_\eta$ and $\tilde{\sigma}_\epsilon$ to obtain fiducial quantity of a new concentration $\mathcal{X}_j$, $j = 1, \dots, r^*$, as follows:

$$\tilde{\mathcal{X}}_j = \frac{\sum_{i=1}^{q^*} \sum_{k=1}^{N_{ij}^*} y_{ijk}^* - \sum_{i=1}^{q^*} \sum_{k=1}^{N_{ij}^*} \tilde{\alpha}_i - \sum_{i=1}^{q^*} \sum_{k=1}^{N_{ij}^*} \tilde{\sigma}_\epsilon z_{ijk}^{\epsilon^*}}{\sum_{i=1}^{q^*} \sum_{k=1}^{N_{ij}^*} \tilde{\beta}_i \exp(\tilde{\sigma}_\eta z_{ijk}^{\eta^*})}. \quad (9)$$

4.3 Fiducial Confidence Interval and Fiducial Point Estimate

So far we have discussed a method of constructing fiducial quantities of the parameters under the model in Subsection 4.1. Furthermore, in subsection 4.2, we have discussed how to obtain fiducial quantities for the unknown true concentrations corresponding to new measurements utilizing the fiducial quantities obtained in subsection 4.1. In this section, we conclude our discussion on fiducial inference for the unknown true concentrations by discussing how to utilize these fiducial quantities to obtain point and interval estimates of the unknown concentrations.

Since fiducial quantities are functions of the observed data and random quantities, whose distributions are free of any unknown parameters, samples from the fiducial quantities can be obtained simply by drawing samples from the corresponding distributions. These samples can then be used to estimate the density of the fiducial distribution for each parameter. Once the density is estimated, it serves as a basis for determining both point and interval estimates for each parameter. By examining the estimated density, we can derive point estimates such as the mode, and the highest density intervals, which offer a measure of uncertainty around these estimates. In this work, we focus on the point and interval estimation of the unknown concentration from new samples $\{y_{ijk}^* : i = 1, \dots, q^*, j = 1, \dots, r^*, k = 1, \dots, N_{ij}^*\}$, where the fiducial quantities of the model parameters have been obtained from the earlier sample $\{y_{ijk} : i = 1, \dots, q, j = 1, \dots, r, k = 1, \dots, N_{ij}\}$, correspond-

ing to known concentrations. For an unknown concentration $\mathcal{X}_j$ (for any $j = 1, \dots, r^*$), one sample from its fiducial distribution can be drawn sequentially based on the discussion in sections 4.1 and 4.2 (we included a summary of the sequential procedure in the supplementary material). One can generate N fiducial quantities of each unknown concentration by replicating these steps N times. Let us denote the density of the fiducial distribution of the unknown concentration $\mathcal{X}_j$ as $d_{\mathcal{X}_j}(\cdot)$ which can be estimated accurately from the N generated fiducial quantities. In the following, we use the estimated 95% fiducial highest density interval as the interval estimate and the fiducial mode as the point estimate of $\mathcal{X}_j$.

100(1 − α)% Fiducial Highest Density Interval (HDI). The 100(1 − α)% fiducial HDI for $\mathcal{X}_j$ is the interval $[L, U]$ such that: $\int_L^U d_{\mathcal{X}_j}(x)dx = 1 - \alpha$ and minimizes the interval length $(U - L)$. The endpoints L and U are determined such that $d_{\mathcal{X}_j}(L) = d_{\mathcal{X}_j}(U)$ ensuring that the interval captures the highest density region. **Fiducial Mode.** The fiducial mode $\hat{\mathcal{X}}_j$ is the value of $\mathcal{X}_j$ that maximizes the fiducial density function; i.e., $\hat{\mathcal{X}}_j = \arg \max_x d_{\mathcal{X}_j}(x)$. This point estimate $\hat{\mathcal{X}}_j$ represents the value of $\mathcal{X}_j$ with the highest probability density under the fiducial distribution.

5 Simulation Studies

To illustrate some key features of the proposed confidence interval for a true unknown concentration, we shall use a few different simulation scenarios. For each simulation scenario, training and test samples (as described in Section 4.4) are generated using two distinct configurations. We label the two configurations for training sample generation as “A” and “B,” and the two configurations for test sample generation as “1” and “2.” This results in a total of four simulation scenarios: 1.A, 2.A, 1.B, and 2.B. The difference between configurations A and B lies in the generation of the training samples. Configuration A includes

3 true concentrations and 3 laboratories, while configuration B uses 9 true concentrations and 10 laboratories. On the other hand, configurations 1 and 2 differ in how the test samples are generated. In configuration 1, multiple replications are drawn from a single laboratory to estimate the unknown concentration, whereas configuration 2 involves a single replication from multiple laboratories. In all scenarios, we compare the performance of our proposed fiducial method for interval estimation with several alternative methods: the parametric bootstrap, the adjusted likelihood ratio test (LRT) statistics proposed by Zhao et al. (2021), Wald-type large-sample method developed by Bhaumik and Gibbons (2005) using Method of Moment Estimator (MME), and Wald-type method based on the MLE.

5.1 Simulation Setup

Details of training samples for simulation Scenarios 1.A and 2.A. Data were simulated for three concentrations (0, 10, 30) $\mu\text{g/L}$ with $N_{ij} = 5$ replicates per concentration nested within $q = 3$ laboratories. Values of the model (1) parameters were set at: $\alpha_i = 1$, $\beta_i = 1$, $\sigma_e = 1$, and $\sigma_\eta = 0.1$. This setup corresponds to a small sample size as only three laboratories and three concentrations per laboratory have been considered.

Details of training samples for simulation Scenarios 1.B and 2.B. Data were simulated with a moderate sample size for nine concentrations (0, 5, 10, 15, 20, 25, 30, 35, 40) $\mu\text{g/L}$ with $N_{ij} = 5$ replicates per concentration nested within $q = 10$ laboratories. Values of the model parameters were set at: $\alpha_i = 0$, $\beta_i = 1$, $\sigma_e = 1$, and $\sigma_\eta = 0.1$.

Details of test samples for simulation Scenarios 1.A and 1.B. Samples for an unknown concentration consist of 5 replications from one laboratory (without loss of generality, say Laboratory 1); with $q^* = 1$, $N_{ij}^* = 5$. For Scenario 1.A, true values of unknown concentrations are (5, 20, 50). In Scenario 1.B, true values of unknown concentrations are

(0, 5, 10, 15, 20, 25, 30, 35, 40, 50) $\mu\text{g/L}$.

Details of test samples for simulation Scenarios 2.A and 2.B. The samples for an unknown concentration consist of a single replication from multiple laboratories. In scenario 2.A, we use three laboratories in the training sample, i.e., $q^* = 3$, $N_{ij}^* = 1$, while, in scenario 2.B, $q^* = 6$, $N_{ij}^* = 1$. The true values of unknown concentrations for 2.A and 2.B are the same as 1.A and 1.B, respectively.

5.2 Simulation Results

The simulation results for each scenario are based on 5,000 simulated datasets. For each dataset, the proposed fiducial method estimates the unknown concentration using 10,000 samples generated from the fiducial quantity. Similarly, 10,000 bootstrap samples are utilized in both the parametric bootstrap and LRT-based methods for comparison. Notably, the LRT-based method is computationally more intensive than the parametric bootstrap method. However, the proposed fiducial method provides significantly more computational efficiency compared to both approaches. The simulation results for interval estimation are presented in Tables 1, 2, 3, 4. For simulation scenarios 1.B and 2.B (Tables 2 and 4), results for concentrations 0, 15, 30, and 50 are presented in this document (Complete tables can be found in Supplementary Material). The performance of the interval estimate is evaluated based on both the average width of the 95% confidence interval and its coverage probability. From the numerical results, we conclude that for interval estimation, the proposed fiducial approach outperforms the other methods considered. As already noted, these methods include (i) the parametric bootstrap, (ii) the modified LRT statistic proposed by Zhao et al. (2021), (iii) Wald-type large-sample method developed by Bhaumik and Gibbons (2005) using MME, and (iv) Wald-type method based on the MLE. When the size of the training

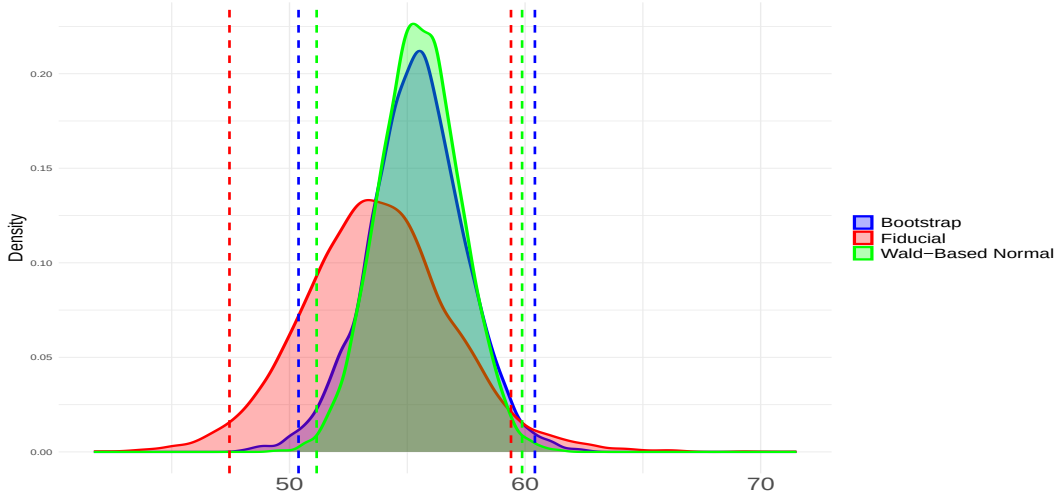


Figure 1: Density plots of the fiducial distribution, bootstrap distribution and Wald-Based normal distribution for an unknown concentration (true value is 50) for a sample from simulation setting 2.A where the MLE is 55.6. Vertical dashed lines shows 95% confidence intervals. The figure highlights the fact that only the fiducial interval includes the true value: a conclusion that agrees with the poor coverage probabilities of the other intervals, as reported in this section.

sample is small (as seen in simulation Scenarios 1.A and 2.A), the coverage probability of the proposed 95% fiducial confidence interval is approximately 0.95. In contrast, other methods such as the bootstrap and Wald, exhibit coverage probabilities that are either significantly smaller than 0.95 or very close to 1 (see Tables 1 and 3), and the expected confidence widths correspondingly reflect the behavior of the coverage probabilities. As the training sample size increases from small to moderate (Scenarios 1.B and 2.B), the performance of the Wald intervals based on the MLE improves, as expected. Nevertheless, even in these scenarios, the coverage probability of the fiducial interval remains the most accurate. In fact we see that the coverage probability associated with the parametric bootstrap solution can be significantly lower than the assumed confidence level. We have provided Figure 1 for illustration for such lower coverage of other intervals.

The simulation results for point estimation are provided in the supplementary material,

Table 1: Estimated coverage probabilities and expected widths of the 95% confidence intervals for the true unknown concentration: Simulation Scenario 1.A

True Conc.	Average 95% Confidence Limits [Average Width]				
	Fiducial	Bootstrap	LRT	Wald MLE	Wald MME
5	(3.8, 6.3) [2.5]	(4.0, 6.4) [2.4]	(3.7, 6.3) [2.6]	(4.0, 6.2) [2.2]	(2.8, 7.2) [4.4]
20	(17.6, 22.7) [5.1]	(17.6, 22.4) [4.8]	(17.8, 22.8) [5.0]	(17.8, 22.3) [4.5]	(15.6, 24.6) [9.0]
50	(44.2, 56.5) [12.3]	(43.3, 55.5) [12.2]	(44.4, 56.5) [12.1]	(44.4, 55.3) [10.9]	(40.1, 60.4) [20.3]
True Conc.	Coverage Probability				
	Fiducial	Bootstrap	LRT	Wald MLE	Wald MME
5	0.949	0.931	0.941	0.931	1.000
20	0.952	0.938	0.889	0.921	1.000
50	0.949	0.918	0.915	0.911	0.994

Table 2: Simulation Scenario 1.B: Estimated coverage probabilities and expected widths of the 95% confidence intervals for the true unknown concentrations 0, 15, 30, 50.

True Conc.	Average 95% Limits [Average Width]				Coverage Probability			
	Fiducial	Bootstrap	Wald MLE	Wald MME	Fiducial	Bootstrap	Wald MLE	Wald MME
0	(0, 1.1) [1.1]	(0, 1.5) [0.9]	(0, 1.3) [1.3]	(0, 1.9) [1.9]	0.948	0.814	0.900	0.996
15	(13.4, 16.7) [3.3]	(13.6, 16.7) [3.1]	(13.4, 16.7) [3.3]	(11.3, 18.6) [7.3]	0.949	0.931	0.943	1.000
30	(27.1, 33.0) [5.9]	(26.8, 32.5) [5.7]	(26.9, 32.7) [5.8]	(23.6, 36.3) [12.7]	0.947	0.914	0.942	1.000
50	(45.3, 55.2) [9.9]	(44.2, 53.6) [9.4]	(44.9, 54.4) [9.5]	(39.8, 60.4) [20.6]	0.949	0.899	0.935	1.000

where we observe that the fiducial point estimates perform similarly to the MLE.

6 Data Analysis

We shall report the results from the analysis of both the cadmium and copper datasets mentioned earlier. One measurement from each laboratory was selected as the test sample for each concentration, while the remaining data were used to train the model. Since we are working with real-world datasets, it is essential to first evaluate the adequacy of the Rocke-Lorenzato model specified in (1) before performing the data analysis to estimate the unknown concentrations. To this end, the training data was used to fit the Rocke-Lorenzato

Table 3: Estimated coverage probabilities and expected widths of the 95% confidence intervals for the true unknown concentration: Simulation Scenario 2.A

True Conc.	Average 95% Limits [Average Width]			Coverage Probability		
	Fiducial	Bootstrap	Wald MLE	Fiducial	Bootstrap	Wald MLE
5	(3.6, 6.4) [2.8]	(3.8, 6.4) [2.6]	(3.8, 6.2) [2.2]	0.949	0.939	0.918
20	(17.3, 22.9) [5.6]	(17.3, 22.6) [5.3]	(17.7, 22.4) [4.7]	0.953	0.938	0.927
50	(43.8, 56.7) [12.9]	(43.2, 55.8) [12.6]	(44.4, 55.6) [11.2]	0.947	0.901	0.895

Table 4: Simulation Scenario 2.B: Estimated coverage probabilities and expected widths of the 95% confidence intervals for the true unknown concentrations 0, 15, 30, 50.

True Conc.	Average 95% Limits [Average Width]			Coverage Probability		
	Fiducial	Bootstrap	Wald MLE	Fiducial	Bootstrap	Wald MLE
0	(0, 0.9) [0.9]	(0, 1.1) [1.1]	(0, 1.0) [1.0]	0.952	0.849	0.936
15	(13.6, 16.5) [2.9]	(13.7, 16.4) [2.7]	(13.6, 16.5) [2.9]	0.948	0.918	0.938
30	(27.5, 32.6) [5.1]	(27.4, 32.2) [4.8]	(27.4, 32.4) [5.0]	0.946	0.932	0.942
50	(46.0, 54.2) [8.2]	(45.4, 53.3) [7.9]	(45.7, 53.7) [8.0]	0.948	0.912	0.932

model (1). Subsequently, 95% upper and lower calibration curves were constructed based on the fitted model, the majority of observed measurements fall within the specified limits, supporting the adequacy of the Rocke-Lorenzato model (1) for these datasets. A detailed goodness-of-fit assessments for the Rocke-Lorenzato calibration model on these datasets is presented in the supplementary material.

After assessing the model fit, the test samples with unknown concentrations are now analyzed, leading to the construction of confidence intervals and point estimates for the unknown concentrations. Table 5 presents point estimates calculated using both the fiducial method and the MLE for each concentration. Accompanying these point estimates are 95% confidence intervals derived from the proposed fiducial method in Section 4, as well as those obtained using the bootstrap and Wald methods based on MLE. It is noteworthy that unlike the fiducial confidence intervals whose coverage probabilities are quite accurate, the

Table 5: Results for Copper and Cadmium data. Each data consists five replicates for each concentration in every laboratory. The first four replicates were used to train the data and the last replicate was used to estimate the concentration.

Cadmium Data						Copper Data					
True	Point Estimate		95% Confidence Interval			True	Point Estimate		95% Confidence Interval		
Conc.	Fiducial	MLE	Fiducial	Bootstrap	MLE	Conc.	Fiducial	MLE	Fiducial	Bootstrap	MLE
0	1.40	1.51	(0, 5.6)	(0, 5.4)	(0, 5.0)	0	0.14	0.20	(0, 2.7)	(0, 2.1)	(0, 2.1)
20	20.5	19.4	(17, 24)	(16, 24)	(16, 23)	2	2.6	2.6	(0.6, 5.0)	(1.1, 4.7)	(0.7, 4.5)
100	103.2	102.8	(91, 107)	(96, 110)	(96, 110)	10	9.1	8.9	(6.5, 11.5)	(6.8, 11.0)	(7.0, 10.8)
						50	47.3	47.1	(42.6, 51.5)	(44.7, 49.8)	(44.7, 49.5)
						200	197.1	195.9	(181, 208)	(189, 205)	(189, 203)

bootstrap and Wald MLE methods result in coverage probabilities that can be significantly below the assumed confidence level. Clearly, this has implications in the application of these methodologies. For example, for the copper data, when the true concentration is 50, both the bootstrap and Wald MLE methods yield upper limits of the 95% interval that fall below 50. This demonstrates the potential unreliability of these methods.

7 Discussion

The fiducial approach developed in this article is unique in that it allows estimates for a single laboratory or multiple laboratories analyzing split samples from the same environmental source to borrow strength from other laboratories for which calibration data using the same analytic methodology are available. The distinct advantage of this approach is that both the point estimate and confidence interval are more precise. The approach presented in this paper addresses a crucial challenge in analytical chemistry, namely the uncertainty in chemical measurements. Uncertainty is frequently underestimated in this field, leading to the erroneous belief that measurement estimates reflect true values. The techniques devised

in this study facilitate the estimation of uncertainty related to the actual true concentration upon which these measurements are predicated. Given that measurements are typically collected from various instruments, analyses, and often across multiple laboratories, it is essential to consider these sources of variance when evaluating the overall uncertainty in the analytical process. Therefore, these variances should be factored into calculations that define the bounds for the underlying true concentration. In practical applications, when assessing a measurement against a regulatory threshold limit, or compliance standard (referred to as RT), we can compute a confidence interval for the true unknown concentration X . An exceedance is declared only if the lower confidence limit (LCL) exceeds the RT, providing $100(1 - \alpha)\%$ confidence that the RT has been surpassed. Thus, determination of the LCL with the prespecified confidence level is extremely important for comparisons to the RT, and our fiducial approach provides an LCL for this purpose.

Two additional noteworthy applications of this method are worth mentioning. First, it can be employed to estimate a detection limit, defined as the concentration threshold where true values exceed zero (see Currie (1968); Danzer and Currie (1998); Gibbons (1996)). In this context, we identify the minimum value of X such that the LCL is greater than zero. Secondly, this method can aid in determining a quantification limit (see Currie (1968); Gibbons et al. (1992)). The goal in this case is to identify the true concentration for which the confidence interval is sufficiently narrow to ensure a reliable quantitative measurement. Currie (1968) defined the “determination limit” as the true concentration at which the percent relative standard deviation (RSD) is 10%. We propose redefining the quantification limit as the value of X where the interval length is 10% of the point estimate of X given Y . Previous methodologies primarily concentrated on single-laboratory confidence bounds on Y given X (see Gibbons (1996); Gibbons and Bhaumik (2001)),

rather than directly computing bounds on X based on single or multiple laboratories, as this paper demonstrates. It is crucial to recognize that achieving a 10% RSD may not always be attainable, contingent upon the variance (σ_η).

SUPPLEMENTARY MATERIAL

Supplementary Material outlines the sequential steps required to generate the fiducial pivot quantities for estimating unknown concentrations, including the procedure for handling new test samples. In addition, the supplement discusses specific cases where traditional methods like MLE and MME become computationally infeasible due to limited sample sizes. It also includes comprehensive simulation results with complete tables and graphical analyses that demonstrate the performance of the fiducial approach covering both interval and point estimation, as well as detailed goodness-of-fit assessments for the Rocke–Lorenzato calibration model on cadmium and copper datasets.

References

- Bebu, I. and Mathew, T. (2008), “Comparing the means and variances of a bivariate lognormal distribution,” *Statistics in Medicine*, 27, 2684–2696.
- Bhaumik, D. K. and Gibbons, R. D. (2005), “Confidence regions for random-effects calibration curves with heteroscedastic errors,” *Technometrics*, 47, 223–231.
- Bhaumik, D. K., Kapur, K., Bhaumik, R., and Reda, D. J. (2013), “Sample size determination and hypothesis testing for the mean of a lognormal distribution,” *Journal of Environmental Statistics*, 5.

- Bhaumik, D. K., Shi, H., Reda, D. J., and Sinha, B. K. (2021), “Generalized Confidence Interval for An Agreement Between Raters,” *Statistics in Medicine*, 40, 2230–2238.
- Cameron, D., Hartig, F., Minnuno, F., Oberpriller, J., Reineking, B., Van Oijen, M., and Dietze, M. (2022), “Issues in calibrating models with multiple unbalanced constraints: the significance of systematic model and data errors,” *Methods in Ecology and Evolution*.
- Currie, L. A. (1968), “Limits for qualitative detection and quantitative determination, application to radiochemistry,” *Analytical Chemistry*, 40, 586–593.
- Danzer, K. and Currie, L. A. (1998), “Guidelines for calibration in analytical chemistry part 1. fundamentals and single component calibration,” *Pure Applied Chemistry*, 70.
- Eberhardt, K. R. and Mee, R. W. (1994), “Constant-width calibration intervals for linear regression,” *Journal of Quality Technology*, 26, 21–29.
- Gibbons, R. D. (1996), “Some conceptual and statistical issues in analysis of groundwater monitoring data,” *Environmetrics*, 7, 125–199.
- Gibbons, R. D. and Bhaumik, D. K. (2001), “Weighted random-effects regression models with application to interlaboratory calibration,” *Technometrics*, 43, 192–198.
- Gibbons, R. D., Grams, N. E., Jarke, F. H., and Stoub, K. P. (1992), “Practical quantification limits,” *Chemometrics and Intelligent Laboratory Systems*, 12, 225–235.
- Hannig, J., Iyer, H., Lai, R. C., and Lee, T. C. (2016), “Generalized Fiducial Inference: A Review and New Results,” *Journal of the American Statistical Association*, 111, 1346–61.
- Johnson, D. and Krishnamoorthy, K. (1996), “Combining independent studies in a calibration problem,” *Journal of the American Statistical Association*, 91, 1707–1715.

- Krishnamoorthy, K. and Mathew, T. (2003), “Inferences on the means of lognormal distributions using generalized p-values and generalized confidence intervals,” *Journal of Statistical Planning and Inference*, 115, 103–121.
- (2009), *Statistical Tolerance Regions: Theory, Applications and Computation*, New York: John Wiley.
- Osborne, C. (1991), “Statistical calibration:a review,” *International Statistical Review*, 59.
- Rocke, D. M. and Lorenzato, S. (1995), “A two-component model for measurement error in analytical chemistry,” *Technometrics*, 37, 176–184.
- Tsui, K. W. and Weerahandi, S. (1989), “Generalized p-values in significance testing of hypotheses in the presence of nuisance parameters,” *Journal of the American Statistical Association*, 84, 602–607.
- Weerahandi, S. (1993), “Generalized Confidence Intervals,” *Journal of the American Statistical Association*, 88, 899–905.
- Xu, L.-W. and Wang, S.-G. (2008), “A new generalized p-value for anova under heteroscedasticity,” *Statistics & Probability Letters*, 78, 963–969.
- Zhao, J., Mathew, T., and Bebu, I. (2021), “Accurate confidence intervals for inter-laboratory calibration and common mean estimation,” *Chemometrics and Intelligent Laboratory Systems*, 208, 104218.

**SUPPLEMENTARY MATERIAL of
‘Fiducial Inference for Random-Effects
Calibration Models: Advancing Reliable
Quantification in Environmental Analytical
Chemistry’**

1 Sequential Steps to Generate Fiducial Pivot Quantity for Unknown Concentrations

We want to generate fiducial pivot quantity for unknown concentrations from new samples $\{y_{ijk}^* : i = 1, \dots, q^*, j = 1, \dots, r^*, k = 1, \dots, N_{ij}^*\}$, where the fiducial quantities of the model parameters have been obtained from the earlier sample $\{y_{ijk} : i = 1, \dots, q, j = 1, \dots, r, k = 1, \dots, N_{ij}\}$, corresponding to known concentrations.

Based on the discussion in sections 4.1 and 4.2 in the paper, for an unknown concentration $\mathcal{X}_j$ (for any $j = 1, \dots, r^*$), one sample from its fiducial distribution can be drawn in the following sequential steps:

1. Generate $z_{ijk}^\eta \overset{i.i.d}{\sim} N(0, 1)$, $z_{ijk}^\epsilon \overset{i.i.d}{\sim} N(0, 1)$, $i = 1, \dots, q$, $j = 1, \dots, r$, $k = 1, \dots, N_{ij}$.
2. Generate initial fiducial quantities for σ_ϵ and α_i , as $\tilde{\sigma}_\epsilon^{init}$ and $\tilde{\alpha}_i^{init}$ respectively, $i = 1, \dots, q$.
3. Generate fiducial quantities for β_i , $i = 1, \dots, q$, using the expression (3) in the main manuscript and obtain updated fiducial quantities for α_i , $i = 1, \dots, q$, using the expression (5) in the main manuscript.
4. Generate the fiducial quantity for σ_η and updated fiducial quantity of σ_ϵ by simultaneously solving the two non-linear structural equations (6) and (7) in the main manuscript.
5. finally, generate $z_{ijk}^\eta \overset{i.i.d}{\sim} N(0, 1)$, $z_{ijk}^\epsilon \overset{i.i.d}{\sim} N(0, 1)$, $i = 1, \dots, q^*$, $k = 1, \dots, N_{ij}^*$, and obtain the fiducial quantity of the unknown concentration $\mathcal{X}_j$ from the expression (9) in the main manuscript for any $j = 1, \dots, r^*$, using the fiducial quantities in the

previous steps. One can generate N fiducial quantities of each unknown concentration by replicating these steps N times.

2 Special Case: MLE and MME Estimates Computationally Infeasible due to Limited Sample Size

As this article focuses on the small sample inference, we highlight some limitations of the traditional methods, such as maximum likelihood estimation (MLE) and method of moment estimation (MME). In the context of a multilaboratory testing problem, estimating unknown concentrations involves two types of samples: (i) Training Sample: it consists of historical data from laboratories with known concentrations, primarily used for estimating the model parameters; (ii) the Test Sample: the new sample for which we aim to estimate unknown concentrations. Note that generally the test sample is much smaller in size than the training sample. The sample size for the training sample is $\sum_{i=1}^q \sum_{j=1}^r N_{ij}$ and the sample size for each unknown concentration in the test sample is $\sum_{i=1}^{q^*} N_{ij}^*$. When the training sample is derived from only a few laboratories and the historical data includes only a limited number of concentrations, parameter estimation using methods such as MLE or MME can become unstable and, at times, computationally infeasible. We illustrate this by an example with the MME estimates proposed in ?. The MME of β_i , $i = 1, \dots, q$, and σ_η are obtained as $\hat{\beta}_i^{MME} = \sqrt{\frac{\mu_{zi}^4}{\sigma_{zi}^2 - \sigma_u^2 + \mu_{zi}^2}}$ and $\hat{\sigma}_\eta^{MME} = \sqrt{\frac{2 \sum_{i=1}^q \ln(\mu_{zi}/\beta_i)}{q}}$, where, $z_{ijk} = (y_{ijk} - \alpha_i)/x_j$ and $u_{ijk} = e_{ijk}/x_j$ for $x_j \neq 0$, $\mu_{zi} = E(z_{ijk})$, $\sigma_{zi}^2 = \frac{\sum_{j=1}^r \text{var}(z_{ijk})}{r}$, and $\sigma_u^2 = \frac{\sigma_e^2 \sum_{j=1}^r 1/x_j^2}{r}$. In the case where σ_η is small (for example, $\sigma_\eta = 0.01$) and σ_e is moderately large (for example, $\sigma_e = 2$), the observed values of z_{ijk} can become negative, and consequently, the average

μ_{zi} can also be negative for positive β_i . In such cases, the numerical value of σ_u^2 may exceed that of $\sigma_{zi}^2 + \mu_{zi}^2$, since the z_i 's are close to zero. As a result, estimates $\hat{\beta}_i^{MME}$ and $\hat{\sigma}_\eta^{MME}$ may be computationally infeasible. Likewise, when the test sample with an unknown concentration comes from a single laboratory with only a few replicates of the measurement, the ML estimates may face similar computational challenges.

The proposed fiducial method is not affected by such computational challenges; however, the proposed method requires a consistent estimate of σ_η . In cases where the MME or MLE of σ_η is computationally infeasible, one can obtain an initial fiducial quantity for σ_η using the data corresponding to high concentrations and substitute that in expression (3) in the main manuscript at the place of $\hat{\sigma}_\epsilon$. When the concentration x_j is high for some $j = j^*$, the model (2) in the main manuscript can be well approximated as $y_{ij^*k} \approx \beta_i x_{j^*} \exp(\sigma_\eta z_{ij^*k}^\eta)$ and the distribution of $\ln(y_{ij^*k}/x_j)$ can be approximated by a normal distribution with mean $\ln(\beta_i)$ and variance σ_η . Thus, an initial fiducial quantity for σ_η can be obtained by solving the structural equation
$$\sum_{i=1}^q \sum_{k=1}^{N_{ij^*}} \left(\ln(y_{ij^*k}/x_j) - \sum_{i=1}^q \sum_{k=1}^{N_{ij^*}} \frac{1}{q^* N_{ij^*}} \ln(y_{ij^*k}/x_j) \right)^2 = \sigma_\eta \sum_{i=1}^q \sum_{k=1}^{N_{ij^*}} \left(z_{ij^*k}^\eta - \sum_{i=1}^q \sum_{k=1}^{N_{ij^*}} \frac{1}{q^* N_{ij^*}} z_{ij^*k}^\eta \right)^2$$
 for σ_η , where the solution is evident.

3 Complete Tables For Interval Estimation Results

4 Simulation Results For Point Estimation

The performance of the point estimate is evaluated based on relative bias, relative absolute bias, and relative root mean squared error, which are defined as follows. For an unknown concentration $\mathcal{X}$ and its point estimate $\hat{\mathcal{X}}_i$ computed from the i th simulated dataset, $i =$

Table C.1: Estimated coverage probabilities and expected widths of the 95% confidence intervals
for the true unknown concentration: Simulation Scenario 1.B

True Conc.	Average 95% Limits [Average Width]				Coverage Probability			
	Fiducial	Bootstrap	Wald MLE	Wald MME	Fiducial	Bootstrap	Wald MLE	Wald MME
0	(0, 1.1) [1.1]	(0, 1.5) [0.9]	(0, 1.3) [1.3]	(0, 1.9) [1.9]	0.948	0.814	0.900	0.996
5	(3.8, 6.1) [2.3]	(4.4, 6.5) [2.1]	(4.1, 6.3) [2.1]	(2.7, 7.1) [4.4]	0.947	0.845	0.921	0.998
10	(8.7, 11.3) [2.6]	(9.0, 11.5) [2.5]	(8.8, 11.4) [2.6]	(7.1, 12.7) [5.6]	0.945	0.918	0.932	0.997
15	(13.4, 16.7) [3.3]	(13.6, 16.7) [3.1]	(13.4, 16.7) [3.3]	(11.3, 18.6) [7.3]	0.949	0.931	0.943	1.000
20	(17.9, 22.1) [4.2]	(18.0, 22.0) [4.0]	(18.0, 22.0) [4.0]	(15.5, 24.5) [9.0]	0.948	0.936	0.941	1.000
25	(22.5, 27.5) [5.0]	(22.5, 27.3) [4.8]	(22.5, 27.3) [4.8]	(19.6, 30.4) [10.8]	0.950	0.914	0.943	1.000
30	(27.1, 33.0) [5.9]	(26.8, 32.5) [5.7]	(26.9, 32.7) [5.8]	(23.6, 36.3) [12.7]	0.947	0.914	0.942	1.000
35	(31.7, 38.6) [6.9]	(31.1, 37.7) [6.6]	(31.5, 38.1) [6.6]	(27.7, 42.4) [14.7]	0.952	0.916	0.945	1.000
40	(36.2, 44.1) [7.9]	(35.6, 43.2) [7.6]	(35.9, 43.5) [7.6]	(31.7, 48.3) [16.6]	0.953	0.920	0.944	1.000
50	(45.3, 55.2) [9.9]	(44.2, 53.6) [9.4]	(44.9, 54.4) [9.5]	(39.8, 60.4) [20.6]	0.949	0.899	0.935	1.000

1, ..., 5000, relative bias is defined as $\frac{1}{5000} \sum_{i=1}^{5000} \frac{(\hat{\mathcal{X}}_i - \mathcal{X})}{\mathcal{X}}$, relative absolute bias is defined as $\frac{1}{5000} \sum_{i=1}^{5000} \frac{(|\hat{\mathcal{X}}_i - \mathcal{X}|)}{\mathcal{X}}$, relative root mean squared error is defined as $\frac{1}{5000} \sqrt{\sum_{i=1}^{5000} \frac{(\hat{\mathcal{X}}_i - \mathcal{X})^2}{\mathcal{X}^2}}$, provided $\mathcal{X} \neq 0$. When $\mathcal{X} = 0$, we use the bias, absolute bias and the root mean squared error.

We compare the proposed fiducial method to the MLE. The simulation results for point estimation are displayed in Tables D.3 and D.4. The numerical results indicate that the fiducial method demonstrates less bias than the MLE across all simulation scenarios. This feature is particularly pronounced at low concentrations, such as 0 and 5, where the relative bias of the MLE is significantly larger than that of the fiducial method (see Tables D.3 and D.4). Furthermore, aside from exhibiting less bias, the fiducial method shows comparable performance with the MLE in terms of both absolute relative bias and relative root mean

Table C.2: Estimated coverage probabilities and expected widths of the 95% confidence intervals for the true unknown concentration: Simulation Scenario 2.B

True Conc.	Average 95% Limits [Average Width]			Coverage Probability		
	Fiducial	Bootstrap	Wald MLE	Fiducial	Bootstrap	Wald MLE
0	(0, 0.9) [0.9]	(0, 1.1) [1.1]	(0, 1.0) [1.0]	0.952	0.849	0.936
5	(4.1, 6.0) [2.1]	(4.4, 6.1) [1.7]	(4.3, 6.0) [1.7]	0.951	0.888	0.921
10	(8.9, 11.2) [2.3]	(9.1, 11.3) [2.2]	(9.0, 11.2) [2.2]	0.952	0.916	0.932
15	(13.6, 16.5) [2.9]	(13.7, 16.4) [2.7]	(13.6, 16.5) [2.9]	0.948	0.918	0.938
20	(18.3, 21.9) [3.6]	(18.2, 21.6) [3.4]	(18.3, 21.8) [3.5]	0.953	0.916	0.941
25	(22.9, 27.2) [4.3]	(22.8, 26.9) [4.1]	(22.8, 27.1) [4.3]	0.953	0.937	0.941
30	(27.5, 32.6) [5.1]	(27.4, 32.2) [4.8]	(27.4, 32.4) [5.0]	0.946	0.932	0.942
35	(32.1, 38.0) [5.9]	(31.7, 37.3) [5.6]	(32.0, 37.7) [5.7]	0.947	0.903	0.940
40	(36.8, 43.5) [6.7]	(36.4, 42.8) [6.4]	(36.6, 43.1) [6.5]	0.954	0.923	0.942
50	(46.0, 54.2) [8.2]	(45.4, 53.3) [7.9]	(45.7, 53.7) [8.0]	0.948	0.912	0.932

squared error.

5 Checking Goodness of Fit of Rocke and Lorenzato Model in Cadmium and Copper Data

Since we are working with real-world datasets, it is essential to first evaluate the adequacy of the Rocke-Lorenzato model specified in the model (1) in the main manuscript before performing the data analysis to estimate the unknown concentrations. To this end, the training data for each dataset were used to fit the Rocke-Lorenzato model (1) in the main manuscript. Subsequently, 95% upper and lower calibration curves were constructed based on the fitted model, as described below.

Table D.3: Comparison of point estimates by the fiducial method and maximum likelihood for simulation Scenarios 1.A and 1.B. Abbreviations: Fid: Fiducial Mode, MLE: Maximum Likelihood Estimator, Rel. Bias: Relative Bias, Rel. Abs. Bias: Relative Absolute Bias, Rel. RMSE: Relative Root Mean Square Error.(For zero concentration Bias, Abs. Bias, RMSE are computed)

Scenario 1.A							Scenario 1.B						
True Conc.	Rel. Bias ($\times 10^2$)		Rel. Abs. Bias ($\times 10^2$)		Rel. RMSE ($\times 10^2$)		True Conc.	Rel. Bias ($\times 10^2$)		Rel. Abs. Bias ($\times 10^2$)		Rel. RMSE ($\times 10^2$)	
	Fid	MLE	Fid	MLE	Fid	MLE		Fid	MLE	Fid	MLE	Fid	MLE
5	-0.22	2.53	9.8	9.7	12	12	0	-1.54	24	46	51	59	63
20	-0.01	0.20	5.0	5.0	6.4	6.3	5	-0.58	4.56	9.0	9.5	11	12
50	0.01	-0.31	4.8	4.7	5.9	5.9	10	0.50	1.46	5.8	5.7	7.2	7.1
							15	-0.21	0.41	4.4	4.3	5.6	5.5
							20	0.10	-0.17	4.3	4.1	5.2	5.2
							25	0.08	-0.24	4.2	4.1	5.3	5.1
							30	0.13	-0.50	4.0	4.0	5.0	5.0
							35	-0.21	0.94	4.3	4.3	5.4	5.3
							40	-0.18	0.86	3.9	3.8	4.9	4.8
							50	0.08	0.91	4.3	4.1	5.3	5.1

Upper and Lower Fitted Calibration Curves Based on Observed data: For every known concentration $\mathcal{X}$ and laboratory i , the fiducial quantity for a random observation following the Rocke-Lorenzato model was computed as $\tilde{\alpha}_i + \tilde{\beta}_i \mathcal{X} \exp(\tilde{\sigma}_\eta z^\eta) + \tilde{\sigma}_\epsilon z^\epsilon$, where $z^\eta, z^\epsilon \sim N(0, 1)$ and $\tilde{\alpha}_i, \tilde{\beta}_i, \tilde{\sigma}_\eta, \tilde{\sigma}_\epsilon$ are the fiducial quantities of $\alpha_i, \beta_i, \sigma_\eta, \sigma_\epsilon$ respectively. Details of the derivation of these fiducial quantities are provided in Section 4 of the main manuscript. Using samples generated from the above fiducial quantity corresponding to a given concentration $\mathcal{X}$, a 95% highest density interval (HDI) was computed. Most of the observed data are expected to fall inside the interval if the data truly followed the Rocke-

Table D.4: Comparison of point estimates by the fiducial method and maximum likelihood for simulation Scenarios 2.A and 2.B. Abbreviations: Fid: Fiducial Mode, MLE: Maximum Likelihood Estimator, Rel. Bias: Relative Bias, Rel. Abs. Bias: Relative Absolute Bias, Rel. RMSE: Relative Root Mean Square Error.(For zero concentration Bias, Abs. Bias, RMSE are computed)

Scenario 2.A							Scenario 2.B						
True Conc.	Rel. Bias ($\times 10^2$)		Rel. Abs. Bias ($\times 10^2$)		Rel. RMSE ($\times 10^2$)		True Conc.	Rel. Bias ($\times 10^2$)		Rel. Abs. Bias ($\times 10^2$)		Rel. RMSE ($\times 10^2$)	
	Fid	MLE	Fid	MLE	Fid	MLE		Fid	MLE	Fid	MLE	Fid	MLE
5	-0.84	1.12	10	10	13	13	0	-0.76	15	34	37	42	45
20	-0.06	0.19	5.3	5.3	6.7	6.7	5	-0.12	3.27	7.5	7.8	9.5	10
50	0.00	-0.10	5.5	5.5	6.9	6.8	10	-0.09	1.27	4.5	4.5	5.7	5.7
							15	-0.22	0.32	4.1	4.0	5.2	5.1
							20	0.06	0.18	3.6	3.6	4.6	4.5
							25	0.00	-0.18	3.5	3.6	4.4	4.3
							30	-0.00	-0.35	3.5	3.4	4.4	4.3
							35	0.05	-0.52	3.5	3.5	4.3	4.3
							40	0.26	-0.31	3.4	3.4	4.2	4.1
							50	0.11	-0.61	3.5	3.4	4.3	4.3

Lorenzato model in (1) in the main manuscript. The upper and lower calibration curves of the fitted model are constructed by joining the upper and lower limits of these 95% intervals, respectively, for every known concentration. For clarification, it is important to distinguish the 95% fiducial HDI so constructed from the 95% fiducial HDI for the unknown concentration discussed in Section 4 in the main manuscript. The 95% fiducial HDI around a known concentration represents the range of random variation in the observed data and is used solely to assess the goodness of fit of the model with the training sample. In contrast, the 95% fiducial HDI for an unknown concentration reflects the uncertainty in

the estimation of that concentration and is determined based on the test sample. These intervals serve distinct purposes and should not be confused.

Figures 1 and 2 display the upper and lower fitted calibration curves for the cadmium and copper datasets, respectively. In both cases, the majority of observed measurements fall within the specified limits, supporting the adequacy of the Rocke-Lorenzato model (1) in the main manuscript for these datasets.

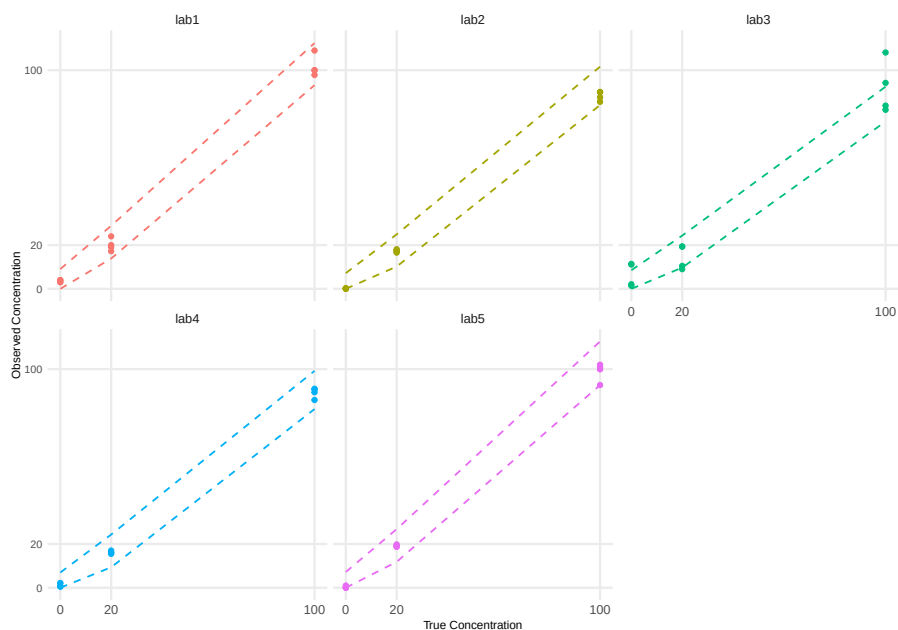


Figure 1: Goodness of fit of cadmium data: Lower and Upper Fitted Calibration Curves (see details in Section 6 of the main manuscript) for each laboratory (dashed lines) and observed measurements (points)

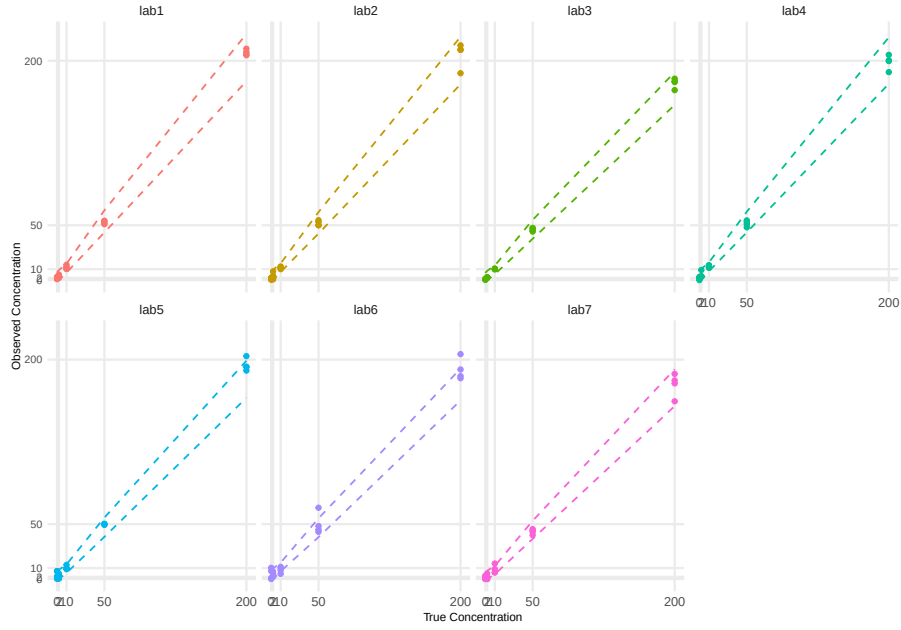


Figure 2: Goodness of fit of copper data: Lower and Upper Fitted Calibration Curves (see details in Section 6 of the main manuscript) for each laboratory (dashed lines) and observed measurements (points)

6 Plots Regarding Interval Estimation of Unknown Concentration for Cadmium and Copper Data

Figures 3 and 4 provide a visual summary of the 95% confidence intervals alongside the upper and lower calibration curves for the unknown concentrations in both cadmium and copper datasets. The upper and lower calibration curves are constructed by joining the upper and lower limits of 95% intervals, respectively, for every unknown concentration.

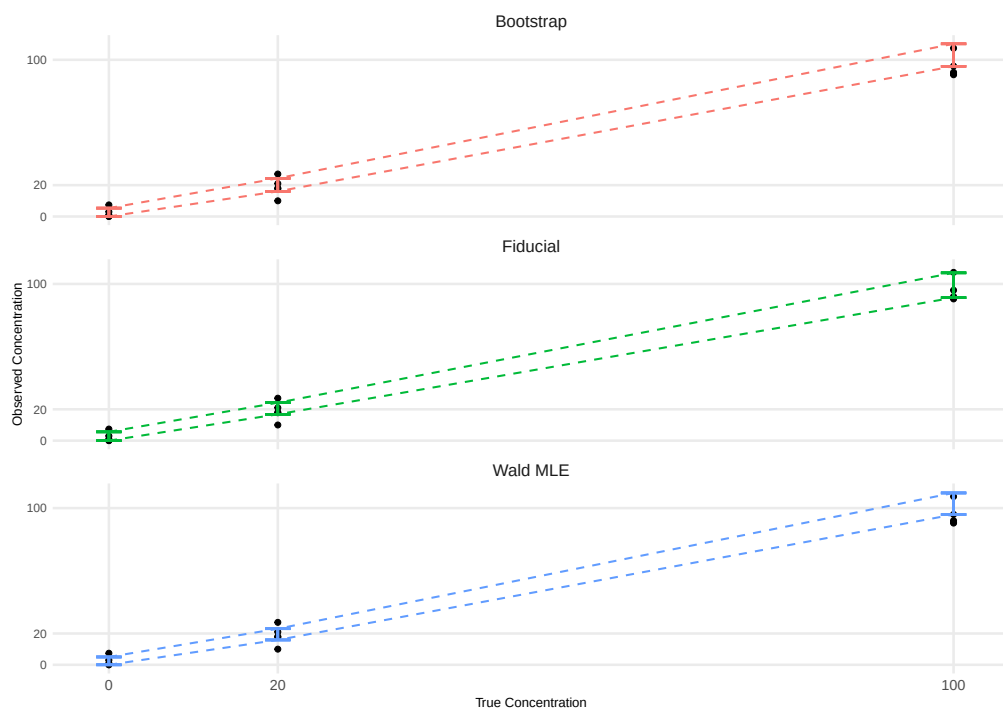


Figure 3: Interval estimation for cadmium data: 95% Confidence intervals of unknown concentrations (True values: (0,20,100)) by bootstrap, fiducial, and the MLE-based Wald method

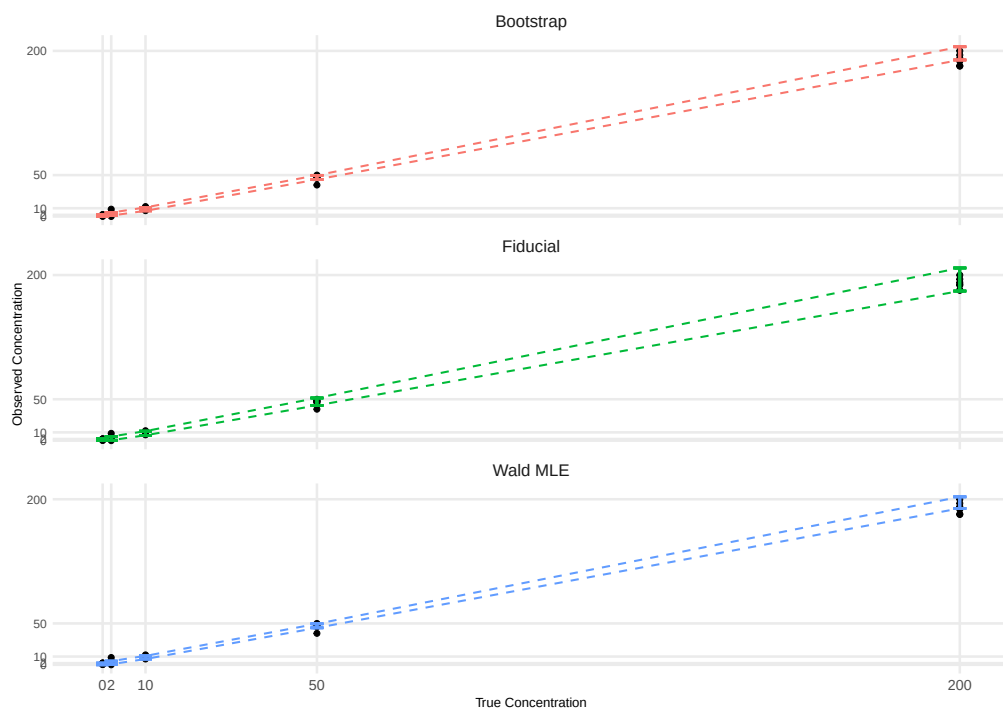


Figure 4: Interval estimation for copper data: 95% Confidence intervals of unknown concentrations (True values: (0,2,10,50,200)) by bootstrap, fiducial, and the MLE-based Wald method.