
GEM: Empowering MLLM for Grounded ECG Understanding with Time Series and Images

Xiang Lan¹ Feng Wu¹ Kai He¹ Qinghao Zhao² Shenda Hong^{3,4} Mengling Feng¹

Abstract

While recent multimodal large language models (MLLMs) have advanced automated ECG interpretation, they still face two key limitations: (1) insufficient multimodal synergy between time series signals and visual ECG representations, and (2) limited explainability in linking diagnoses to granular waveform evidence. We introduce GEM, the first MLLM unifying ECG time series, 12-lead ECG images and text for grounded and clinician-aligned ECG interpretation. GEM enables feature-grounded analysis, evidence-driven reasoning, and a clinician-like diagnostic process through three core innovations: a dual-encoder framework extracting complementary time series and image features, cross-modal alignment for effective multimodal understanding, and knowledge-guided instruction generation for generating high-granularity grounding data (ECG-Grounding) linking diagnoses to measurable parameters (*e.g.*, QRS/PR Intervals). Additionally, we propose the Grounded ECG Understanding task, a clinically motivated benchmark designed to comprehensively assess the MLLM’s capability in grounded ECG understanding. Experimental results on both existing and our proposed benchmarks show GEM significantly improves predictive performance (CSN 7.4% \uparrow), explainability (22.7% \uparrow), and grounding (24.8% \uparrow), making it more suitable for real-world clinical applications. GitHub repository: <https://github.com/lanxiang1017/GEM.git>

1. Introduction

Electrocardiography (ECG), a cornerstone of cardiac diagnostics, captures the heart’s electrical activity through body-surface electrodes, enabling non-invasive assessment of cardiac physiology and pathology (Berkaya et al., 2018; Hannun et al., 2019). Clinical ECG interpretation synergizes computational and clinical expertise: automated algorithms process raw signals to detect fiducial points (*e.g.*, R-wave peaks) and generate diagnostic hypotheses, while clinicians validate these findings through 12-lead waveform analysis (Smulyan, 2019). By contextualizing algorithmic outputs with patient-specific factors, clinicians resolve ambiguities, detect subtle anomalies, and formulate diagnoses. This synergy between computational precision and expert judgment ensures reliable and holistic diagnoses in clinical practice.

Deep learning models have achieved promising results in cardiac anomalies detection (Hong et al., 2020; Zhu et al., 2021; Kiyasseh et al., 2021; Lan et al., 2022; Li et al., 2024; Lan et al., 2024; McKeen et al., 2024) yet lack conversational capability and model explainability. While recent MLLMs like PULSE (Liu et al., 2024c) have advanced language-based ECG interpretation via large-scale instruction tuning (*e.g.*, ECG-Instruct’s 1M+ samples), they primarily focus on static image inputs and predefined diagnostic tasks, leaving two critical gaps: **Limited Modality Exploration**. Current models typically process only one non-text modality, overlooking the benefits of a synergistic approach that combining time series and image-based modeling. For example, time series models capture dynamic changes but may overlook spatial patterns, while image-based models detect global structures but may miss subtle temporal details. This limits their ability to replicate clinicians’ holistic reasoning, which integrates both machine-measured temporal signals and waveform patterns from 12-lead plots. **Insufficient Explainability and Grounding**. Most existing models designed for ECG provides limited explainability, failing to explicitly connect diagnoses to granular waveform evidence and provide insight into their diagnostic reasoning. A trustworthy ECG model should not only predict cardiac conditions but also explicitly highlight which ECG features led to those conclusions. Such grounded explanations enhance transparency and make model outputs

¹Saw Swee Hock School of Public Health and Institute of Data Science, National University of Singapore, Singapore ²Department of Cardiology, Peking University People’s Hospital, Beijing, China ³National Institute of Health Data Science, Peking University, Beijing, China ⁴Institute for Artificial Intelligence, Peking University, Beijing, China. Correspondence to: Shenda Hong <hongshenda@pku.edu.cn>, Mengling Feng <ephfm@nus.edu.sg>.

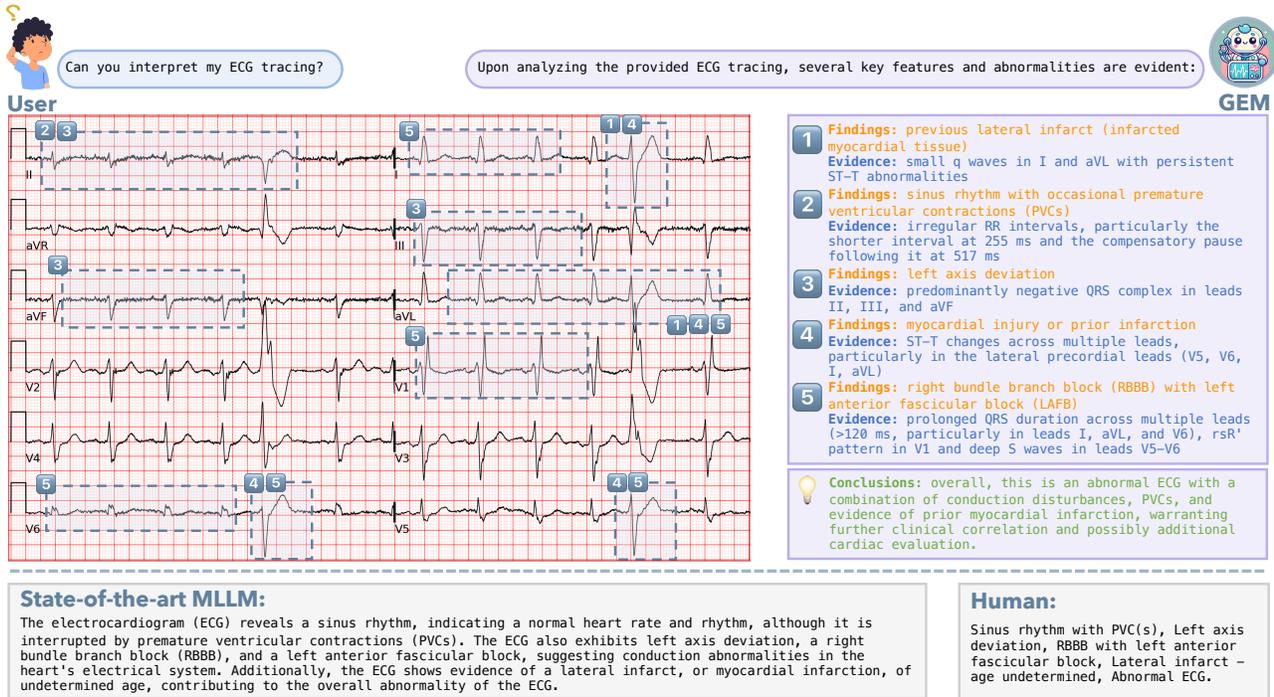


Figure 1. GEM offers superior granularity in ECG interpretation compared to state-of-the-art models and human-written reports. GEM’s core capabilities: 1.Feature-Grounded Analysis: findings are precisely linked to measurable ECG features (e.g., QRS/RR intervals). 2.Evidence-Driven Diagnosis: conclusions are supported by clear and logical reasoning directly linked to ECG findings. 3.Realistic Interpretation Process: mimics how a clinician analyzes ECGs and arrive at a diagnosis.

more reliable and actionable for clinical decision-making.

To address these challenges, we introduce GEM, a multi-modal large language model designed for grounded ECG understanding by integrating time series and image data. As a conversational cardiology AI assistant, GEM differs from other MLLMs through three key features. First, it provides **feature-grounded analysis**, ensuring that its findings are explicitly tied to detailed ECG features. Second, it offers **evidence-driven diagnosis**, where its conclusions are supported by clear and logical reasoning directly linked to ECG findings. Lastly, GEM simulates a **realistic interpretation process**, mimicking how a clinician analyzes ECGs and arrive at a diagnosis.

These features are enabled by three major components of GEM. **Multimodal Encoding:** GEM adopts a dual-encoder architecture with one encoder for raw ECG time series and another for their transformed images. Each encoder specialized in extracting modality-specific features using established models from the time series and vision domains, allowing GEM to effectively leverage the strengths from both modalities. **Cross-modal Alignment Learning:** To enable the LLM to interpret ECG time series and images, GEM first projects time series representations to the image representations dimensionality, then applies a shared

projector to transform both into language-like embeddings that are comprehensible to LLM. These aligned representations are fused with textual instruction embeddings for next-token prediction training, enabling the LLM to effectively understand multimodal inputs for interactive ECG interpretation. **Knowledge-Guided Instruction Generation:** GEM is trained on high-granularity instruction data annotated with heartbeat-level physiological features. These data are generated using a knowledge-guided instruction strategy that incorporates grounding feature extraction and a diagnosis guider. Grounding feature extraction derives precise heartbeat-level physiological features from ECG time series, ensuring explicit grounding for ECG analysis. The cardiology-specific diagnosis guider processes these features into structured prompts, leveraging GPT-4o’s latent medical knowledge to generate clinically detailed and feature-grounded ECG instruction data.

The main contributions of this work are three-folds:

1.First Unified Multimodal ECG Model. We introduce GEM, the first multimodal framework to synergistically integrate raw ECG time series signals, 12-lead plots, and textual instructions, leveraging their complementary strengths to advance grounded ECG understanding.

2. *First High-granularity ECG Grounding Dataset.* We curate ECG-Grounding, a novel dataset comprising 30,000 instruction pairs annotated with heartbeat-level physiological features. This is the first high-granularity ECG grounding dataset, enabling evidence-based diagnosis and improving the trustworthiness of medical AI.

3. *Clinically Oriented Diagnostic System.* We propose the Grounded ECG Understanding task, a clinically motivated benchmark designed to comprehensively assess a model’s ECG interpretation capability. Experimental results demonstrate that GEM not only excels in predictive performance but also significantly enhances explainability and grounding, making it more applicable for real-world clinical settings while fostering greater trust among medical professionals.

2. Related Work

2.1. Multimodal Large Language Models

Large Language Models (LLMs), such as GPTs (Brown et al., 2020; Achiam et al., 2023), LLaMA (Touvron et al., 2023), have made significant advancements in artificial intelligence. Despite their superior performance on numerous natural language processing tasks, LLMs are inherently limited to the text modality, making them “blind” to other modalities such as images, audio and video. To mitigate this constraint, Multimodal Large Language Models (Liu et al., 2024b;a; Zhu et al., 2023; Gao et al., 2023; Gong et al., 2023; Li et al., 2023b; Pi et al., 2023; Lyu et al., 2023; Chen et al., 2023) have been recently developed to extend the ability of LLMs in comprehending multiple modalities. By integrating LLMs with various data sources, MLLMs enable the handling of diverse information beyond text. For example, LLaVA (Liu et al., 2024b) enables LLMs to comprehend visual inputs by adopting a learnable projector to map image features into the word embedding space. Video-LLaMA (Zhang et al., 2023) further enhances LLMs by enabling video perception and understanding. Qwen-Audio (Chu et al., 2023) introduces an audio-language model capable of processing various audio types, including human speech, natural sounds, and music. Medical data, by nature, are inherently multimodal, encompassing diverse formats such as X-ray images, physiological time series (*e.g.*, ECG signals), and textual reports. These modalities collectively form the foundation for clinical decision-making, driving the advancement of MLLMs in medical applications (He et al., 2025; Li et al., 2023a; Thirunavukarasu et al., 2023; Kim et al., 2024; Moor et al., 2023; Liu et al., 2023; Fu et al., 2024; Zhu et al., 2024). Different from these works, our focus is on empowering MLLMs with the grounded ECG understanding capability.

2.2. Language-based ECG Analysis

Language-based ECG diagnosis and interpretation is still in its early stages of development. Only a few recent studies have explored LLM-based approaches for ECG analysis. For instance, Yu et al. (2023) proposed a zero-shot retrieval-augmented diagnosis technique, embedding domain knowledge from textbooks and research papers into a vector database to improve zero-shot diagnostic accuracy. Liu et al. (2024c) introduced PULSE, an LLM-based framework designed to enhance ECG image understanding for diagnosis and report generation. PULSE synthesizes realistic ECG images from raw ECG signals, enabling better utilization of image-based LLaVA models. Zhao et al. (2024) developed ECG-CoCa, an ECG encoder trained on ECG-text pairs, alongside ECG-Chat, a modified LLaVA model capable of processing ECG time series. Chan et al. (2024) proposed an analytical framework integrating time series data with LLMs, combining physiological signal analysis with contextual textual information. Wan et al. (2024) designed an instruction-tuning framework for automated ECG report generation, converting ECG-text pairs into chatbot-style instructions and fine-tuning the LLM’s linear layers.

GEM distinguishes itself from existing language-based ECG models two key aspects. First, in model architecture, unlike existing methods that limit analysis to isolated modalities (either ECG signals or transformed ECG images) with text, GEM introduces a unified architecture that integrates both time series data and 12-lead images. This approach mirrors a clinician’s natural workflow, where dynamic signal trends and spatial waveform patterns are jointly analyzed for a more comprehensive interpretation. Second, in ECG understanding, GEM establishes a new paradigm for evidence-driven diagnosis. While current models often lack grounded understanding, GEM enables heartbeat-level interpretability by directly linking each diagnostic conclusion to quantifiable physiological evidence, enhancing explainability and clinical reliability. By combining the two paradigm-shifting innovations, GEM addresses fundamental limitations in existing models and advances language-based ECG analysis, setting a new standard for AI-assisted cardiac diagnostics.

3. Method

3.1. Architecture

GEM mainly comprises three modules: the multimodal encoding module, the cross-modal alignment learning module, and the knowledge-guided instruction generation module. Figure 2 provides an overview of our proposed GEM model. We will elaborate on each module in the following sections.

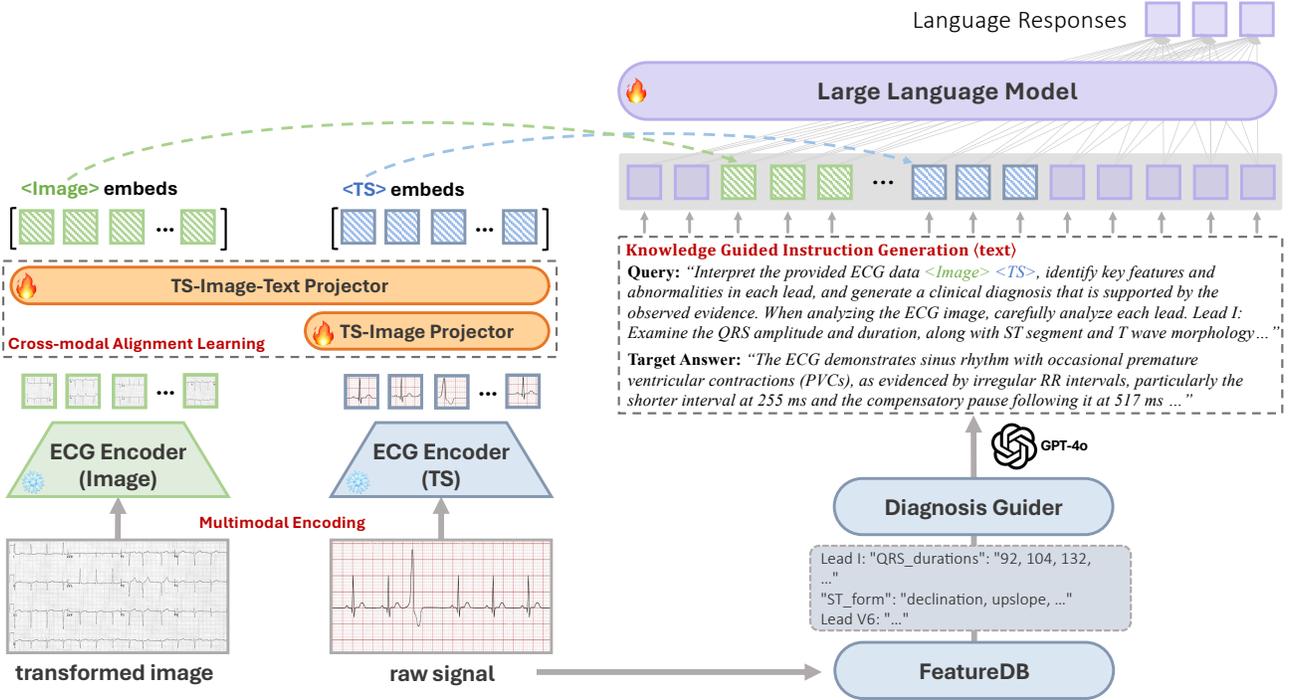


Figure 2. GEM’s Architecture. *Multimodal Encoding*: Separate encoders process ECG time series and images to generate modality-specific representations, enabling a holistic analysis of ECG data. *Cross-modal Alignment Learning*: Time series and image representations are first aligned and then mapped to a textual space using a shared projector, ensuring coherent understanding for the LLM. *Knowledge Guided Instruction Generation*: Physiological features extracted from all 12 leads are sequenced and structured using a diagnosis guider, which prompts GPT-4o with domain-specific instructions to generate high-granularity instructional data.

3.2. Multimodal Encoding

We represent the ECG time series as $\mathbf{x}_{ts} \in \mathbb{R}^{C \times L}$, where C is the number of leads in multi-lead ECG and L is the length of the signal. The transformed ECG image derived from the time series is denoted as $\mathbf{x}_{img} \in \mathbb{R}^{H \times W \times 3}$, where H and W are the height and width of the image, respectively. For ECG time series encoder $E_{ts}(\theta_{ts})$, we adopt the pre-trained ECG-CoCa model (Zhao et al., 2024), which has been extensively trained on a large number of ECG-Text pairs with contrastive learning to effectively capture the intricate patterns within the ECG time series:

$$\mathbf{e}_{ts} \in \mathbb{R}^{n_s \times d_s} = E_{ts}(\mathbf{x}_{ts} | \theta_{E_{ts}}), \quad (1)$$

where n_s is the number of time series feature patches and d_s is the dimension of time series features. For ECG image encoder $E_{img}(\theta_{img})$, we utilize the pre-trained CLIP encoder from LLaVA (Liu et al., 2024b). This model is adept at understanding and processing visual information, making it suitable for extracting features from the ECG images:

$$\mathbf{e}_{img} \in \mathbb{R}^{n_m \times d_m} = E_{img}(\mathbf{x}_{img} | \theta_{E_{img}}), \quad (2)$$

where n_m is the number of image feature patches and d_m is the dimension of image features. These two encoders enable

the separate extraction of features from time series and images. This dual-encoder approach allows us to harness the distinct advantages of each data type, enhancing the overall interpretative power of our model.

3.3. Cross-modal Alignment Learning

Considering that the time series and image encoders are trained independently, their generated representations often show inconsistencies within the representation space. Meanwhile, to ensure that the LLM can interpret ECG time series and images effectively, it is essential that these multimodal inputs are rendered as comprehensible as textual data. Therefore, aligning these diverse modality inputs within a unified representation space becomes crucial.

In our approach, the ECG time series representation is first mapped to the same dimensionality as the ECG image representation. This is accomplished using a multi-layer perceptron (MLP) projector:

$$\hat{\mathbf{e}}_{ts} \in \mathbb{R}^{n_s \times d_m} = MLP_{ts}(\mathbf{e}_{ts} | \theta_{M_{ts}}). \quad (3)$$

Subsequently, we employ an additional projector to map both the time series and image representations into a consis-

tent textual space:

$$\mathbf{h}_{ts} \in \mathbb{R}^{n_s \times d_t} = MLP(\hat{\mathbf{e}}_{ts} | \theta_M), \quad (4)$$

$$\mathbf{h}_{img} \in \mathbb{R}^{n_m \times d_t} = MLP(\mathbf{e}_{img} | \theta_M), \quad (5)$$

where d_t is the dimension of the text embeddings. This step is for ensuring that the features extracted from both modalities are not only aligned dimensionally but are also interpretable in a LLM-friendly format.

Once we have obtained the features from both the time series and image modalities, we integrate these with the embeddings of the textual query \mathbf{x}_q :

$$\mathbf{x} = \text{Concatenate}(\mathbf{h}_{ts}, \mathbf{h}_{img}, \text{Embed}(\mathbf{x}_q)). \quad (6)$$

The integration is crucial for creating a cohesive representation that encapsulates the full spectrum of information from the multimodal inputs.

3.4. Knowledge-Guided Instruction Generation

Instruction data forms the foundation of multimodal training, directly shaping how the MLLM generates responses for given queries. This is because the language response $\theta_{LLM}(\mathbf{x})$ is optimized to match the target answer \mathbf{y} for each multimodal input $(\mathbf{x}_{ts}, \mathbf{x}_{img}, \mathbf{x}_q)$. Here, we introduce two key mechanisms to guarantee that \mathbf{y} aligns with feature-grounded analysis, evidence-driven diagnosis, and a realistic interpretation process, thereby empowering GEM with a grounded understanding of ECG data.

Grounding Features Extraction. To enable the feature-grounded analysis, we propose to further excavate more detailed information from raw ECG time series. We begin by extracting universal elements, including waveforms, amplitudes of fiducial points, and intervals, from each heartbeat in each lead. These elements are then structured into a time-ordered sequence for further analysis. For instance, an ECG time series with ten visible heartbeats allows us to construct a QRS duration sequence as $[QRS_1, QRS_2, \dots, QRS_{10}]$, where each element precisely represents the QRS duration for a specific heartbeat. These sequences provide a fine grade description of the physiological features of the ECG, enabling the model to assess heart conditions at the level of individual heartbeats. In our implementation, we incorporate 14 feature sequences for each of the 12 ECG leads, covering: Heart Rate, RR Interval 1, RR Interval 2, P Amplitude, P Duration, PR Interval, QRS Amplitude, QRS Duration, T Amplitude, T Duration, ST Duration, ST Form, QT Interval, and QTc Interval. This comprehensive set of sequences captures the temporal evolution of key physiological features, enabling granular analysis of cardiac activity. The feature extraction process is formulated as:

$$\mathbf{x}_{fs} = \text{FeatureDB}(\mathbf{x}_{ts}), \quad (7)$$

where \mathbf{x}_{fs} is a dictionary in which the keys represent feature names and the values correspond to their respective feature values. $\text{FeatureDB}(\mathbf{x}_{ts})$ represents the function that extracts structured physiological features from the ECG time series \mathbf{x}_{ts} . Note that there is no trainable parameters in FeatureDB.

Diagnosis Guider. With the feature sequences \mathbf{x}_{fs} extracted, the challenge lies in generating high-granularity \mathbf{y} without relying on costly human-expert annotation. To address this, we design a diagnosis guider that constructs a prompt \mathbf{x}_p to effectively guide GPT-4o in generating accurate and clinically grounded responses \mathbf{y} for each sample:

$$\mathbf{x}_p = \text{DiagnosisGuider}(\mathbf{x}_{fs}). \quad (8)$$

The diagnosis guider provides cardiology-specific instructions for detailed analysis of each aspects of ECG data (*e.g.*, instruct GPT-4o to assess the P wave amplitude and duration in Lead II to evaluate atrial enlargement) and incorporates guidance reflecting real-world clinical diagnostic processes. As a result, each sample receives a unique \mathbf{x}_p tailored to its \mathbf{x}_{fs} , ensuring precise activation of GPT-4o’s latent medical knowledge for accurate and personalized analysis. See Appendix A for more details of the diagnosis guider.

ECG-Grounding Data. Using the knowledge-guided instruction generation method, we employ GPT-4o to curate approximately 30,000 fine-grained instruction-response pair $(\mathbf{x}_q, \mathbf{y})$ from the MIMIV-IV-ECG (Gow et al., 2023) database:

$$\mathbf{y} = \text{GPT-4o}(\mathbf{x}_p). \quad (9)$$

A comparison illustrating the advantages of our instruction data is shown in Figure 3. Instruction data from PULSE, primarily derived from the original report, occasionally provides incorrect explanations due to hallucinations, which can arise when LLM has limited reference information. In contrast, our ECG-Grounding provides more accurate, holistic, and evidence-driven interpretations with diagnoses grounded in measurable ECG features. Training with these fine-grained instructional data, GEM significantly enhances its explainability and grounding capabilities, making it more suitable for clinical applications and fostering greater trust among medical professionals.

3.5. Training

Unlike conventional multi-step training pipelines for most MLLMs, which first train a linear projector using brief image captions for cross-modal alignment and subsequently fine-tune the LLM with instruction data, we adopt an one-step training approach for GEM. In this approach, we freeze θ_{ts} and θ_{img} for feature extraction, while jointly training $\theta_{M_{ts}}$, θ_M , and θ_{LLM} . This one-step process enhances the consistency in training multiple modalities with limited data and improves the efficiency of the training phase.

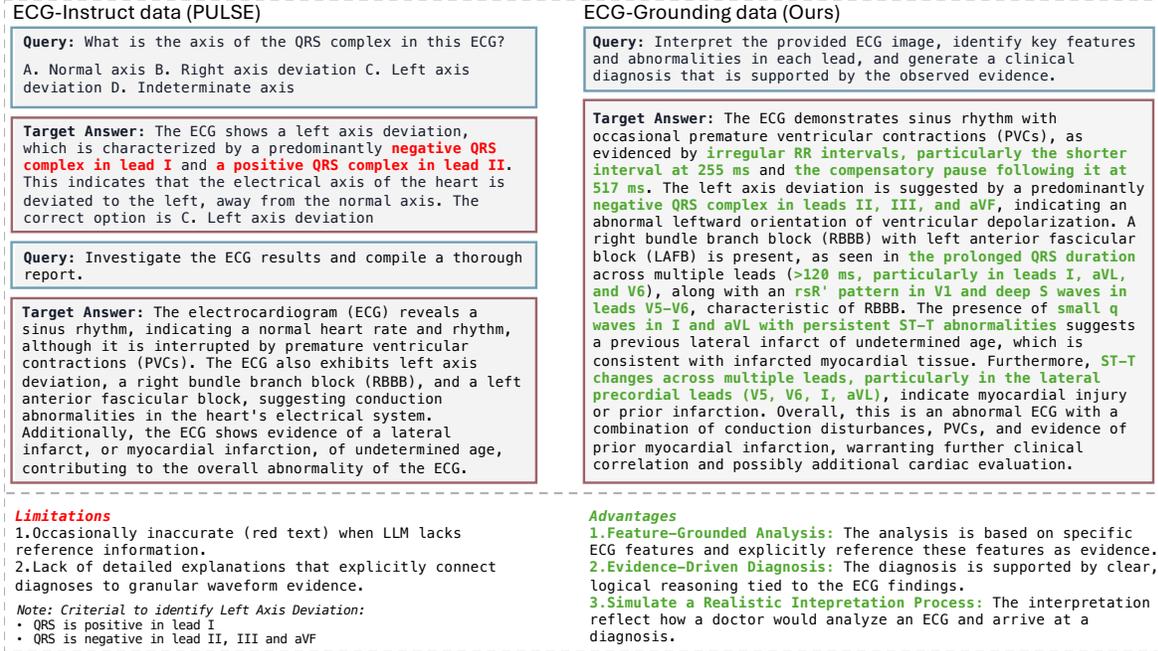


Figure 3. Comparison of ECG-Instruct and our ECG-Grounding.

The training objective is formulated as minimizing the negative log-likelihood of the target answer y and the LLM response $\theta_{LLM}(\mathbf{x})$, given the embedding of multimodal input \mathbf{x} :

$$L(\mathbf{y}; \theta_{LLM}(\mathbf{x})) = - \sum_{i=1}^N \log P(y_i | \mathbf{x}, \theta_{LLM}), \quad (10)$$

where N represents the number of tokens in \mathbf{y} and y_j is the j -th token in \mathbf{y} . We adopt two foundational LLMs for the GEM framework. The first, LLaVA-7B, is the unmodified version of the original LLaVA model, which has not undergone training on any ECG data. The second, referred to as PULSE-7B, represents state-of-the-art LLMs specifically trained on millions of ECG images. For both models, we implement supervised fine-tuning (SFT) for a single epoch. We use 8 NVIDIA A100 GPUs for the training process.

4. Experiments

4.1. Training Data

The training of GEM involves two datasets: the ECG-Instruct data from PULSE (Liu et al., 2024c), which includes 1,156,110 conversations, and our own ECG-Grounding data, comprising 30,000 conversations. The ECG-Grounding data is sampled from MIMIC-IV-ECG, selecting only samples that have not been used in training models such as PULSE. We use the ECG-image-kit (Shivashankara et al., 2024) for the generation of ECG images from the original ECG signal,

and use the FeatureDB (Hong et al., 2017; 2019) to extract ECG features. To facilitate further research and benefit the community, we will publicly release ECG-Grounding data.

4.2. Evaluation Tasks and Metrics

Grounded ECG Understanding. To comprehensively evaluate whether the model achieves clinically grounded ECG interpretation capabilities comparable to real cardiologists, we introduce the Grounded ECG Understanding task. This task evaluates the MLLM’s ability to identify detailed clues in ECG analysis, requiring it to provide specific details and relevant domain knowledge to support its interpretation.

We utilize GPT-4o to score the responses of the MLLM using a predefined set of metrics that measures the accuracy and comprehensiveness of the details provided. Specifically, these metrics include: *DiagnosisAccuracy* evaluates whether the generated diagnosis is correct, specific, and supported by ECG findings. Results are expressed as a percentage, indicating the average accuracy across identified key diagnoses. *AnalysisCompleteness* checks if all key ECG components (e.g., rhythm, intervals, waveforms, and lead-specific findings) are discussed. Results are provided in absolute terms, indicating the average number of correctly addressed key ECG features for each sample. *AnalysisRelevance* assesses whether each explanation directly supports the diagnosis, with results showing on average how many points support the diagnosis with clear ECG evidence for each sample. *LeadAssessmentCoverage* evaluates

Datasets	MIMIC-IV-ECG				
	Diagnosis Accuracy	Analysis Completeness	Analysis Relevance	LeadAssessment Coverage	LeadAssessment Accuracy
PULSE	75.94	2.37	2.39	5.84	2.00
GEM (Ours)					
SFT LLaVA-7B	86.34	4.41	5.01	51.60	33.27
SFT PULSE-7B	86.21	4.43	4.91	51.60	33.07

Table 1. Grounded ECG Understanding results (Part 1).

Datasets	MIMIC-IV-ECG			
	ECG Feature Grounding	Evidence-Based Reasoning	Clinical Diagnostic Fidelity	Average
PULSE	50.18	52.40	51.63	51.40
GEM (Ours)				
SFT LLaVA-7B	74.48	75.09	75.28	75.28
SFT PULSE-7B	74.95	74.70	74.87	74.84

Table 2. Grounded ECG Understanding results (Part 2).

how many of the 12 ECG leads are analyzed. Results indicate the average number of leads analyzed per sample, providing insight into the comprehensiveness of the ECG review. *LeadAssessmentAccuracy* verifies the accuracy of described lead findings (*e.g.*, QRS, ST, T waves, amplitude, intervals, ST segments) against the ground truth interpretation. *ECGFeatureGrounding* determines if the interpretation references actual ECG features (*e.g.*, QRS amplitude, PR interval) instead of generic terms. Results are scaled from 0 to 100. *EvidenceBasedReasoning* evaluates whether the diagnosis follows logical, evidence-supported steps. Results range from 0 to 100. and *ClinicalDiagnosticFidelity* assesses if the model mimics how a clinician interprets ECG data, considering all relevant factors. Results are scaled from 0 to 100. Further details of the scoring method for each criteria can be found in the Appendix B. The entire test dataset comprises 2,381 samples from the MIMIC-IV-ECG database that were not included in the training set.

ECG-Bench. In addition, we use the ECG-Bench (Liu et al., 2024c) to assess our model’s capability in cardiac abnormality detection and report generation. ECG-Bench is a comprehensive benchmark for evaluating MLLMs on ECG image interpretation. It incorporates several datasets including the PTB and PTB-XL dataset (Wagner et al., 2020), CPSC2018 dataset (Liu et al., 2018), G12EC dataset, CODE-15% dataset (Ribeiro et al., 2021), and CSN dataset (Zheng et al., 2020). ECG-Bench contains two main tasks: abnormality detection and report generation. For the abnormality detection task, it uses AUC, F1, and Hamming Loss (HL) as metrics for multi-label datasets, and accuracy for others.

In the report generation task, it employs GPT-4o to evaluate the reports based on their accuracy in rhythms, waveform descriptions, and diagnoses, with a maximum score of 100.

4.3. Results

In the Grounded ECG Understanding task, Table 1 shows the robust performance of the GEM models. The proposed GEM models, including both SFT LLaVA-7B and SFT PULSE-7B, demonstrate notable advancements in diagnosis accuracy, achieving scores over 86% compared to 75.94% by PULSE. This highlights GEM’s superior capability to accurately identify and interpret complex ECG data. In terms of the analysis completeness and relevance, the results reveal significant performance advantages of GEM models over PULSE. Specifically, the GEM models achieve an 87% improvement in analysis completeness and a 110% enhancement in analysis relevance compared to the baseline. This substantial performance gap indicates that GEM not only identifies and interprets a greater number of critical ECG components (with an average two-fold increase in feature coverage) but also maintains stronger clinical relevance by effectively connecting these observations to the diagnostic reasoning process. Moreover, the substantial improvements in lead assessment coverage, increasing from 5.84 to over 51, and in accuracy, rising from 2 to over 33, underscore a thorough analysis across multiple ECG leads by GEM. This enhancement reflects the effectiveness of our ECG-Grounding data in providing a more comprehensive and precise evaluation of the cardiac activity captured in ECG recordings.

Datasets	PTB-XL Super			CODE-15%			CPSC 2018			CSN	G12EC
Metric	AUC	F1	HL	AUC	F1	HL	AUC	F1	HL	Accuracy	Accuracy
Random	50.3	33.2	50.1	48.8	15.0	32.1	51.2	15.1	28.8	11.6	12.1
GPT-4o	55.6	28.3	26.2	59.9	24.9	15.7	50.9	10.6	18.2	57.5	49.2
PULSE	82.4	74.8	11.0	90.7	85.4	5.0	76.9	57.6	8.6	85.2	78.2
GEM (Ours)											
SFT LLaVA-7B	81.8	73.6	11.6	90.5	84.8	5.1	74.1	52.0	9.0	92.6	81.8
SFT PULSE-7B	83.4	75.8	11.0	91.5	86.4	4.7	79.1	61.1	8.1	86.2	80.5
Ablations											
TS only	81.2	72.5	11.9	90.8	84.9	5.0	76.3	54.0	8.5	91.6	81.4
TS+IMG	82.7	74.8	11.1	91.3	86.3	4.6	74.4	51.5	8.8	90.1	81.1

Table 3. ECG-Bench abnormality detection results.

Datasets	PTB-XL Report	ECG-QA
Metric	Report Score	Accuracy
Random	0	16.2
GPT-4o	50.2	35.2
PULSE	61.3	73.8
GEM (Ours)		
SFT LLaVA-7B	65.0	71.0
SFT PULSE-7B	67.1	73.6

Table 4. ECG-Bench report generation and QA results.

Table 2 expands the performance evaluation by incorporating ECG feature grounding and finer aspects of clinical reasoning, including evidence-based reasoning and clinical diagnostic fidelity, both essential for real-world clinical decision-making. ECG feature grounding assesses the extent to which interpretations reference specific ECG features (*e.g.*, QRS amplitude, PR interval) rather than relying on generic descriptions. A perfect score of 100 in this metric would indicate that all ECG findings are comprehensively documented, clearly linked to diagnoses, and cover all pertinent ECG features. In this metric, GEM achieves a score near 75, demonstrating that most ECG findings are explicitly cited and connected to diagnostic conclusions. Furthermore, the metrics for evidence-based reasoning and clinical diagnostic fidelity yield impressive scores above 75. These results demonstrate that the GEM’s findings logically progress toward diagnoses with clear clinical justifications and that its analytical approach aligns closely with structured clinical workflows. Such capabilities are essential for fostering trust and ensuring practical utility in clinical environments.

In the ECG-Bench task, Table 3 offers a detailed examination of the models’ performance on cardiac abnormality detection across a variety of datasets. The results clearly show that our GEM (SFT PULSE-7B) model consistently outper-

forms the state-of-the-art PULSE model. Furthermore, the GEM (SFT LLaVA-7B) version, which has not been previously trained on ECG data, still manages to achieve comparable, and in some cases superior (*e.g.*, 7.4% in CSN and 3.6% in G12EC), performance across most datasets, despite being trained for only a single epoch. These results highlight the robustness and efficacy of the GEM framework, demonstrating its capacity to deliver substantial performance gains even with minimal domain-specific training. Ablation studies further highlight the critical role of multimodal inputs in achieving efficient training and superior classification performance. We assess two model variants: TS-only, which is trained exclusively on ECG time series, and TS+IMG, which incorporates both time series and image data. Notably, the TS+IMG model surpasses PULSE on PTB-XL Super and CODE-15% datasets, despite being trained for only one epoch. In Table 4, GEM showcases exceptional report generation capabilities, achieving a substantial 5.8% improvement over PULSE in the PTB-XL Report while maintaining comparable performance in ECG-QA. These results further evident GEM’s capability in delivering holistic, accurate, and clinically-aligned ECG interpretations.

Collectively, GEM demonstrates superior performance in both Grounded ECG Understanding and ECG Bench tasks, establishing its effectiveness in ECG interpretation across multiple dimensions.

4.4. Cardiologist Evaluation

We conducted cardiologist evaluations to further assess GEM’s applicability in real-world clinical scenarios. Specifically, cardiologists review GEM generated interpretations along with corresponding target answers produced by GPT-4o. We ask them to highlight two types of findings: (1) those exceeding their expectations for a clinical AI assistant in real-world settings, and (2) those where they hold differing opinions. We present six representative examples involving

complex cardiac conditions, as shown in Appendix C. Overall, GEM demonstrated its capability to generate clinically insightful findings, often surpassing expert expectations, suggesting its potential for real-world clinical applications.

It is noteworthy that cardiologists also highlights certain cases where their interpretations diverged from GEM or GPT-4o. Although our knowledge-guided instruction generation approach avoids costly expert annotations while producing high-quality target answers, GPT-4o still occasionally generates target answers misaligned with cardiologist interpretations. For example, in Figure 5, the GPT-4o suggest no evidence ischemia or infarction, while cardiologist suspects there are ischemia in the precordial leads. These deviations highlight opportunities for future refinement of model reasoning to better align with cardiologist level clinical judgment.

5. Conclusion

In this work, we present GEM, the first MLLM for grounded ECG interpretation, integrating ECG time series, 12-lead ECG images, and textual instructions. GEM achieves feature-grounded analysis, evidence-driven diagnosis, and clinician-style diagnostic workflows through three core technical innovations: (1) multimodal encoding, (2) cross-modal alignment learning, and (3) knowledge-guided instruction generation. To address the granularity limitations of existing instruction data for ECG understanding, we develop ECG-Grounding, a novel dataset comprising 30,000 fine-grained instruction pairs annotated with heartbeat-level physiological features. This establishes the first high-resolution benchmark for grounded ECG understanding. Furthermore, we introduce the Grounded ECG Understanding task, a clinically motivated evaluation benchmark that comprehensively assesses diagnostic fidelity, interpretative completeness, and grounded ECG understanding through multi-dimensional metrics. Experimental results show that GEM significantly improves predictive performance, explainability, and grounding, outperforming state-of-the-art models even with limited training epochs.

Future works. To further improve GEM, we plan the following enhancements: First, we plan to expand the ECG-Grounding data to include diverse demographic and pathological cases to enhance generalizability. We will gradually release the updated version of our ECG-Grounding data. Second, we plan to incorporate real-world clinical feedback from cardiologists to further refine the model’s diagnostic reasoning capability. Third, we aim to further enhance multimodal data alignment to improve representation learning, ensuring that extracted features capture more fine-grained, clinically relevant information closely aligned with the textual instructions.

By bridging computational precision with clinician-level reasoning, GEM represents a step toward more reliable, explainable, and clinically applicable AI-assisted ECG analysis, with potential for broader integration into real-world healthcare workflows.

Acknowledgment

This research is supported by A*STAR, CISCO Systems (USA) Pte. Ltd, and National University of Singapore under its Cisco-NUS Accelerated Digital Economy Corporate Laboratory (Award I21001E0002) and the National Research Foundation, Singapore under its AI Singapore Programme (AISG Award No: AISG2-TC-2022-004).

References

- Achiam, J., Adler, S., Agarwal, S., Ahmad, L., Akkaya, I., Aleman, F. L., Almeida, D., Altenschmidt, J., Altman, S., Anadkat, S., et al. Gpt-4 technical report. *arXiv preprint arXiv:2303.08774*, 2023.
- Berkaya, S. K., Uysal, A. K., Gunal, E. S., Ergin, S., Gunal, S., and Gulmezoglu, M. B. A survey on ecg analysis. *Biomedical Signal Processing and Control*, 43:216–235, 2018.
- Brown, T., Mann, B., Ryder, N., Subbiah, M., Kaplan, J. D., Dhariwal, P., Neelakantan, A., Shyam, P., Sastry, G., Askell, A., et al. Language models are few-shot learners. *Advances in neural information processing systems*, 33: 1877–1901, 2020.
- Chan, N., Parker, F., Bennett, W., Wu, T., Jia, M. Y., Fackler, J., and Ghobadi, K. Leveraging llms for multimodal medical time series analysis. In *Machine Learning for Healthcare Conference*. PMLR, 2024.
- Chen, K., Zhang, Z., Zeng, W., Zhang, R., Zhu, F., and Zhao, R. Shikra: Unleashing multimodal llm’s referential dialogue magic. *arXiv preprint arXiv:2306.15195*, 2023.
- Chu, Y., Xu, J., Zhou, X., Yang, Q., Zhang, S., Yan, Z., Zhou, C., and Zhou, J. Qwen-audio: Advancing universal audio understanding via unified large-scale audio-language models. *arXiv preprint arXiv:2311.07919*, 2023.
- Fu, G., Zheng, J., Abudayyeh, I., Ani, C., Rakovski, C., Ehwerhemuepha, L., Lu, H., Guo, Y., Liu, S., Chu, H., et al. Cardiogpt: An ecg interpretation generation model. *IEEE Access*, 2024.
- Gao, P., Han, J., Zhang, R., Lin, Z., Geng, S., Zhou, A., Zhang, W., Lu, P., He, C., Yue, X., et al. Llama-adapter v2: Parameter-efficient visual instruction model. *arXiv preprint arXiv:2304.15010*, 2023.

- Gong, T., Lyu, C., Zhang, S., Wang, Y., Zheng, M., Zhao, Q., Liu, K., Zhang, W., Luo, P., and Chen, K. Multimodal-gpt: A vision and language model for dialogue with humans. *arXiv preprint arXiv:2305.04790*, 2023.
- Gow, B., Pollard, T., Nathanson, L. A., Johnson, A., Moody, B., Fernandes, C., Greenbaum, N., Waks, J. W., Eslami, P., Carbonati, T., et al. Mimic-iv-ecg: Diagnostic electrocardiogram matched subset. *Type: dataset*, 6:13–14, 2023.
- Hannun, A. Y., Rajpurkar, P., Haghpanahi, M., Tison, G. H., Bourn, C., Turakhia, M. P., and Ng, A. Y. Cardiologist-level arrhythmia detection and classification in ambulatory electrocardiograms using a deep neural network. *Nature medicine*, 25(1):65–69, 2019.
- He, K., Mao, R., Lin, Q., Ruan, Y., Lan, X., Feng, M., and Cambria, E. A survey of large language models for healthcare: from data, technology, and applications to accountability and ethics. *Information Fusion*, 118:102963, 2025. ISSN 1566-2535. doi: <https://doi.org/10.1016/j.inffus.2025.102963>. URL <https://www.sciencedirect.com/science/article/pii/S1566253525000363>.
- Hong, S., Wu, M., Zhou, Y., Wang, Q., Shang, J., Li, H., and Xie, J. ENCASE: an ensemble classifier for ECG classification using expert features and deep neural networks. In *CinC*, 2017. doi: 10.22489/CinC.2017.178-245. URL <https://doi.org/10.22489/CinC.2017.178-245>.
- Hong, S., Zhou, Y., Wu, M., Shang, J., Wang, Q., Li, H., and Xie, J. Combining deep neural networks and engineered features for cardiac arrhythmia detection from ECG recordings. *Physiological Measurement*, 40(5):054009, Jun 2019. doi: 10.1088/1361-6579/ab15a2. URL <https://doi.org/10.1088%2F1361-6579%2Fab15a2>.
- Hong, S., Zhou, Y., Shang, J., Xiao, C., and Sun, J. Opportunities and challenges of deep learning methods for electrocardiogram data: A systematic review. *Computers in Biology and Medicine*, pp. 103801, 2020.
- Kim, Y., Xu, X., McDuff, D., Breazeal, C., and Park, H. W. Health-llm: Large language models for health prediction via wearable sensor data. *arXiv preprint arXiv:2401.06866*, 2024.
- Kiyasseh, D., Zhu, T., and Clifton, D. A. Clocs: Contrastive learning of cardiac signals across space, time, and patients. In Meila, M. and Zhang, T. (eds.), *Proceedings of the 38th International Conference on Machine Learning*, volume 139 of *Proceedings of Machine Learning Research*, pp. 5606–5615. PMLR, 18–24 Jul 2021. URL <https://proceedings.mlr.press/v139/kiyasseh21a.html>.
- Lan, X., Ng, D., Hong, S., and Feng, M. Intra-inter subject self-supervised learning for multivariate cardiac signals. In *Proceedings of the AAAI Conference on Artificial Intelligence*, volume 36, pp. 4532–4540, 2022.
- Lan, X., Yan, H., Hong, S., and Feng, M. Towards enhancing time series contrastive learning: A dynamic bad pair mining approach. In *The Twelfth International Conference on Learning Representations*, 2024. URL <https://openreview.net/forum?id=K2c04ulKXn>.
- Li, C., Wong, C., Zhang, S., Usuyama, N., Liu, H., Yang, J., Naumann, T., Poon, H., and Gao, J. Llava-med: Training a large language-and-vision assistant for biomedicine in one day. *arXiv preprint arXiv:2306.00890*, 2023a.
- Li, J., Liu, C., Cheng, S., Arcucci, R., and Hong, S. Frozen language model helps ecg zero-shot learning. In Oguz, I., Noble, J., Li, X., Styner, M., Baumgartner, C., Rusu, M., Heinmann, T., Kontos, D., Landman, B., and Dawant, B. (eds.), *Medical Imaging with Deep Learning*, volume 227 of *Proceedings of Machine Learning Research*, pp. 402–415. PMLR, 10–12 Jul 2024. URL <https://proceedings.mlr.press/v227/li24a.html>.
- Li, K., He, Y., Wang, Y., Li, Y., Wang, W., Luo, P., Wang, Y., Wang, L., and Qiao, Y. Videochat: Chat-centric video understanding. *arXiv preprint arXiv:2305.06355*, 2023b.
- Liu, F., Liu, C., Zhao, L., Zhang, X., Wu, X., Xu, X., Liu, Y., Ma, C., Wei, S., He, Z., et al. An open access database for evaluating the algorithms of electrocardiogram rhythm and morphology abnormality detection. *Journal of Medical Imaging and Health Informatics*, 8(7):1368–1373, 2018.
- Liu, F., Zhu, T., Wu, X., Yang, B., You, C., Wang, C., Lu, L., Liu, Z., Zheng, Y., Sun, X., et al. A medical multimodal large language model for future pandemics. *NPJ Digital Medicine*, 6(1):226, 2023.
- Liu, H., Li, C., Li, Y., and Lee, Y. J. Improved baselines with visual instruction tuning. In *Proceedings of the IEEE/CVF Conference on Computer Vision and Pattern Recognition*, pp. 26296–26306, 2024a.
- Liu, H., Li, C., Wu, Q., and Lee, Y. J. Visual instruction tuning. *Advances in neural information processing systems*, 36, 2024b.
- Liu, R., Bai, Y., Yue, X., and Zhang, P. Teach multimodal llms to comprehend electrocardiographic images. *arXiv preprint arXiv:2410.19008*, 2024c.

- Lyu, C., Wu, M., Wang, L., Huang, X., Liu, B., Du, Z., Shi, S., and Tu, Z. Macaw-llm: Multi-modal language modeling with image, audio, video, and text integration. *arXiv preprint arXiv:2306.09093*, 2023.
- McKeen, K., Oliva, L., Masood, S., Toma, A., Rubin, B., and Wang, B. Ecg-fm: An open electrocardiogram foundation model. *arXiv preprint arXiv:2408.05178*, 2024.
- Moor, M., Banerjee, O., Abad, Z. S. H., Krumholz, H. M., Leskovec, J., Topol, E. J., and Rajpurkar, P. Foundation models for generalist medical artificial intelligence. *Nature*, 616(7956):259–265, 2023.
- Pi, R., Gao, J., Diao, S., Pan, R., Dong, H., Zhang, J., Yao, L., Han, J., Xu, H., and Zhang, L. K. T. Detgpt: Detect what you need via reasoning. *arXiv preprint arXiv:2305.14167*, 2023.
- Ribeiro, A. H., Paixao, G., Lima, E. M., Ribeiro, M. H., Pinto Filho, M. M., Gomes, P. R., Oliveira, D. M., Meira Jr, W., Schon, T. B., and Ribeiro, A. L. P. Code-15%: A large scale annotated dataset of 12-lead ecgs. *Zenodo*, Jun, 9, 2021.
- Shivashankara, K. K., Shervedani, A. M., Clifford, G. D., Reyna, M. A., Sameni, R., et al. Ecg-image-kit: a synthetic image generation toolbox to facilitate deep learning-based electrocardiogram digitization. *Physiological Measurement*, 45(5):055019, 2024.
- Smulyan, H. The computerized ecg: Friend and foe. *The American Journal of Medicine*, 132(2):153–160, 2019. ISSN 0002-9343. doi: <https://doi.org/10.1016/j.amjmed.2018.08.025>. URL <https://www.sciencedirect.com/science/article/pii/S0002934318308532>.
- Thirunavukarasu, A. J., Ting, D. S. J., Elangovan, K., Gutierrez, L., Tan, T. F., and Ting, D. S. W. Large language models in medicine. *Nature medicine*, 29(8):1930–1940, 2023.
- Touvron, H., Lavril, T., Izacard, G., Martinet, X., Lachaux, M.-A., Lacroix, T., Rozière, B., Goyal, N., Hambro, E., Azhar, F., et al. Llama: Open and efficient foundation language models. *arXiv preprint arXiv:2302.13971*, 2023.
- Wagner, P., Strodthoff, N., Bousseljot, R.-D., Kreiseler, D., Lunze, F. I., Samek, W., and Schaeffter, T. Ptb-xl, a large publicly available electrocardiography dataset. *Scientific data*, 7(1):1–15, 2020.
- Wan, Z., Liu, C., Wang, X., Tao, C., Shen, H., Peng, Z., Fu, J., Arcucci, R., Yao, H., and Zhang, M. Electrocardiogram instruction tuning for report generation. *arXiv preprint arXiv:2403.04945*, 2024.
- Yu, H., Guo, P., and Sano, A. Zero-shot ecg diagnosis with large language models and retrieval-augmented generation. In *Machine Learning for Health (ML4H)*, pp. 650–663. PMLR, 2023.
- Zhang, H., Li, X., and Bing, L. Video-llama: An instruction-tuned audio-visual language model for video understanding. *arXiv preprint arXiv:2306.02858*, 2023.
- Zhao, Y., Zhang, T., Wang, X., Han, P., Chen, T., Huang, L., Jin, Y., and Kang, J. Ecg-chat: A large ecg-language model for cardiac disease diagnosis. *arXiv preprint arXiv:2408.08849*, 2024.
- Zheng, J., Zhang, J., Danioko, S., Yao, H., Guo, H., and Rakovski, C. A 12-lead electrocardiogram database for arrhythmia research covering more than 10,000 patients. *Scientific data*, 7(1):48, 2020.
- Zhu, D., Chen, J., Shen, X., Li, X., and Elhoseiny, M. Minigt-4: Enhancing vision-language understanding with advanced large language models. *arXiv preprint arXiv:2304.10592*, 2023.
- Zhu, L., Mou, W., Wu, K., Lai, Y., Lin, A., Yang, T., Zhang, J., and Luo, P. Multimodal chatgpt-4v for electrocardiogram interpretation: Promise and limitations. *Journal of Medical Internet Research*, 26:e54607, 2024.
- Zhu, Z., Lan, X., Zhao, T., Guo, Y., Kojodjojo, P., Xu, Z., Liu, Z., Liu, S., Wang, H., Sun, X., et al. Identification of 27 abnormalities from multi-lead ecg signals: An ensemble se_resnet framework with sign loss function. *Physiological Measurement*, 42(6):065008, 2021.

A. Diagnosis Guider Prompts

```
# Your task: Interpret the provided ECG image, identify key features and abnormalities in
↳ each lead, and generate a clinical diagnosis that is supported by the observed
↳ evidence.

## Key objectives:
1. Simulate a Realistic Diagnostic Process: The interpretation should reflect how a
↳ doctor would analyze an ECG, ask clarifying questions, and arrive at a diagnosis.
2. Grounded ECG Understanding: The analysis should be based on specific ECG features
↳ and explicitly reference these features as evidence.
3. Evidence-Based Reasoning: The diagnosis should be supported by clear, logical
↳ reasoning tied to the ECG findings.

## Guidelines for the ECG analysis:
1. Data:
  ECG image: an image that display the 12-lead ECG tracings. Make the task centered
  ↳ on the ECG image, assuming direct ECG image analysis.
  Machine measurements: A time-ordered list of ECG features computed for each
  ↳ heartbeat in every lead. Each entry in the list corresponds to the features
  ↳ calculated for a single heartbeat.

2. Act as a cardiologist and use medical knowledge to analyze the provided ECG image
↳ step-by-step:
  Initial Analysis: Analyze the provided ECG image to identify key features such as
  ↳ rhythm, intervals, and any apparent abnormalities.
  Detailed Reasoning: Explain your thought process step-by-step, referencing
  ↳ specific ECG features (e.g., "The ST segment is elevated in leads V1-V4, which
  ↳ suggests anterior myocardial infarction").
  Evidence-Based Diagnosis: Propose a diagnosis or differential diagnoses,
  ↳ justifying your conclusions with explicit ECG data.

3. When analyzing the ECG image, carefully analyze each lead:
  Lead I: Examine the QRS amplitude and duration, along with ST segment and T wave
  ↳ morphology. Abnormalities may indicate lateral wall issues such as left
  ↳ ventricular hypertrophy, bundle branch block, or lateral ischemia/infarction.
  Lead II: Look at the P wave amplitude and duration to assess right or left atrial
  ↳ enlargement; the PR interval can reveal conduction delays. ST and T wave
  ↳ changes here suggest inferior wall ischemia or infarction.
  Leads III and aVF: Primarily reflect inferior wall status. Abnormal Q waves, along
  ↳ with ST segment and T wave changes, point toward inferior infarction or
  ↳ ischemia.
  Lead aVL: Focuses on the high lateral region; QRS, ST, and T wave abnormalities
  ↳ here suggest high lateral ischemia or infarction.
  Lead aVR: ST elevation may indicate left main or multivessel disease, and T wave
  ↳ inversion can be associated with ventricular arrhythmia.
  Lead V1: An increased R wave, a characteristic rsR pattern, and ST-T changes help
  ↳ identify right ventricular hypertrophy, right bundle branch block, or
  ↳ ischemia.
  Leads V2-V4: Assessing the anterior or anteroseptal regions. The presence of Q
  ↳ waves, along with ST segment and T wave deviations, suggests anterior wall
  ↳ infarction or ischemia.
  Leads V5-V6: Focus on the lateral wall, where similar QRS, ST, and T wave changes
  ↳ can indicate lateral ischemia or infarction.

4. When analyzing the machine measurements, you should aware that:
  a. If any abnormalities appear in the computed measurements that are not mentioned
  ↳ in the report, you must strictly follow and trust the report.
  b. Evaluate and interpret the machine measurements as if you had computed them
  ↳ yourself. In your analysis, refer to these values as your own computed
  ↳ measurements rather than using phrases like "machine measurements provided".
```

Guidelines for the response generation:

1. Synthesize your findings to deduce a likely diagnosis or set of diagnoses. Clearly
→ explain how the evidence supports your conclusion.
2. Ensure your diagnosis is comprehensive and strictly based on the report. Do not
→ include diagnosis that not mentioned in the report.
3. Make sure your diagnosis are grounded in the given ECG image and machine
→ measurements, and you should explicitly reference (e.g., specify lead and the
→ position of the abnormal heartbeat).
4. Strictly follow the output format and requirements specified in your task
→ instructions.
5. The given report only served as the ground truth for you to analyze the ECG image.
→ The generated text must not show that you are aware of the existence of the
→ report.
6. Never make up explanations.

ECG Report:

{{report}}

ECG Machine Measurements:

{{machine_measurements}}

Generataion rule

The generated text must not show that you are aware of the existence of the report.

Do not include phrases like "Based the report", or "Given the ECG report".

The primary objective is to analyze the ECG and identify evidence that supports the

→ results. The analysis should focus solely on the ECG itself, never analyze the report.

Present your work in this format:

****Response:**** [Comprehensive response following the task's guidelines, strictly based on

→ the report. Using a complete paragraph with more natural expression. Do not use a list

→ format. Limit your responses within 300 words.]

B. Evaluation Prompts

Your task: Evaluate the alignment and quality of a generated ECG interpretation by
↪ comparing it to a ground truth clinician's interpretation.

Evaluation Criteria:

1. **DiagnosisAccuracy**: Evaluates whether the generated diagnosis is correct, specific, and supported by ECG findings.
↪
 - Scoring
 - +2 per diagnosis: Each correctly identified key diagnosis with supporting ECG features.
 - +1 per diagnosis: Each mostly correct diagnosis but lacking key supporting details.
 - +0 per diagnosis: Each incorrect or vague diagnosis not supported by ECG features.
2. **AnalysisCompleteness**: Checks if all key ECG components (rhythm, intervals, waveforms, and lead-specific findings) are discussed.
↪
 - Scoring
 - +1 per feature: For each correctly addressed key ECG feature (e.g., rhythm, PR interval, QRS duration, ST segment, T wave morphology).
 - +0 per missing feature: For each key feature omitted or inaccurately described.
3. **AnalysisRelevance**: Assesses whether each provided explanation directly supports the diagnosis.
↪
 - Scoring
 - +2 per feature or per lead: Each point that strongly supports the diagnosis with clear ECG evidence.
 - +1 per feature or per lead: Some points are relevant but not fully justified.
 - +0: Includes unrelated or misleading explanations.
4. **LeadAssessmentCoverage**: Evaluates how many of the 12 ECG leads are analyzed.
↪
 - Scoring
 - +1 per lead: For each lead adequately assessed.
 - +0 per missing lead: For each lead omitted or inaccurately described.
5. **LeadAssessmentAccuracy**: Checks if the described lead findings (e.g., QRS, ST, T waves, amplitude, intervals, ST segments) match standard ECG interpretation.
↪
 - Scoring
 - +2 per lead: Findings closely match expected values.
 - +1 per lead: Findings are somewhat accurate but have minor inconsistencies.
 - +0 per lead: Findings contradict ECG norms.
6. **ECGFeatureGrounding**: Determines if the interpretation references actual ECG features (e.g., QRS amplitude, PR interval) instead of generic terms.
↪
 - Scoring (0-100)
 - 100: ECG findings are comprehensively cited, linked to diagnoses, and cover all relevant ECG features.
 - 80: ECG findings are explicitly cited and linked to diagnoses.
 - 50: Some ECG references exist but are incomplete.
 - 0: Lacks specific waveform references.
7. **EvidenceBasedReasoning**: Evaluates whether the diagnosis follows logical, evidence-supported steps.
↪
 - Scoring (0-100)
 - 100: Findings logically progress to diagnosis with thorough and clear justifications covering all necessary steps.
 - 80: Findings logically progress to diagnosis with clear justifications.
 - 50: Some reasoning exists but lacks complete step-by-step analysis.
 - 0: Reasoning is unclear or not derived from ECG findings.

- 8. ClinicalDiagnosticFidelity: Assesses if the model mimics how a clinician interprets
↳ an ECG, considering all relevant factors.
 - Scoring (0-100)
 - 100: The analysis follows a structured clinical approach and considers all
↳ relevant clinical factors.
 - 80: The analysis follows a structured clinical approach.
 - 50: Some clinical reasoning is present but incomplete.
 - 0: The approach lacks structured clinical reasoning.


NOTE: Each score must be calculated based on strict criteria to ensure objective
↳ evaluation.

Generated ECG Interpretation:
{{generated}}

Ground Truth Clinician's Interpretation:
{{groundtruth}}

****Response:**** [Please organize your output in a JSON format for each criterion, including
↳ a brief explanation for each aspect. Strictly follow this JSON format that records
↳ every scoring for each criterion: {'DiagnosisAccuracy': [{'Score': , 'Explanation':},
↳ ...], 'AnalysisCompleteness': [{'Score': , 'Explanation':}, ...],
↳ 'AnalysisRelevance': [{'Score': , 'Explanation':}, ...],
↳ 'LeadAssessmentCoverage': [{'Score': , 'Explanation':}, ...],
↳ 'LeadAssessmentAccuracy': [{'Score': , 'Explanation':}, ...],
↳ 'ECGFeatureGrounding': [{'Score': , 'Explanation':}, ...],
↳ 'EvidenceBasedReasoning': [{'Score': , 'Explanation':}, ...],
↳ 'ClinicalDiagnosticFidelity': [{'Score': , 'Explanation':}, ...]]

C. Case Studies by Cardiologist

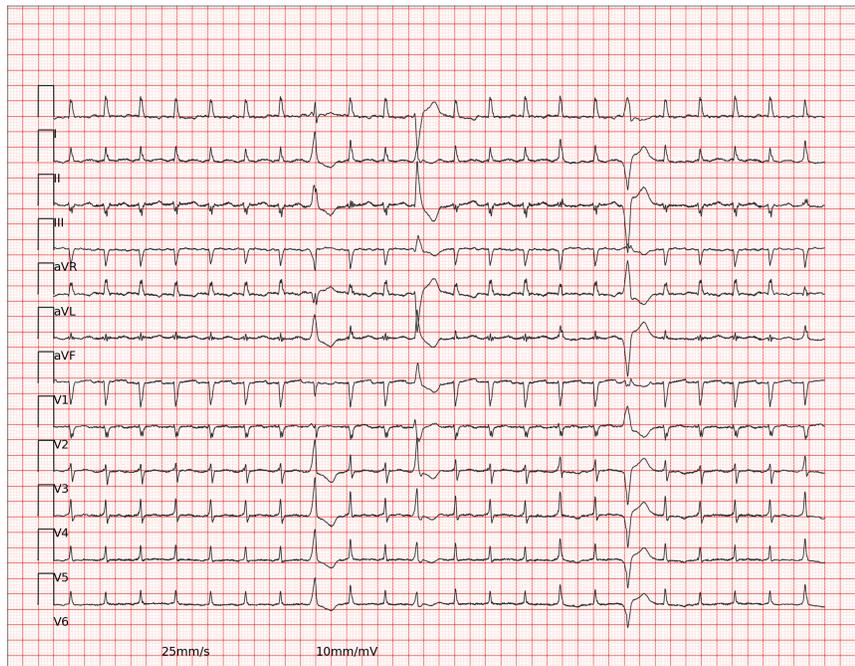


GEM Interpretation: Upon analyzing the provided ECG image and computed measurements, several key features and abnormalities are evident. The heart rate is elevated at 121 bpm, indicating a rapid ventricular response. The rhythm appears irregular, consistent with atrial fibrillation, as seen by the absence of distinct P waves and the presence of irregular RR intervals. In the precordial leads, particularly V1 through V4, there are low QRS voltages, which could suggest conditions such as pericardial effusion or obesity. The QRS complexes in these leads are notably diminished, which aligns with the low voltage finding. Furthermore, the presence of premature ventricular contractions (PVCs) is suggested by the irregular RR intervals and the occasional absence of preceding P waves. The ST-T changes observed in the anterolateral leads, particularly V2 to V6, are nonspecific but could indicate ischemic changes or other cardiac conditions. The QTc intervals are prolonged in several leads, which may suggest a risk for arrhythmias. The combination of these findings, including the rapid ventricular response, irregular rhythm, low QRS voltages, and nonspecific ST-T changes, supports a diagnosis of atrial fibrillation with rapid ventricular response and PVCs. The low QRS voltages in the precordial leads add another layer of complexity, potentially indicating underlying cardiac conditions. The ECG is classified as abnormal due to these findings, which warrant further clinical correlation and possibly additional diagnostic testing to determine the underlying cause of these abnormalities.

Target Answer (GPT-4o): Upon examining the provided ECG, we identify several critical features and abnormalities. The first notable observation is the presence of atrial fibrillation, characterized by an irregular rhythm and lack of distinct P waves in the ECG tracings. This is supported by the irregular RR intervals, which vary significantly between heartbeats. The heart rate is elevated at 132 bpm, indicative of a rapid ventricular response. Additionally, the ECG reveals frequent premature ventricular contractions, or multifocal PVCs, with varying QRS complex morphologies, suggesting aberrant ventricular conduction or possible ventricular ectopy. Lead I and precordial leads demonstrate low QRS voltages, which can correlate with various conditions, including myocardial infarction or myopathies. In terms of ST-T changes, the anterior leads (V1-V4) may suggest possible anterior myocardial infarction, as evidenced by QRS amplitude and morphology abnormalities. Moreover, inferior and lateral leads, particularly lead aVL, display nonspecific ST-T changes, which could indicate ischemia or other pathological processes. A comprehensive review of the QT interval reveals prolonged QTc intervals particularly in the precordial leads, reflecting altered repolarization processes. Furthermore, no distinct high, lateral, or inferior lead pattern suggests definitive infarction, but the changes are significant enough to warrant further clinical investigation and correlation with clinical symptoms and history. Conclusively, the ECG findings strongly support a diagnosis of atrial fibrillation with rapid ventricular response, accompanied by frequent multifocal PVCs or aberrant ventricular conduction. The anterior region's ST-T changes and precordial low QRS voltages emphasize the need for further assessment, possibly indicating a previous or ongoing ischemic event. However, additional clinical context and investigations are essential to corroborate these findings and formulate a complete treatment strategy.

Ground-truth Report: Atrial fibrillation with rapid ventricular response with frequent multifocal PVCs or aberrant ventricular conduction; Possible anterior infarct; Inferior/lateral ST-T changes are nonspecific; Low QRS voltages in precordial leads; Abnormal ECG.

PULSE Interpretation: The electrocardiogram (ECG) reveals atrial fibrillation with a rapid ventricular response, accompanied by premature ventricular contractions (PVCs) or aberrant ventricular conduction, indicating irregular heart rhythm. Additionally, the ECG suggests possible anterior myocardial infarction, although the age of the infarct cannot be determined. The inferior and lateral ST-T changes are nonspecific, and the QRS voltages are low in the precordial leads, contributing to an overall abnormal ECG pattern.

Figure 4. Cardiologist Evaluations. Blue: Findings exceeding expert expectations. Yellow: Findings with differing expert opinions.

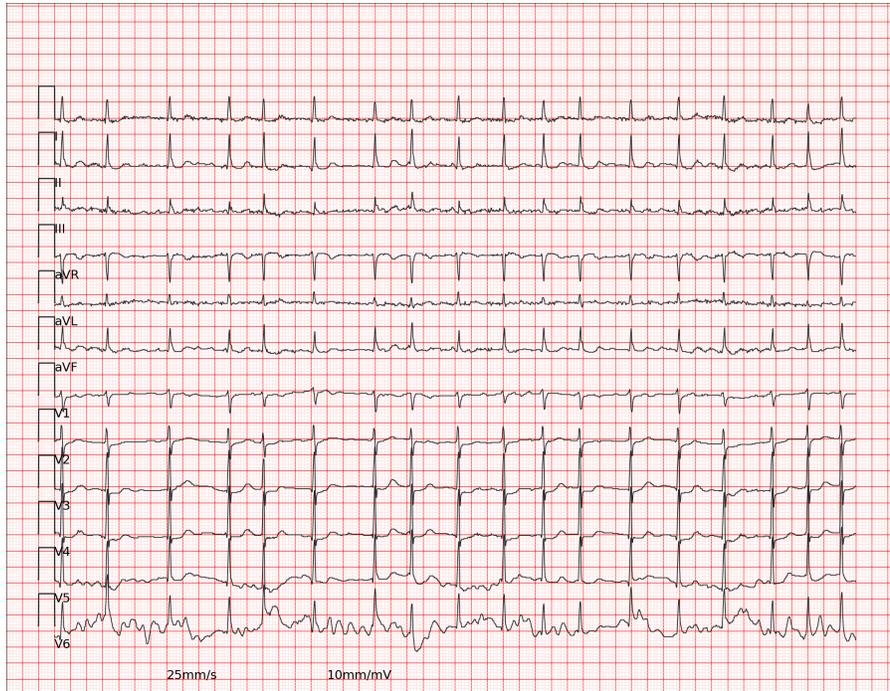


GEM Interpretation: Upon analyzing the provided ECG image and computed measurements, several key features and abnormalities are evident. The rhythm is irregular, with a heart rate of 102 bpm, which is consistent with atrial fibrillation. This is supported by the absence of distinct P waves across multiple leads, such as Lead I and Lead II, and the presence of irregular RR intervals, as seen in the measurements. The QRS complexes are generally narrow, indicating a supraventricular origin, which aligns with atrial fibrillation. In Lead I, the QRS amplitude and duration are within normal limits, but the irregular rhythm is evident. Lead II shows a similar pattern with absent P waves and irregular intervals, further supporting atrial fibrillation. Leads III and aVF, which reflect the inferior wall, do not show significant Q waves or ST segment changes that would suggest acute ischemia or infarction. Lead aVL, focusing on the high lateral region, also does not exhibit significant abnormalities in the QRS, ST, or T wave morphology. Lead aVR, often used to assess global ischemic changes, does not show significant ST elevation or T wave inversion, which might have indicated more severe conditions like left main coronary artery disease. In the precordial leads V1 through V6, the QRS complexes are narrow, and there are no significant ST segment elevations or depressions that would suggest acute myocardial infarction. Overall, the ECG findings are consistent with atrial fibrillation, characterized by the absence of distinct P waves and irregular RR intervals. The lack of significant ST segment changes or pathological Q waves across the leads suggests that there is no acute ischemic event present. This analysis supports the diagnosis of atrial fibrillation as the primary abnormality, with no evidence of acute myocardial infarction or other significant ischemic changes.

Target Answer (GPT-4o): Upon examining the provided ECG image, several key features indicate the presence of atrial fibrillation (AF). The irregularly irregular rhythm, as evidenced by varying RR intervals, supports this observation. Notably, the RR intervals range widely, which is a hallmark of AF due to inconsistent atrial impulses reaching the ventricles. Additionally, the absence of distinct P waves on leads with traditionally observable P waves is another pivotal indicator; this suggests a loss of organized atrial electrical activity, consistent with AF. The heart rate approximates 107 bpm, indicative of a moderate ventricular rate, aligning with a presentation of AF with a controlled rate. No consistent pattern in the PR interval is observed, further corroborating the diagnosis of AF, as atrial activity is chaotic and not leading to a predictable conduction pattern through the AV node. Analyzing leads such as lead II, where P wave characteristics are more observable under normal conditions, we see minimal or absent P wave formations. The variability and inconsistency of PR intervals signify typical conduction during AF episodes, contributing to its erratic rhythm. In terms of additional abnormalities, minor QRS and T wave deviations across various leads were noted, but they did not form a pathological pattern significant enough to suggest additional diagnoses such as ischemia or infarction. The ST segment and QTc intervals are within normal limits given the context of AF. This ECG analysis supports a diagnosis of atrial fibrillation, underpinning the clinical picture with evidenced irregular RR intervals and absent P waves. My interpretation remains focused on the features aligning with AF without assuming further conditions not supported by the ECG findings.

Ground-truth Report: Atrial fibrillation.

PULSE Interpretation: The electrocardiogram (ECG) reveals atrial fibrillation, characterized by an irregular heart rhythm with no discernible P waves, indicating a lack of coordinated atrial contraction. Additionally, the ECG shows a repolarization abnormality, suggesting a potential issue with the heart's recovery phase after contraction. The probable rate-related abnormality implies that the repolarization abnormality may be related to the atrial fibrillation, which is causing an irregular heart rate.

Figure 5. Cardiologist Evaluations. Blue: Findings exceeding expert expectations. Yellow: Findings with differing expert opinions.

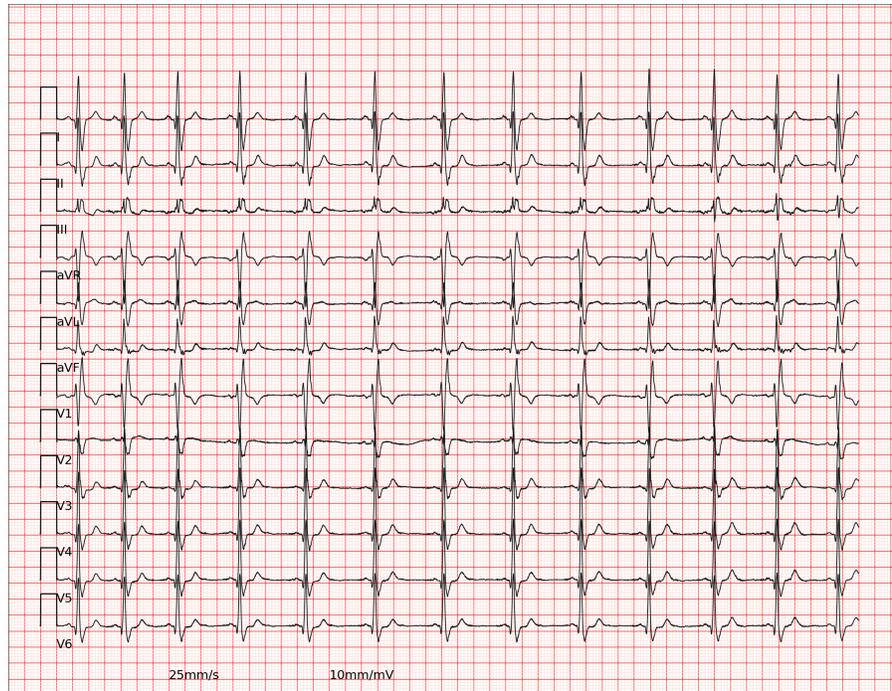


GEM Interpretation: Upon analyzing the provided ECG image and computed measurements, several key features and abnormalities are evident. The rhythm is sinus, as indicated by the consistent P wave morphology across leads and a **regular heart rate of 75 bpm**. In Lead I, the QRS complex shows a prolonged duration, which is consistent with a right bundle branch block (RBBB). **This is further supported by the characteristic rsR' pattern in Lead V1, a hallmark of RBBB.** The QRS duration in Lead I is notably prolonged, ranging from 134 to 150 ms, which is indicative of this conduction delay. In Lead II, the PR interval is within normal limits, but the QRS duration is slightly prolonged, further supporting the presence of RBBB. The ST segment and T wave morphology in Lead II do not show significant deviations, suggesting no acute ischemic changes. Leads III and aVF, which reflect the inferior wall, do not exhibit pathological Q waves or significant ST segment changes, ruling out inferior wall ischemia or infarction. Lead aVL shows no significant abnormalities in the QRS complex or ST segment, which suggests no high lateral ischemia. Lead aVR does not show ST elevation, which would have indicated left main coronary artery disease. In the precordial leads V2-V4, the QRS complexes are widened, consistent with RBBB, and there are no significant ST segment elevations or depressions to suggest anterior wall ischemia. Overall, the ECG findings are consistent with a right bundle branch block, as evidenced by the prolonged QRS duration and characteristic patterns in the precordial leads. The sinus rhythm and absence of acute ischemic changes suggest that the conduction abnormality is the primary finding. This analysis aligns with the diagnosis of an abnormal ECG due to the presence of RBBB.

Target Answer (GPT-4o): The 12-lead ECG image reveals key characteristics supporting the diagnosis of a right bundle branch block (RBBB) in the context of sinus rhythm. In Lead I, the QRS duration extends to approximately 134 milliseconds, **coupled with a notable upslope in the ST segment across multiple heartbeats.** The amplitude of the QRS complex peaks at approximately 2.463 mV, which might suggest normal conduction patterns as the lead does not display criteria typical of a RBBB. Lead V1, however, showcases the classic features of RBBB with an RsR' pattern evidenced by a marginally increased QRS duration and the presence of a characteristic "rabbit ears" morphology, as indicated by significant amplitudes and uniform morphology, compatible with RBBB diagnostics. Furthermore, **different elevations and variations, including notable upslope ST segments observed in leads V1 to V3, can further affirm the presence of bundle branch blocks without overt ischemic patterns.** Additional leads such as II, III, aVR, and aVF do not exhibit significant anomalies, aside from mild ST variations, which are consistent across these leads without inferring significant ischemia or infarction. The computed average heart rate is normal, around 73 bpm, **with consistent RR intervals suggesting regular sinus rhythm.** Although the overall ECG demonstrates sinus rhythm, the broad QRS complexes and rightward axis deviation are hallmarks of a right bundle branch block. These observations, together with consistent measurements, fully support the diagnosis of RBBB. No acute ischemic changes are discerned based on ST segment and T wave morphology across most of the ECG leads, correlating with the absence of acute myocardial infarction signs. The ECG also appears abnormal due to the prolonged QRS complex duration, supporting the finding of a RBBB when considered in conjunction with the specific electrocardiographic patterns observed.

Ground-truth Report: Sinus rhythm; Right bundle branch block; Abnormal ECG.

PULSE Interpretation: The electrocardiogram (ECG) reveals a sinus tachycardia, characterized by a rapid heart rate with a normal P-wave axis, indicating that the sinus node is functioning normally but at an accelerated pace. Additionally, the ECG shows a right bundle branch block (RBBB), which is a conduction abnormality affecting the right ventricle, resulting in a delayed activation of the right ventricle. Overall, the ECG suggests a diagnosis of sinus tachycardia with right bundle branch block.

Figure 6. Cardiologist Evaluations. Blue: Findings exceeding expert expectations. Yellow: Findings with differing expert opinions.

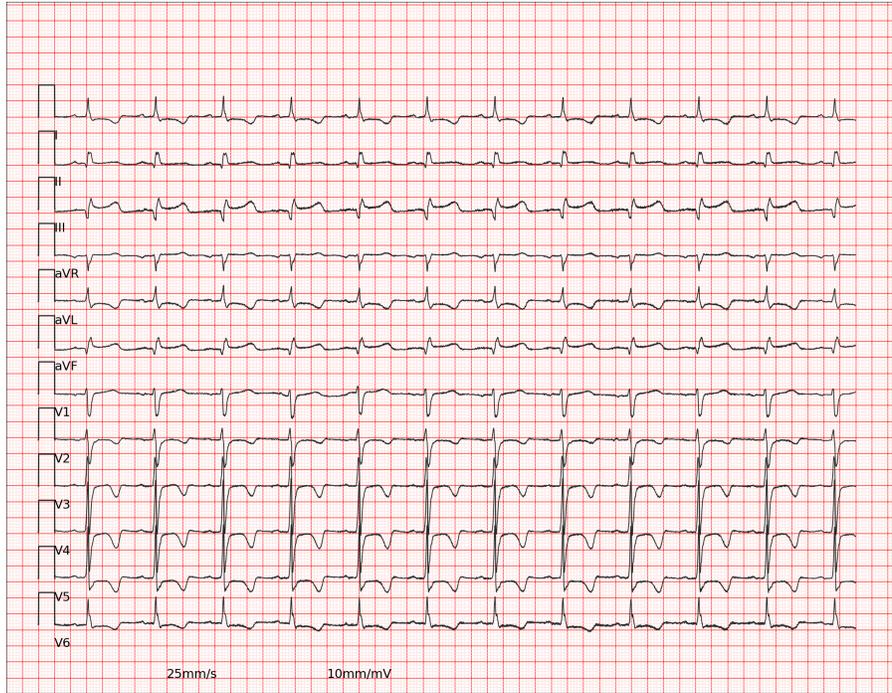


GEM Interpretation: Upon analyzing the provided ECG image and computed measurements, several key features and abnormalities are evident. The rhythm is sinus, as indicated by the regular P waves preceding each QRS complex. However, the QT interval is prolonged across multiple leads, with QTc values exceeding 500 ms in several instances, suggesting a risk for arrhythmias. In leads II, III, and aVF, there are notable changes consistent with an inferior infarct, including abnormal Q waves and ST segment deviations. Additionally, leads V2 to V4 show ST segment elevation and T wave inversion, which are indicative of anterior wall ischemia or infarction. The lateral leads, particularly I, aVL, V5, and V6, exhibit ST-T changes, such as ST segment depression and T wave inversion, which may suggest lateral myocardial ischemia. The presence of these changes in multiple leads supports the diagnosis of myocardial ischemia. The combination of these findings, including the prolonged QT interval and the ST-T changes, suggests a complex cardiac condition involving both ischemic and infarct patterns. The ECG is abnormal, with evidence pointing towards an inferior infarct, anterior and lateral ischemia, and a prolonged QT interval, all of which require further clinical correlation and management.

Target Answer (GPT-4o): Analyzing the provided 12-lead ECG shows evidence of an inferior myocardial infarction, as evidenced by significant deviations in the leads primarily reflective of the inferior wall, such as leads II, III, and aVF, which exhibit abnormalities in ST segments and T wave morphology. The QRS complexes and ST segment changes suggest an old or acute inferior infarct. Additionally, ST-T changes are notable in the anterior leads (V2-V4), indicating anterior or anteroseptal myocardial ischemia. The presence of significant ST elevation or depression, coupled with T wave inversion in these leads, corroborates the potential for ischemia in the anterior regions of the heart. In the lateral leads (I, aVL, V5, and V6), the presence of ST segment and T wave changes could indicate lateral wall ischemia or infarction. Specifically, notable changes in ST form and T wave morphology support this suspicion, which aligns with the observed abnormalities. The heart rate, computed at an average of 70 beats per minute, appears normal as it shows regular RR intervals, supporting a sinus rhythm without apparent rhythm disturbances. The computed PR intervals and QRS durations in different leads suggest no significant atrioventricular or bundle branch blocks, except for possible localized conduction delays that contribute to the differential ECG patterns observed. Collectively, these ECG findings demonstrate a combination of inferior infarct with the potential for ischemia notable in both anterior and lateral regions. This pattern of changes, particularly the ST-T variations across multiple leads, supports the likelihood of broadened myocardial involvement. The diagnosis aligns with the classification of an abnormal ECG, confirming myocardial ischemia and infarction.

Ground-truth Report: Sinus rhythm; Inferior infarct; Ant/septal and lateral ST-T changes may be due to myocardial ischemia; Abnormal ECG.

PULSE Interpretation: The electrocardiogram (ECG) reveals a sinus rhythm, indicating a normal heart rate and rhythm. However, it also shows an intraventricular conduction defect, which suggests a delay in the electrical conduction system of the heart. Furthermore, the ECG indicates evidence of a previous inferior myocardial infarction, although the age of the infarct cannot be determined. Additionally, there are ST-T changes observed in the anterior, septal, and lateral leads, which may be indicative of myocardial ischemia. Overall, the ECG is abnormal, suggesting cardiac disease. Clinical Diagnosis: Abnormal ECG with evidence of previous inferior myocardial infarction and possible myocardial ischemia.

Figure 7. Cardiologist Evaluations. Blue: Findings exceeding expert expectations. Yellow: Findings with differing expert opinions.

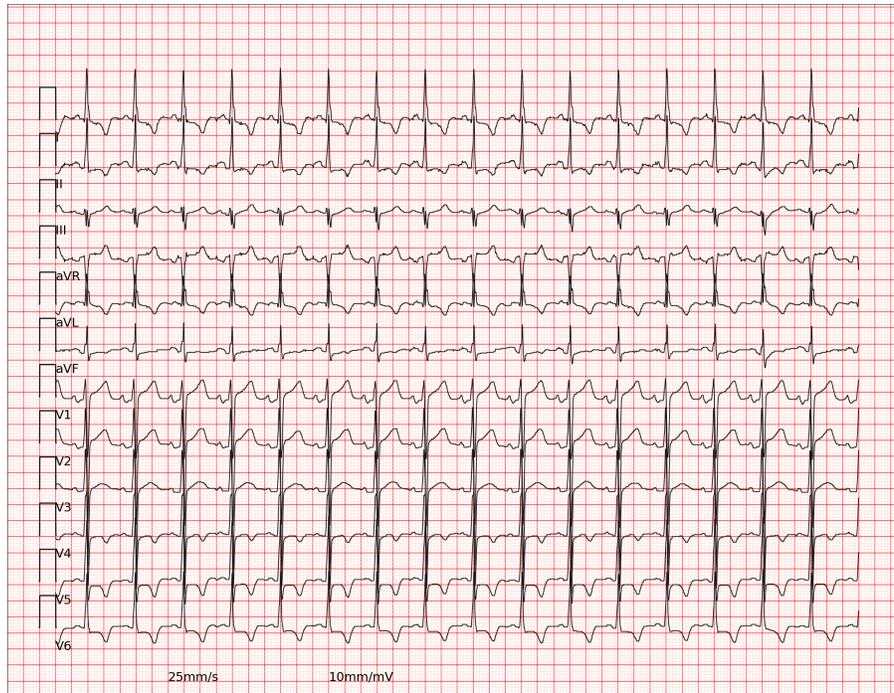


GEM Interpretation: Upon analyzing the provided ECG image and computed measurements, several key features and abnormalities are evident. The heart rate is elevated at 101 bpm, indicating sinus tachycardia. In Lead I, the QRS amplitude is notably high, suggesting left ventricular hypertrophy (LVH), which is further supported by the increased QRS amplitude in Leads V5 and V6. The presence of secondary repolarization abnormalities is suggested by the ST segment changes and T wave morphology in these leads. In Lead II, the P wave amplitude and duration are within normal limits, but the PR interval is slightly prolonged, which could indicate a first-degree AV block. The QRS complexes in Leads III and aVF are diminished, which might suggest an inferior wall issue, but the lack of significant Q waves or ST segment changes makes this less likely. Lead aVL shows a high QRS amplitude, consistent with LVH, and the T wave morphology is abnormal, which could indicate lateral ischemia. Lead aVR does not show significant ST elevation, reducing the likelihood of left main coronary artery disease. In Lead V1, the absence of an rsR' pattern and the lack of significant ST-T changes suggest that right ventricular hypertrophy or right bundle branch block is unlikely. Leads V2 to V4 show no significant Q waves or ST segment elevation, reducing the likelihood of an anterior myocardial infarction. Overall, the ECG findings are consistent with sinus tachycardia, probable left atrial enlargement, and left ventricular hypertrophy with secondary repolarization abnormalities. The absence of significant ST segment elevation or depression in the anterior leads suggests that acute ischemia or infarction is unlikely. The findings are consistent with chronic changes rather than acute pathology.

Target Answer (GPT-4o): In analyzing the ECG image, several key features and abnormalities are identified that support the clinical findings. The rhythm across leads is sinus, as observed from the consistent P wave morphology, notably prominent in Lead II, indicating a regular atrial rhythm. In Lead II, the P wave amplitude varies, but the durations range around 92 to 104 ms, suggesting left atrial enlargement. This enlargement is further evidenced by the increased P wave amplitude in several leads, including V1, where the amplitude is significantly high. The QRS complexes show increased amplitudes across multiple leads, particularly in V3 and V4, with corresponding secondary ST-T changes. This pattern of high voltage QRS complexes and secondary repolarization abnormalities is indicative of left ventricular hypertrophy (LVH). The computer measurements confirm increased QRS amplitude in these leads, supporting the presence of LVH. Additionally, there is anterior ST segment elevation seen in leads V2-V4, which can be attributed to the strain pattern associated with LVH rather than acute myocardial infarction, as it correlates with the presence of LVH. There is no significant Q wave development or evidence of acute ischemic injury that suggests an infarction. Overall, the ECG findings collectively indicate sinus rhythm with left atrial enlargement and left ventricular hypertrophy with secondary repolarization abnormalities. The anterior ST elevation observed is likely a consequence of LVH, consistent with a LVH strain pattern. These interpretations align with the clinical manifestations of a patient with left heart strain and structural changes, corroborated by the data synthesis from both visual ECG interpretation and detailed interval measurements.

Ground-truth Report: Sinus rhythm; Left atrial enlargement; LVH with secondary repolarization abnormality; Anterior ST elevation, probably due to LVH.

PULSE Interpretation: The electrocardiogram (ECG) reveals a sinus tachycardia, characterized by a rapid heart rate with a normal P-wave axis, indicating that the sinus node is functioning normally but at an accelerated pace. Additionally, the ECG shows evidence of biatrial enlargement, suggesting that both the left and right atria are enlarged, which can be a sign of cardiac chamber dilation or hypertrophy. Furthermore, the ECG demonstrates left ventricular hypertrophy (LVH) with secondary repolarization abnormality, indicating that the left ventricle is thickened, leading to abnormal electrical conduction and repolarization patterns.

Figure 8. Cardiologist Evaluations. Blue: Findings exceeding expert expectations. Yellow: Findings with differing expert opinions.

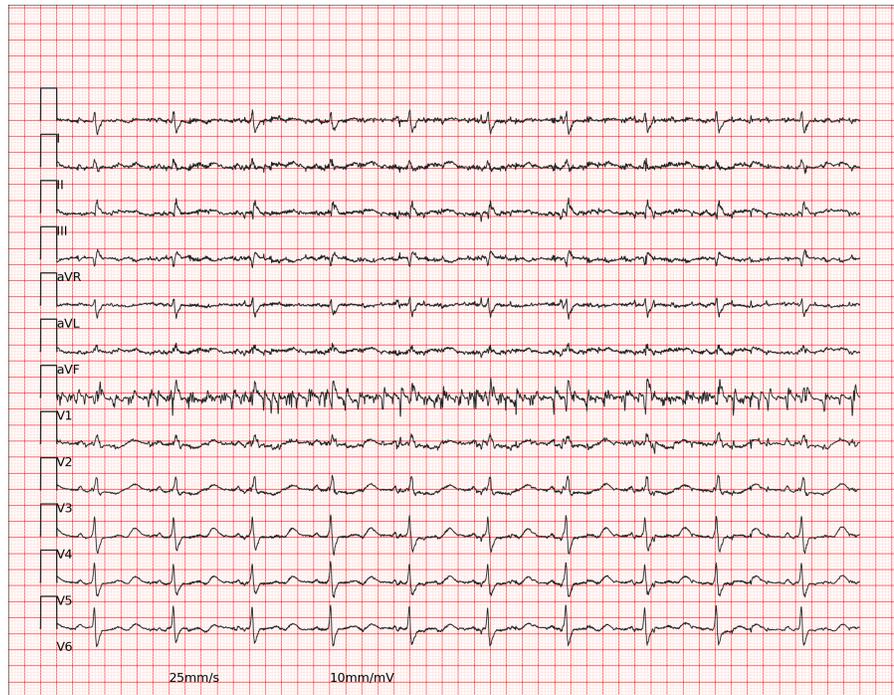

GEM Interpretation: Upon analyzing the provided ECG image and computed measurements, several key features and abnormalities are evident. The rhythm is sinus, as indicated by the consistent P wave morphology across leads, but there is a borderline first-degree A-V block, suggested by the prolonged PR intervals, particularly noticeable in leads I and II. The QRS complexes exhibit a right bundle branch block pattern, characterized by a widened QRS duration, particularly in lead V1, where an rsR' pattern is observed. The axis is rightward, which is consistent with the right bundle branch block. In the inferior leads (II, III, and aVF), there are nonspecific T wave changes, which could indicate ischemia or other non-specific changes. The QRS amplitude in lead III is notably low, which might suggest an inferior wall issue, although the changes are not definitive for infarction. **The lateral leads (I, aVL, V5, and V6) show some T wave changes, but these are nonspecific and do not clearly indicate lateral ischemia or infarction.** The QTc intervals are prolonged in several leads, which could be a concern for potential arrhythmogenic risk. The overall ECG is abnormal, with the combination of these findings suggesting a complex cardiac conduction issue. The rightward axis and right bundle branch block are significant, and the nonspecific T wave changes in the inferior and lateral leads warrant further clinical correlation and possibly additional diagnostic testing to rule out underlying ischemic heart disease or other cardiac conditions.

Target Answer (GPT-4o): Upon reviewing the 12-lead ECG image, the diagnosis reveals a sinus rhythm with a rightward axis and the presence of a right bundle branch block. The rhythm is consistent with a heart rate around 70 bpm, indicated by relatively regular RR intervals. The rightward axis is evident as lead I shows a smaller QRS amplitude compared to lead aVF, which reflects more positive deflections. A right bundle branch block is supported by the widened QRS complexes, particularly evident in leads V1 and V2, where an rsR' pattern could be observed, consistent with this conduction abnormality. Examination of the lateral leads (I, aVL, V5, V6) reveals mild T wave abnormalities, which are noted as nonspecific changes. In these leads, the morphology of the T waves shows slight variations that do not definitively point to ischemia but suggest some lateral wall changes. Despite these T wave changes being nonspecific, they may warrant further clinical correlation given the right bundle branch block already present. In summary, the ECG findings primarily indicate a right bundle branch block with a rightward axis. The T wave changes noted in the lateral leads are nonspecific and should be clinically correlated to assess their significance. These findings characterize the ECG as abnormal, and the conduction delay could suggest underlying pathology that may require additional testing or clinical evaluation to rule out underlying cardiac conditions. This evaluation is critical to ensure comprehensive cardiovascular assessment and appropriate subsequent management.

Ground-truth Report: Sinus rhythm; Rightward axis; Right bundle branch block; Lateral T wave changes are nonspecific; Abnormal ECG.

PULSE Interpretation: The electrocardiogram (ECG) reveals a sinus rhythm, indicating a normal heart rate and rhythm. However, the ECG also exhibits several abnormalities, including a prolonged QT interval, which may indicate a prolonged repolarization phase, and a right axis deviation, suggesting an abnormal orientation of the heart's electrical axis. Additionally, the ECG shows signs of right ventricular hypertrophy, low QRS voltages in the limb leads, and extensive ST-T changes, which may be attributed to either hypertrophy or ischemia, leading to an overall abnormal ECG.

Figure 9. Cardiologist Evaluations. Blue: Findings exceeding expert expectations. Yellow: Findings with differing expert opinions.