

Influence of voltmeter input impedance on quantum Hall effect measurements

F. Fischer and M. Grayson

Walter Schottky Institut, Technische Universität München, 85748, Garching, Germany

We report on the influence of voltmeters on measurements of the longitudinal resistance in the quantum Hall effect regime. We show that for typical input resistances for standard digital lock-in amplifiers the longitudinal resistance can show a non-zero minimum which might be mistaken for parallel conduction in the doping layer. This residual resistance can be calculated with $R_{xx}^{res} = R_{xy}^2/R_{in} + j\omega CR_{xy}^2$, where R_{in} is the input resistance of the voltmeter and C the measurement capacitance. In contrast to a real parallel conduction the effect disappears when either the current source and ground contact are swapped or the polarity of the B-field is changed. We discuss the influence of input capacitances and stray capacitances on the measurement. The data demonstrates the influence of the voltmeter input impedance on the longitudinal resistance measurement.

I. INTRODUCTION

The quantum Hall effect (QHE) [1, 2] of a two-dimensional system (2DS) in a perpendicular magnetic field B is characterized by the quantization of the transversal resistance to $R_{xy} = h/\nu e^2$ ($\nu = 1, 2, \dots$) coincident with minima in the longitudinal resistance R_{xx} in the quantum Hall regime. R_{xx} vanishes to $R_{xx}=0$ in high mobility samples and sufficiently low temperatures. A similar behaviour occurs for the fractional quantum Hall effect (FQHE) [3, 4] for odd denominator fractions ($\nu = \frac{1}{3}, \frac{2}{5}, \frac{2}{3}, \dots$). A flat low temperature non-zero minimum is normally attributed to a parallel conducting dopant layer [5, 6, 7].

In the following article we will show measurements where parallel conduction seems to be evident but in reality is due to a measurement error, because of a finite input impedance of the measuring voltmeter.

II. ANALYSIS

The standard method to measure the quantum Hall effect is a four-point measurement with a Hall-bar indicated in Fig. 1. A typical Hall-bar has a current input I_i and a current output I_o . There are 4 voltage contacts V_1, V_2, V_3 and V_4 to measure the electrical potential at fixed points along the sample. For an ideally homogeneous and symmetric sample with $I_i = I_o = I$, R_{xy} is given by $R_{xy} = |(V_1 - V_3)|/I = |(V_2 - V_4)|/I$ and $R_{xx} = |(V_1 - V_2)|/I = |(V_3 - V_4)|/I$. In the common QHE edge state picture [8] the current flows only along a small strip at the edge of the sample as indicated in Fig. (1) for electrons. Following the Landauer-Büttiker formalism for QHE edge channels [9, 10] the total current along the edge is given by

$$\nu \frac{e}{h} \mu_k = \nu \frac{e^2}{h} V_k \quad (1)$$

with ν being the filling factor, μ_k the chemical potential of the edge channel with index k and V_k the corresponding voltage. The chemical potential of each edge channel is determined by the voltage of the injecting contact. In

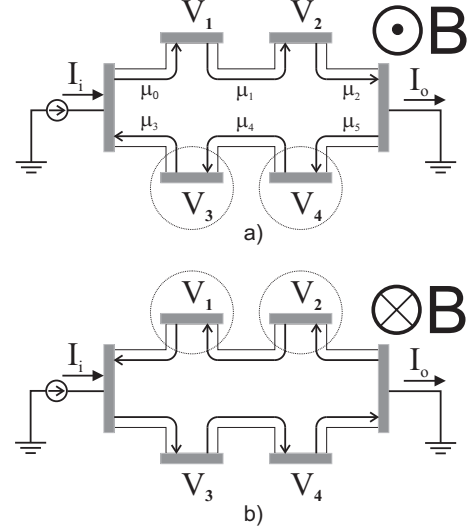


FIG. 1: Typical example of a Hall-bar geometry. I_i and I_o indicate the input and output currents. Figure a: For the indicated B-field pointing out of the plane the arrow-lines represent electrons in the QHE regime travelling clockwise along the edge having a certain chemical potential μ_i ($i = 1, \dots, 6$). For hole edge channels one can simply flip the indicated B-field polarity. Figure b: For the B-field pointing into the plane electrons travel counterclockwise.

The circled contacts show proper $R_{xx} = 0$ minima, whereas the voltage contacts opposite show residual R_{xx}^{res} in the minima described in Eqs. (9a) and (9b)

contrast to the conventional current, the direction of the edge channel flow equals the direction of the carrier movement and is determined by the sign of the charge and the polarity of the magnetic field. The following discussion assumes edge currents running clockwise, as in the case of electrons and a magnetic field pointing out of the Hall-bar (see Fig. 1).

If a current is driven through the Hall-bar but the voltage contacts are not connected then $\mu_0 = \mu_1 = \mu_2$ is the chemical potential of the current source and $\mu_3 = \mu_4 = \mu_5$ is the chemical potential of the ground. The same is applicable if the measurement devices connected to the voltage contacts V_i do not draw any current, and the

edge current flowing into a voltage contact equals the edge current flowing out. In the QHE regime R_{xx} is then given by

$$R_{xx}^{1-2} = \frac{V_1 - V_2}{I} = e \frac{\mu_1 - \mu_2}{I} = 0 \quad \text{and} \quad (2a)$$

$$R_{xx}^{3-4} = \frac{V_3 - V_4}{I} = e \frac{\mu_3 - \mu_4}{I} = 0. \quad (2b)$$

When the voltage contacts 1 and 2 are connected to a voltmeter with a finite input-impedance $Z_{in} = (R_{in}^{-1} + j\omega C)^{-1}$ a current I^{lost} flows through the voltmeter to ground (see Fig. 2). Since the potential of a contact is

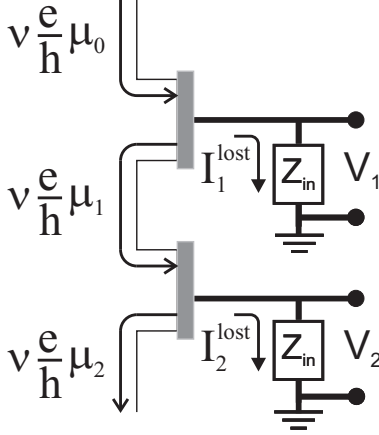


FIG. 2: Circuit diagram for a measurement of the longitudinal resistance. A voltmeter with input impedance Z_{in} measures the differential voltage $V_{xx} = V_1 - V_2$.

given by the chemical potential of the reduced outflowing edge current, this leads to a change of the chemical potential of the edge current at every voltage contact. One would then measure a finite residual resistance $R_{xx}^{res} \neq 0$ in a QHE minimum of R_{xx} . The measured voltages V_{xx}^{1-2} and V_{xx}^{3-4} are given by

$$V_{xx}^{1-2} = e(\mu_1 - \mu_2) = \frac{h}{\nu e^2} I_2^{lost} = \frac{h}{\nu e^2} \frac{1}{Z_{in}} e \mu_2 \quad (3a)$$

$$V_{xx}^{3-4} = \frac{h}{\nu e^2} \frac{1}{Z_{in}} e \mu_3 \quad (3b)$$

with $(\nu e/h)\mu_1 = (\nu e/h)\mu_2 + I_2^{lost}$ required by current conservation and $I_2^{lost} = V_2/Z_{in}$.

Measurements of the QHE are normally performed with one of the current contacts set to ground and driving a constant DC or constant AC sinusoidal current through the sample. In the QHE regime the voltage of the non-grounded current contact is given by $V_{xy} = R_{xy}I = h/\nu e^2 I$ with I being the driving current. By eliminating μ_1 in the set of equations:

$$e\mu_0 = I_1^{lost} R_{xy} + e\mu_1 \quad (4a)$$

$$e\mu_1 = I_2^{lost} R_{xy} + e\mu_2 \quad (4b)$$

one deduces

$$\mu_2 = \mu_0 \left(\frac{R_{xy}^2}{Z_{in}^2} + 2 \frac{R_{xy}}{Z_{in}} + 1 \right)^{-1} \quad (5)$$

$$(6)$$

respectively

$$\mu_3 = \mu_5 \left(\frac{R_{xy}^2}{Z_{in}^2} + 2 \frac{R_{xy}}{Z_{in}} + 1 \right)^{-1} \quad (7)$$

In standard voltmeters $R_{xy} \ll Z_{in}$, since $R_{xy} \sim h/e^2 = 25 \text{ k}\Omega$, and $Z_{in} \geq 10 \text{ M}\Omega$. The above expression reduces to $\mu_2 \approx \mu_0 = V_{xy}/e$ and $\mu_3 = \mu_5 = 0$, thus

$$V_{xx}^{1-2} = R_{xy}^2 \frac{1}{Z_{in}} I = \left(\frac{h}{\nu e^2} \right)^2 \frac{1}{Z_{in}} I$$

$$= R_{xy}^2 \left(\frac{1}{R_{in}} + j\omega C \right) I \quad (8a)$$

$$V_{xx}^{3-4} = 0 \quad (8b)$$

the real and imaginary components of V_{xx}^{1-2} are then given by

$$Re(V_{xx}^{1-2}) = \frac{R_{xy}^2}{R_{in}} I$$

$$\Rightarrow Re(R_{xx}^{1-2}) = \frac{Re(V_{xx})}{I} = \frac{R_{xy}^2}{R_{in}} \quad (9a)$$

$$Im(V_{xx}^{1-2}) = R_{xy}^2 \omega C I$$

$$\Rightarrow Im(R_{xx}^{1-2}) = \frac{Im(V_{xx})}{I} = R_{xy}^2 \omega C \quad (9b)$$

The results (8a) and (8b) show that on one side of the Hall-bar there is a spurious measured resistance and on the other side there is none. Whether the residual resistance $R_{xx}^{res} = R_{xy}^2/R_{in} + j\omega C R_{xy}^2$ is measured is determined by the chemical potential sourcing the edge currents. If an effect of the input impedance on the measurement is suspected, one can investigate that by either swapping the current source contact with the ground contact, by changing the polarity of the magnetic field or by using voltage contacts on the other side of sample. It is worth reminding the reader that a parallel conducting layer would also lead to a non-zero R_{xx} minima increasing with the magnetic field [5], but *would be seen at all R_{xx} contact pairs, would not disappear upon reversing the B -field and is normally not quadratic in $1/\nu$.*

III. MEASUREMENT

A voltage measurement can either be done DC or AC, the latter normally performed with lock-in amplifiers. Table III gives an overview over the input impedances

manufacturer	model	mode	R	C	residual $\text{Re}(R_{xx})$
Keithley	1801	DC	$> 1 \text{ G}\Omega$	-	$< 0.66 \text{ }\Omega/\nu^2$
Agilent	34420A	DC	$> 10 \text{ G}\Omega$	$< 3.6 \text{ nF}$	$< 0.06 \text{ }\Omega/\nu^2$
Signal Recovery (EG&G)	7260	digital AC	$10 \text{ M}\Omega$	30 pF	$66.6 \text{ }\Omega/\nu^2$
Signal Recovery (EG&G)	210	analog AC	$100 \text{ M}\Omega$	25 pF	$6.66 \text{ }\Omega/\nu^2$
Stanford Research	810	digital AC	$10 \text{ M}\Omega$	25 pF	$66.6 \text{ }\Omega/\nu^2$
Stanford Research	510	analog AC	$100 \text{ M}\Omega$	25 pF	$6.66 \text{ }\Omega/\nu^2$

TABLE I: Input impedances of different commercially available voltmeters (data taken from the product specifications).

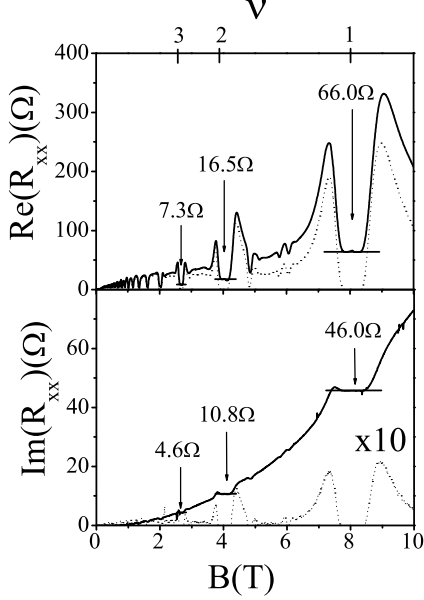


FIG. 3: Real part (top figure) and imaginary part (bottom figure) of R_{xx} corresponding to R_{xx}^{1-2} (solid line) and R_{xx}^{3-4} (dotted line) as measured in Fig. 1a)

$Z_{in} = (R_{in}^{-1} + j\omega C)^{-1}$ of commercially available voltmeters. The data in Fig. III shows the real and imaginary parts of R_{xx}^{1-2} (solid line) and R_{xx}^{3-4} (dotted line) for a 2DEG in an AlGaAs/GaAs heterostructure measured as in Fig. 1a) with B coming out of the sample. The data was taken with an EG&G 7260 lock-in amplifier in a vander Pauw geometry instead of a Hall-Bar geometry at an excitation frequency of $\omega = 2\pi f = 17 \text{ Hz}$. From the finite R_{xx}^{1-2} resistance minima one can directly deduce the input impedance of the lock-in amplifier using eq. (8a) $Z_{in} = 10.1 \text{ M}\Omega + j\omega 4 \text{ nF}$. While the real part is nearly exactly the value given in Table III, the capacitance from the imaginary part is about a factor of $\times 20$ larger than the lock-in specifications. This is because the total capacitance $C_{total} = 4 \text{ nF}$ is $C_{total} = C_{lock-in} + C_{setup}$, where $C_{lock-in} = 30 \text{ pF}$ is the capacitance of the lock-in given in Table III and $C_{setup} = 3970 \text{ pF}$ the stray capac-

itance of the rest of the measurement setup principally the probe wiring. Since $\text{Im}(R_{xx}) \sim \omega$ scales linearly with frequency, any structures in R_{xx} that change frequency can be assumed to derive from the total measurement capacitance.

Standard lock-in amplifiers are equipped with optional internal line filters at frequencies 50 Hz and 100Hz (or 60 Hz and 120 Hz depending on the line standard). If activated, these filters can even at low lock-in frequencies (1-10 Hz), induce a phase shift of the residual resistance signal $R_{xx}^{res} = R_{xx}^0 e^{i\phi}$ by $\Delta\phi$ to $R_{xx} = R_{xx}^0 e^{i(\phi+\Delta\phi)}$ and rotating some of the imaginary part into the real. Therefore a large capacitive signal can be rotated into the real R_{xx} trace and vice versa. Because this can lead to a deviation of the expected residual resistance values it is recommended not to use internal filters unless they are needed, and then to carefully calibrate out the induced phase shift.

IV. CONCLUSION

We conclude that special care has to be taken while setting up a measurement of the longitudinal resistance. We showed that although the value of the voltmeter input resistance of $R = 10 \text{ M}\Omega$ seems to be sufficiently high, it can highly influence a measurement of R_{xx} by showing a residual longitudinal resistance of $R_{xx}^{res} = R_{xy}^2/R_{in} + j\omega C R_{xy}^2$ and lead to the false assumption of a parallel channel. If a finite R_{xx} minimum drops to zero with any of the following 3 tests, the input resistance is responsible for the residual R_{xx} .

1. Switching the polarity of the magnetic field.
2. Swapping the current source contact.
3. Using voltage contacts on the other side of the sample.
4. For phasing errors, one can check if changing the excitation frequency by a certain factor, changes the R_{xx} signal at the minimum.

In contrast to AC lock-in amplifiers, commercially available DC-voltmeters usually have sufficiently high input impedances of $R > 1 \text{ G}\Omega$ giving a residual $R_{xx}^{res} = 0.66 \text{ }\Omega/\nu^2$ (see Table III), a value small enough that it

does not disturb a measurement. The analog lock-in amplifiers have a residual $Re(R_{xx}^{res}) = 6.6 \, \Omega/\nu^2$, a value which could already be seen in measurements especially in the FQHE-regime. The digital lock-in amplifiers have residual input resistances of $Re(R_{xx}^{res}) = 66.6 \, \Omega/\nu^2$, a value which can clearly be seen in measurements. It is advisable to make use of pre-amplifiers with higher input impedances. The total capacitance of the measurement setup C_{setup} (typically several hundred pF or more) must be considered when checking for $Im(R_{xx}^{res})$ effects. Built-in line filters in dilution refrigerator systems to filter out radio-frequencies have a larger contribution to C_{setup} .

In summary, the ideal setup has a very high input impedance preamplifier to the lock-in with a low input

capacitance and small stray capacitances. As a rule of thumb, one should *always* know the chirality of the edge states and measure R_{xx} at the ground-potential side of the sample.

Acknowledgments

This work was supported financially by Deutsche Forschungsgemeinschaft via Schwerpunktprogramm Quantum-Hall-Systeme. M. Grayson would like to thank A. v. Humboldt Foundation for support.

-
- [1] K. v. Klitzing, G. Dorda, and M. Pepper, Phys. Rev. Lett. **45**, 494 (1980).
 - [2] K. von Klitzing, Rev. Mod. Phys. **58**, 519 (1986).
 - [3] H. L. Stormer, A. Chang, D. C. Tsui, J. C. M. Hwang, A. C. Gossard, and W. Wiegmann, Phys. Rev. Lett. **50**, 1953 (1983).
 - [4] H. L. Stormer, D. C. Tsui, and A. C. Gossard, Rev. Mod. Phys. **71**, S298 (1999).
 - [5] M. Grayson and F. Fischer, submitted to J. Appl. Phys. in tandem with this paper
 - [6] M. van der Burgt, V. C. Karavolas, F. M. Peeters, J. Singleton, R. J. Nicholas, F. Herlach, J. J. Harris, M. V. Hove, and G. Borghs, Phys. Rev. B **52**, 12218 (1995).
 - [7] S. Contreras, W. Knap, E. Frayssinet, M. L. Sadowski, M. Goiran, and M. Shur, J. of Appl. Phys. **89**, 1251 (2001).
 - [8] B. I. Halperin, Phys. Rev. B **25**, 2185 (1982).
 - [9] J. K. Jain and S. A. Kivelson, Phys. Rev. B **37**, 4276 (1988).
 - [10] M. Buttiker, Phys. Rev. B **38**, 9375 (1988).